

AR 226 - 3682

INTERIM REPORT #23 Analysis of Sediment Samples

STUDY TITLE

Analysis of Perfluorooctanoic Acid (PFOA) in Water, Soil, Sediment, Fish, Clams,
Vegetation, Small Mammal Liver and Small Mammal Serum Using LC/MS/MS for the
3M Decatur Monitoring Program

DATA REQUIREMENTS

EPA TSCA Good Laboratory Practice Standards 40 CFR 792

STUDY DIRECTOR

Jaisimha Kesari P.E., DEE
Weston Solutions, Inc.
1400 Weston Way
West Chester, PA 19380
Phone: 610-701-3761

INTERIM REPORT COMPLETION DATE

May 2, 2006

PERFORMING LABORATORY

Exygen Research
3058 Research Drive
State College, PA 16801
Phone: 814-272-1039

STUDY SPONSOR

3M Company
3M Building 0236-01-B-10
St. Paul, MN 55144
Phone: 651-733-6374

PROJECT

Protocol Number: P0000760
Exygen Study Number: P0000760

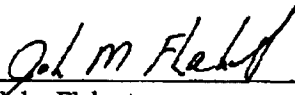
Total Pages: 110

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

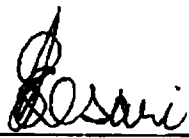
GOOD LABORATORY PRACTICE COMPLIANCE STATEMENT

Exygen Study Number P0000760, entitled "Analysis of Perfluorooctanoic Acid (PFOA) in Water, Soil, Sediment, Fish, Clams, Vegetation, Small Mammal Liver and Small Mammal Serum Using LC/MS/MS for the 3M Decatur Monitoring Program," conducted for 3M Company, is being performed in compliance with EPA TSCA Good Laboratory Practice Standards 40 CFR 792 by Exygen Research, with one exception. The laboratory control and control spikes for the direct injection and alternate SPE methods and the spiking procedures for the surrogate were not recorded promptly at the time of extraction.



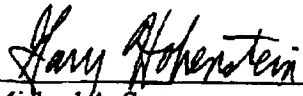
John Flaherty
Principal Investigator
Exygen Research

5/2/06
Date



Jaisimha Kesari P.E., DEE
Study Director
Weston Solutions, Inc.

5/3/06
Date



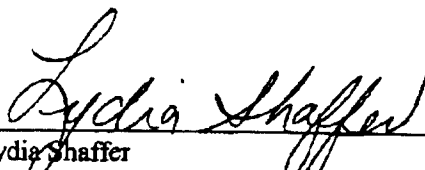
Michael A. Santoro
Sponsor Representative
3M Company

5/4/06
Date

QUALITY ASSURANCE STATEMENT

Exygen Research's Quality Assurance Unit reviewed Exygen Study Number P0000760, entitled, "Analysis of Perfluorooctanoic Acid (PFOA) in Water, Soil, Sediment, Fish, Clams, Vegetation, Small Mammal Liver and Small Mammal Serum Using LC/MS/MS for the 3M Decatur Monitoring Program". All reviewed phases¹ were inspected for conduct according to Exygen Research's Standard Operating Procedures, the Study Protocol, and all applicable Good Laboratory Practice Standards. All findings were reported to the Exygen Principal Investigator and Management and to the Study Director.

<u>Phase</u>	<u>Date Inspected</u>	<u>Date Reported to Principal Investigator</u>	<u>Date Reported to Exygen Management</u>	<u>Date Reported to Study Director</u>
31. Final Raw Data Review and Interim Analytical Report Review	5/2/06	5/2/06	5/2/06	5/2/06


 Lydia Shaffer
 Technical Lead, Quality Assurance Unit

05/02/06
 Date

¹Note: All in-lab inspections will be documented in the QA statement for the final analytical report at the conclusion of the study. This QA statement involves only the review of the interim report and associated raw data.

Interim Report #23 – Analysis of Sediment Samples

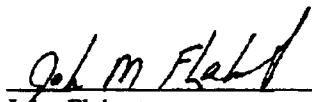
Exygen Study No.: P0000760

CERTIFICATION OF AUTHENTICITY

This interim report, for Exygen Study Number P0000760, is a true and complete representation of the raw data.

Submitted by: Exygen Research
3058 Research Drive
State College, PA 16801
(814) 272-1039

Principal Investigator, Exygen:

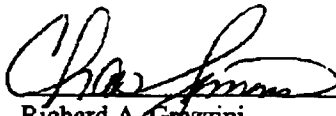


John Flaherty
Vice President
Exygen Research

5/2/06

Date

Exygen Research Facility Management:

 ^{FDR}

Richard A. Grazzini
President
Exygen Research

5/2/06

Date

Study Director, Weston Solutions, Inc.

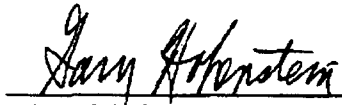


Jaisimha Kesari P.E., DEE
Weston Solutions, Inc.

5/3/06

Date

Sponsor Representative, 3M Company:



Michael A. Santoro
Director of Regulatory Affairs

5/4/06

Date

STUDY IDENTIFICATION

Analysis of Perfluorooctanoic Acid (PFOA) in Water, Soil, Sediment, Fish, Clams,
Vegetation, Small Mammal Liver and Small Mammal Serum Using LC/MS/MS for the
3M Decatur Monitoring Program

PROTOCOL NUMBER: P0000760

EXYGEN STUDY NUMBER: P0000760

TYPE OF STUDY: Residue

SAMPLE MATRIX: Sediment

TEST SUBSTANCE: Perfluorooctanoic acid (PFOA)

SPONSOR: 3M Company
3M Building 0236-01-B-10
St. Paul, MN 55144

STUDY DIRECTOR: Jaisimha Kesari P.E., DEE
Weston Solutions, Inc.
1400 Weston Way
West Chester, PA 19380

STUDY MONITOR: Michael A. Santoro
3M Company
3M Building 0236-01-B-10
St. Paul, MN 55144

PERFORMING LABORATORY: Exygen Research
3058 Research Drive
State College, PA 16801

**ANALYTICAL PHASE
TIMETABLE:** Study Initiation Date: 11/05/04
Interim Analytical Start Date: 04/27/06
Interim Analytical Termination Date: 04/30/06
Interim Report Completion Date: 05/02/06

PROJECT PERSONNEL

The Study Director for this project is Jaisimha Kesari at Weston Solutions, Inc. The following personnel from Exygen Research were associated with various phases of this interim portion of the study:

<u>Name</u>	<u>Title</u>
John Flaherty	Vice President
Karen Risha	Scientist
Paul Connolly	Technical Lead-LC/MS
Christine Edwards	Technician
Mark Ammerman	Sample Custodian
Amy Sheehan	Associate Scientist
Brittany Kravets	Technician
Scott Crain	Technician
Frances Crespi	Technician
Devon Ingram	Technician
Christa Gallant	Technician

TABLE OF CONTENTS

	<u>Page</u>
TITLE PAGE	1
GOOD LABORATORY PRACTICE COMPLIANCE STATEMENT	2
QUALITY ASSURANCE STATEMENT	3
CERTIFICATION OF AUTHENTICITY	4
STUDY IDENTIFICATION.....	5
PROJECT PERSONNEL.....	6
TABLE OF CONTENTS.....	7
LIST OF TABLES	8
LIST OF FIGURES.....	9
LIST OF APPENDICES	10
1.0 SUMMARY	11
2.0 OBJECTIVE	11
3.0 INTRODUCTION.....	12
4.0 ANALYTICAL TEST SAMPLES.....	12
5.0 REFERENCE MATERIAL	12
6.0 DESCRIPTION OF ANALYTICAL METHOD	13
6.1 Extraction Procedure For Sediment.....	13
6.2 Direct Inject Extraction Procedure For Sediment.....	14
6.3 Additional SPE Extraction Procedure For Sediment.....	14
6.4 Percent Solids Procedure For Sediment	14
6.5 Preparation of Standards and Fortification Solutions.....	14
6.6 Chromatography	16
6.7 Instrument Sensitivity.....	16
6.8 Description of LC/MS/MS Instrument and Operating Conditions.....	16
6.9 Quantitation and Example Calculation.....	17
7.0 EXPERIMENTAL DESIGN	19
8.0 RESULTS	19
9.0 CONCLUSIONS	20
10.0 RETENTION OF DATA AND SAMPLES	20

LIST OF TABLES

	<u>Page</u>
Table I. Summary of PFOA in Sediment Samples.....	22
Table II. Summary of PFOA in Direct Injection Sediment Samples	23
Table III. Summary of PFOA in SPE Extracted Sediment Samples	24
Table IV. Matrix Spike Recovery of PFOA in Sediment Samples.....	25
Table V. Matrix Spike Recovery of PFOA in Direct Injection Sediment Samples.....	26
Table VI. Matrix Spike Recovery of PFOA in SPE Extracted Sediment Samples.....	27
Table VII. Surrogate Spike Recovery of ¹³ C PFOA in Sediment Samples	28
Table VIII. Surrogate Spike Recovery of ¹³ C PFOA in Direct Injection Sediment Samples	29
Table IX. Surrogate Spike Recovery of ¹³ C PFOA in SPE Extracted Sediment Samples	30
Table X. Total Percent Solids in Sediment Samples	31

LIST OF FIGURES

	<u>Page</u>
Figure 1. Typical Calibration Curve for PFOA in Methanol	33
Figure 2. Extracted Standards of PFOA in Methanol, 25 ng/L and 50 ng/L, Respectively.....	34
Figure 3. PFOA in Reagent Blank, 0.4 ng/mL Fortified Reagent Spk A, and 4.0 ng/mL Fortified Reagent Spk B, Respectively.....	35
Figure 4. Chromatogram Representing a Sediment Sample Analyzed for PFOA (Exygen ID: C0172892, Data Set: 042806D).....	36
Figure 5. Typical Calibration Curve for ¹³ C PFOA in Methanol.....	37
Figure 6. Extracted Standards of ¹³ C PFOA in Methanol, 25 ng/L and 50 ng/L, Respectively.....	38
Figure 7. ¹³ C PFOA in Reagent Blank, 0.4 ng/mL Fortified Reagent Spk A, and 4.0 ng/mL Fortified Reagent Spk B, Respectively.....	39
Figure 8. Chromatogram Representing a Sediment Sample Analyzed for ¹³ C PFOA (Exygen ID: C0172892, Data Set: 042806D)	40

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

LIST OF APPENDICES**Page**

Appendix A	Study Protocol P0000760 (Exygen Study No. P0000760) with Analytical Method V0001782: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS", and Protocol Amendments	41
-------------------	---	----

1.0 SUMMARY

Exygen Research extracted and analyzed sediment samples for the determination of perfluorooctanoic acid (PFOA) according to Exygen Method V0001782 (Appendix A).

The limit of quantitation for PFOA in sediment was 2.0 ng/g.

The standard sediment method (V0001782) was originally used to extract all of the samples. Use of this method yielded poor spiking recoveries for PFOA. Because of this, a direct injection method (Section 6.2) was used, which also yielded poor recoveries. Finally, an alternative SPE column extraction (Section 6.3) was used to analyze for PFOA and ¹³C PFOA.

Analytical results for the analysis of PFOA in sediment samples extracted according to the method based on dry weight are summarized in Table I. Analytical results for the analysis of PFOA in sediment samples by direct injection, based on dry weight, are summarized in Table II. Analytical results for the analysis of PFOA in sediment samples by alternative SPE extraction, based on dry weight, are summarized in Table III. The overall average percent recovery $\pm$ standard deviation for PFOA in sediment samples based on wet weight was $73 \pm 16\%$. Fortification recoveries for the analysis of PFOA in sediment samples are summarized in Table IV. Fortification recoveries for the analysis of PFOA in direct injection sediment samples are summarized in Table V. Fortification recoveries for the analysis of PFOA in SPE extracted sediment samples are summarized in Table VI. The average percent recovery $\pm$ standard deviation for ¹³C PFOA in sediment samples based on wet weight was $97 \pm 27\%$. Fortification recoveries for ¹³C PFOA in sediment samples based on wet weight are summarized in Table VII. Fortification recoveries for ¹³C PFOA in direct injection sediment samples based on wet weight are summarized in Table VIII. Fortification recoveries for ¹³C PFOA in SPE extracted sediment samples based on wet weight are summarized in Table IX.

The total percent solids for sediment samples are detailed in Table X.

2.0 OBJECTIVE

The objective of the analytical part of this study was to determine levels of perfluorooctanoic acid (PFOA) in sediment according to Protocol P0000760 (Appendix A).

3.0 INTRODUCTION

This report details the results of the analysis for the determination of PFOA in sediment using the analytical methods entitled, "V0001782: Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS."

The study was initiated on November 05, 2004, when the study director signed protocol number P0000760. The analytical start date for this interim report was April 27, 2006 and the analytical termination date for this interim report was April 30, 2006.

4.0 ANALYTICAL TEST SAMPLES

Three sediment samples (Exygen ID: C0172892–C0172894) were received on wet ice on April 22, 2006 from Tim Frinak at Weston Solutions, Inc. The samples were logged in by Exygen personnel and placed in refrigerated storage.

Sample login and chain of custody information is located in the raw data package associated with this interim report. Storage records will be kept at Exygen Research.

5.0 REFERENCE MATERIAL

The analytical standard, PFOA, was purchased from Oakwood Products, Inc. and was received at Exygen on March 17, 2005. The surrogate spiking standard, ¹³C labeled perfluorooctanoic acid (¹³C PFOA SP0004184), was received at Exygen on April 15, 2004 from the 3M Company.

The available information for the reference materials is listed below. PFOA was stored ambient and ¹³C PFOA was stored frozen.

<u>Compound</u>	<u>Exygen Inventory No.</u>	<u>Lot #</u>	<u>Purity (%)</u>	<u>Expiration Date</u>
PFOA	SP0005444	203	99	03/31/07
¹³ C PFOA	SP0004184	3507-195	97	03/29/09

The molecular structure of PFOA is given below:

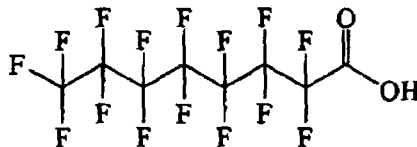
PFOA

Chemical Name: Perfluorooctanoic acid

Molecular Weight: 414

Transitions Monitored: 413 → 369 (for quantification) and
413 → 219 (for confirmation)

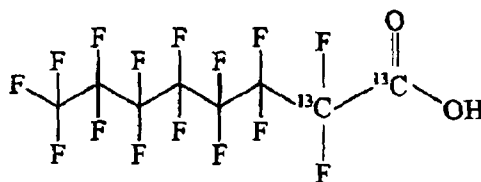
Structure:

**¹³C PFOA**Chemical Name: 1,2-¹³C perfluorooctanoic acid

Molecular Weight: 416

Transition Monitored: 415 → 370

Structure:

**6.0 DESCRIPTION OF ANALYTICAL METHOD**

The analytical method "V0001782: Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS" was used for the sediment samples in this study.

6.1 Extraction Procedure For Sediment

Before the samples were weighed for the extraction, they were mixed thoroughly by vigorously shaking the container. A 5 gram portion of sediment was weighed into a fifty milliliter centrifuge tube for the extraction. After fortification of appropriate samples, 35 mL of 1% acetic acid in water was added to the samples. The samples were vortexed and allowed to shake on a wrist action shaker for ~60 minutes. The samples were centrifuged for ~20 minutes at ~3000 rpm. The supernatant was then loaded onto a C₁₈ SPE cartridge conditioned with 10 mL of methanol and 20 mL of water. The eluate was discarded. Twenty milliliters of methanol was added to the sediment samples left in the centrifuge tube. The samples were vortexed and allowed to shake on a wrist action shaker for another 30 minutes. The samples were centrifuged again for ~20 minutes at ~3000 rpm. The supernatant was then loaded onto the same C₁₈ SPE cartridge. The eluate was collected into a 500 mL Nalgene Bottle. The column was washed with 4 mL of methanol. The wash was collected in the same bottle as the eluate. Approximately two

hundred milliliters of water was added to the bottles. The samples were mixed by shaking and loaded onto another C₁₈ SPE cartridge conditioned with 10 mL of methanol and 20 mL of water. The eluate was discarded. Approximately five milliliters of methanol was added to the cartridge. Five milliliters of eluate was collected into a graduated 15 mL polypropylene centrifuge tube. Each sample was analyzed by LC/MS/MS electrospray.

6.2 Direct Inject Extraction Procedure For Sediment

Before the samples were weighed for the extraction, they were mixed thoroughly by vigorously shaking the container. A one-gram portion of sediment was weighed into a 15-milliliter centrifuge tube for the extraction. Ten milliliters of 1% acetic acid in methanol was added to each sample. The samples were then shaken by hand, vortexed, and sonicated for thirty minutes. The samples were then centrifuged for ~10 minutes at ~3000 rpm. Each sample was analyzed by LC/MS/MS electrospray.

6.3 Additional SPE Extraction Procedure For Sediment

The samples that were prepared in 1% acetic acid for the direct injection method were used for this extraction. Five milliliters of each sample was aliquoted into a 50-mL polypropylene centrifuge tube and the volume was taken to 40 mL with water. The samples were then centrifuged for ~10 minutes at ~3000 rpm. The supernatant was then loaded onto a C₁₈ SPE cartridge conditioned with 10 mL of methanol and 5 mL of water. The eluate was discarded. Approximately five milliliters of methanol was added to the cartridge. Five milliliters of eluate was collected into a graduated 15 mL polypropylene centrifuge tube. Each sample was analyzed by LC/MS/MS electrospray.

6.4 Percent Solids Procedure For Sediment

Percent solids were determined using the procedure indicated in Exygen method V0000427. Approximately 20 grams of sample was weighed into a pan. The weight of the sample plus the pan was recorded. The samples were then dried in an oven overnight at 104 ± 2 °C. Then the samples were transferred to a dessicator and allowed to cool for ~15 minutes. Each sample was then weighed again, including the weight of the pan. The percent solid for each sample was then calculated. See Table X for percent solids results.

6.5 Preparation of Standards and Fortification Solutions

A mixed stock standard solution of PFOA and ¹³C PFOA was prepared at a concentration of 1000 µg/mL by dissolving 100 mg of each of the standards (corrected for purity and

salt content) in 100 mL of methanol. From this solution, the following fortification standards were prepared:

Conc. of Stock or Fort. Solution ($\mu\text{g/mL}$) ¹	Fort Volume (mL)	Final Volume (mL)	Final Conc. of Fortification Std. ($\mu\text{g/mL}$)
1000	10	100	100
100	10	100	10
10	10	100	1.0
1.0	10	100	0.1
0.1	10	100	0.01

¹ of PFOA and ¹³C PFOA

A stock solution of only ¹³C PFOA was also prepared at a concentration of 100 $\mu\text{g/mL}$ by dissolving 10 mg the standard (corrected for purity and salt content) in 100 mL of methanol. From this solution, the following fortification standards were prepared:

Conc. of Stock or Fort. Solution ($\mu\text{g/mL}$) ¹	Fort Volume (mL)	Final Volume (mL)	Final Conc. of Fortification Std. ($\mu\text{g/mL}$)
100	10	100	10
10	10	100	1.0
1.0	10	100	0.1

¹ of ¹³C PFOA

A set of non-extracted calibration standards containing PFOA and ¹³C PFOA was prepared in methanol, as specified in Exygen method V0001782. The following concentrations were prepared

Conc. of Fort or Calibration Solution (ng/mL) ¹	Fort Volume (mL)	Final Volume (mL)	Final Conc. of Calibration Std. (ng/mL)
100	1.0	10	10
100	0.5	10	5.0
100	0.2	10	2.0
10	1.0	10	1.0
5.0	1.0	10	0.5
2.0	1.0	10	0.2

¹ of PFOA and ¹³C PFOA

The standards at final concentrations of 10, 5.0, 2.0, 1.0, 0.5 and 0.2 ng/mL were used for calibration. Extracted calibration standards were prepared following method V0001780. The stock standard solution and all fortification and calibration standard solutions were stored in a refrigerator ($4^{\circ} \pm 2^{\circ}\text{C}$) when not in use. Documentation of standard preparation is located in the raw data package associated with this interim report.

6.6 Chromatography

Quantification of PFOA was accomplished by LC/MS/MS electrospray. The retention time of PFOA was ~11.1 min. Peaks above the LOD were not detected in any of the reagent blank samples corresponding to the analyte retention time.

6.7 Instrument Sensitivity

The smallest standard amount injected during the chromatographic run had a concentration of 0.2 ng/mL of PFOA.

6.8 Description of LC/MS/MS Instrument and Operating Conditions

Instrument: API 4000 Biomolecular Mass Analyzer
Interface: Turbo Ion Spray Liquid Introduction Interface
Computer: DELL OptiPlex GX400
Software: Windows NT, Analyst 1.4.1

HPLC: Hewlett Packard (HP) Series 1100
HP Quat Pump
HP Vacuum Degasser
HP Autosampler
HP Column Oven

HPLC Column: Thermo Fluophase RP, 50 mm x 2.1 mm
Column Temp.: ~30° C
Injection Vol.: 15 µL
Mobile Phase (A): 2 mM Ammonium Acetate in water
Mobile Phase (B): Methanol

<u>Time (min)</u>	<u>% A</u>	<u>% B</u>
0.0	65	35
1.0	65	35
8.0	25	75
10.0	25	75
11.0	65	35
18.0	65	35

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Total run time: ~18 min

Flow Rate: 0.3 mL/min

Ions monitored:

<u>Analyte</u>	<u>Mode</u>	<u>Transition Monitored</u>	<u>Approximate Retention Time (min)</u>
PFOA	negative	413 → 369	~11.1 min.
PFOA Confirm Ion	negative	413 → 219	~11.1 min.
¹³ C PFOA	negative	415 → 370	~11.1 min.

6.9 Quantitation and Example Calculation

PFOA concentrations are quantitated in nanogram of PFOA per milliliter of extract and then converted to nanogram per gram. The PFOA concentrations were quantitated in nanogram of PFOA per liter for Set #3 (Set 042806D). Fifteen microliters of sample or calibration standard were injected into the LC/MS/MS. The peak area was measured and the standard curve was generated (using 1/x fit weighted linear regression) by Analyst software using seven concentrations of standards. The concentration was determined from the equations below.

Please note that the equations in this section are presented in a different order than the equations in the method (V0001782), but this has no effect on the final calculations.

Equation 1 calculated the amount of analyte found (in ng/L, based on peak area) using the standard curve (linear regression parameters) generated by the Analyst software program.

Equation 1:

$$\text{Analyte found (ng/L)} = \frac{(\text{Peak area} - \text{intercept}) \times \text{DF}}{\text{Slope}}$$

Where: DF = Dilution Factor, factor by which the final volume was diluted, if necessary.

For samples fortified with known amounts of PFOA and ¹³C PFOA prior to extraction, Equation 2 was used to convert the amount of PFOA found in ng/L to ng/g (ppb).

Equation 2:

$$\text{PFOA found (ppb)} = \frac{[\text{PFOA found (ng/L)} \times \text{SPE volume extracted}]}{[\text{sample aliquot volume} \times (\text{sample weight} / \text{direct inject sample volume})]}$$

Equation 3 was used to calculate the percent recovery.

Equation 3:

Recovery (%) =

$$\frac{(\text{analyte found (ng/g, wet weight)} - \text{analyte in control (ng/g, wet weight)})}{\text{amount added (ng/g, wet weight)}} \times 100\%$$

Note: For the PFOA recovery calculation, the "control" is the unspiked aliquot of the primary field sample.

Equation 4 was then used to calculate the amount of PFOA found in ppb based on dry weight.

Equation 4:

$$\text{PFOA found (ppb) dry weight} = \text{PFOA found (ppb)} \times [100\% / \text{total solids}(\%)]$$

An example of a calculation using an actual sample follows:

Sediment sample Exygen ID C0172892 Spk D (Set: 042806D), fortified at 4 ng/mL with PFOA where:

peak area	=	236588
intercept	=	4630
slope	=	420
dilution factor	=	100
ng/g PFOA added (fort level)	=	4000
amt in corresponding sample (ng/g)	=	586
SPE volume extracted (L)	=	0.04
sample aliquot volume (mL)	=	5
sample weight (g)	=	1
direct inject sample volume (mL)	=	10
total solids (%)	=	66.54

From equation 1:

$$\begin{aligned} \text{Analyte found (ng/L)} &= \frac{[236588 - 4630]}{420} \times 100 \\ &= 55200 \text{ ng/L} \end{aligned}$$

From equation 3:

$$\begin{aligned} \text{PFOA found (ppb)} &= \frac{(55200 \text{ ng/L} \times 0.04 \text{ L})}{[5 \text{ mL} \times (1 \text{ g} / 10 \text{ mL})]} \\ &= 4420 \text{ ppb} \end{aligned}$$

From equation 3:

$$\begin{aligned} \% \text{ Recovery} &= \frac{(4420 \text{ ng/g} - 586 \text{ ng/g})}{4000 \text{ ng/g}} \times 100\% \\ &= 96\% \end{aligned}$$

From equation 4:

$$\begin{aligned}\text{PFOA found (ppb) dry weight} &= 4416 \text{ ppb} \times (100\% / 66.54\%) \\ &= 6640 \text{ ppb}\end{aligned}$$

NOTE: This value may be slightly different than that of the raw data due to rounding.

7.0 EXPERIMENTAL DESIGN

¹³C PFOA was used as a surrogate for all the samples. ¹³C PFOA was added to the sediment samples in the laboratory after the samples were weighed for extraction. For sediment samples designated as laboratory matrix spikes, PFOA was also added at a known concentration to the samples in the laboratory after the samples were weighed for extraction.

The sediment samples were extracted in one set. The set included one reagent blank and two reagent blanks fortified at known concentrations. The set contained three sediment samples. For each sample, two laboratory matrix spikes and a duplicate were extracted.

8.0 RESULTS

Analytical results for the analysis of PFOA in sediment samples based on dry weight are summarized in Table I. Analytical results for the analysis of PFOA in sediment samples by direct injection, based on dry weight, are summarized in Table II. Analytical results for the analysis of PFOA in sediment samples by alternative SPE extraction, based on dry weight, are summarized in Table III. The overall average percent recovery $\pm$ standard deviation for PFOA in sediment samples based on wet weight was $73 \pm 16\%$. Fortification recoveries for the analysis of PFOA in sediment samples are summarized in Table IV. Fortification recoveries for the analysis of PFOA in direct injection sediment samples are summarized in Table V. Fortification recoveries for the analysis of PFOA in SPE extracted sediment samples are summarized in Table VI. The average percent recovery $\pm$ standard deviation for ¹³C PFOA in sediment samples based on wet weight was $97 \pm 27\%$. Fortification recoveries for ¹³C PFOA in sediment samples based on wet weight are summarized in Table VII. Fortification recoveries for ¹³C PFOA in direct injection sediment samples based on wet weight are summarized in Table VIII. Fortification recoveries for ¹³C PFOA in SPE extracted sediment samples based on wet weight are summarized in Table IX.

The total percent solids for sediment samples are detailed in Table X.

9.0 CONCLUSIONS

The sediment samples were successfully extracted and analyzed for PFOA using an alternative SPE method, which yielded acceptable QC data.

10.0 RETENTION OF DATA AND SAMPLES

All original paper data generated by Exygen Research that pertains to this interim report will be shipped to the study director. This does not include facility-specific raw data such as instrument or temperature logs. Exact copies of all raw data, as well as a signed copy of the interim analytical report and all original facility-specific raw data, will be retained in the Exygen Research archives for the period of time specified in EPA TSCA Good Laboratory Practice Standards 40 CFR 792.

TABLES

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Table I. Summary of PFOA in Sediment Samples

Exygen ID	Client Sample ID	Analyte Found (ppb, ng/g) Dry Weight		Assessed Accuracy (+/- %)
		C8 Acid PFOA	Perfluorooctanoic Acid	
C0172892	DAL-SD-BPP01-0-0000	NR	NR	NR
C0172892 Rep	DAL-SD-BPP01-0-0000*	NR	NR	NR
C0172893	DAL-SD-BPP02-0-0000	NR	NR	NR
C0172893 Rep	DAL-SD-BPP02-0-0000*	NR	NR	NR
C0172894	DAL-SD-BPP03-0-0000	NR	NR	NR
C0172894 Rep	DAL-SD-BPP03-0-0000*	NR	NR	NR

*Laboratory Duplicate

NR = Not reported due to quality control result failures. See Table II for re-analysis results.

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Table II. Summary of PFOA in Direct Injection Sediment Samples

Exygen ID	Client Sample ID	Analyte Found (ppb, ng/g) Dry Weight		Assessed Accuracy (+/- %)
		C8 Acid PFOA	Perfluorooctanoic Acid	
C0172892	DAL-SD-BPP01-0-0000	NR		NR
C0172892 Rep	DAL-SD-BPP01-0-0000*	NR		NR
C0172893	DAL-SD-BPP02-0-0000	NR		NR
C0172893 Rep	DAL-SD-BPP02-0-0000*	NR		NR
C0172894	DAL-SD-BPP03-0-0000	NR		NR
C0172894 Rep	DAL-SD-BPP03-0-0000*	NR		NR

*Laboratory Duplicate

NR = Not reported due to quality control result failures. See Table III for re-analysis results.

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Table III. Summary of PFOA in SPE Extracted Sediment Samples

Exygen ID	Client Sample ID	Analyte Found (ppb, ng/g) Dry Weight		Assessed Accuracy (+/- %)
		C8 Acid PFOA	Perfluorooctanoic Acid	
C0172882	DAL-SD-BPP01-0-0000	881		30
C0172892 Rep	DAL-SD-BPP01-0-0000*	945		30
C0172893	DAL-SD-BPP02-0-0000	927		30
C0172893 Rep	DAL-SD-BPP02-0-0000*	823		30
C0172894	DAL-SD-BPP03-0-0000	2390		30
C0172894 Rep	DAL-SD-BPP03-0-0000*	2380		30

*Laboratory Duplicate

Table IV. Matrix Spike Recovery of PFOA in Sediment Samples

Sample Description	Amount Spiked (ng/g)	C8 Acid PFOA		
		Amt Found in Sample (ng/g) wet weight	Amount Recovered (ng/g) wet weight	Recovery (%)
DAL-SD-BPP01-0-0000 (C0172882 Spk D, 400 ppb Spike)	400	NR	NR	NR
DAL-SD-BPP01-0-0000 (C0172882 Spk E, 4000 ppb Spike)	4000	NR	NR	NR
DAL-SD-BPP02-0-0000 (C0172883 Spk G, 400 ppb Spike)	400	NR	NR	NR
DAL-SD-BPP02-0-0000 (C0172883 Spk H, 4000 ppb Spike)	4000	NR	NR	NR
DAL-SD-BPP03-0-0000 (C0172884 Spk J, 400 ppb Spike)	400	NR	NR	NR
DAL-SD-BPP03-0-0000 (C0172884 Spk K, 4000 ppb Spike)	4000	NR	NR	NR

Average: NR
Standard Deviation: NR

*Sample residue exceeds the spiking level significantly (3x spiking level); therefore, an accurate recovery value cannot be calculated.

NR = Not reported due to quality control result failures. See Table V for re-analysis results.

Note: Since this summary table shows rounded results, recovery values may vary slightly from the values in the raw data.

Table V. Matrix Spike Recovery of PFOA in Direct Injection Sediment Samples

Sample Description	Amount Spiked (ng/g)	CB Acid PFOA		
		Amt Found in Sample (ng/g) wet weight	Amount Recovered (ng/g) wet weight	Recovery (%)
DAL-SD-BPP01-0-0000 (C0172882 Spk D, 4000 ppb Spike)	4000	NR	NR	NR
DAL-SD-BPP01-0-0000 (C0172882 Spk E, 40000 ppb Spike)	40000	NR	NR	NR
DAL-SD-BPP02-0-0000 (C0172883 Spk G, 4000 ppb Spike)	4000	NR	NR	NR
DAL-SD-BPP02-0-0000 (C0172883 Spk H, 40000 ppb Spike)	40000	NR	NR	NR
DAL-SD-BPP03-0-0000 (C0172884 Spk J, 4000 ppb Spike)	4000	NR	NR	NR
DAL-SD-BPP03-0-0000 (C0172884 Spk K, 40000 ppb Spike)	40000	NR	NR	NR

Average: NR
Standard Deviation: NR

*Sample residue exceeds the spiking level significantly (3x spiking level); therefore, an accurate recovery value cannot be calculated.

NR = Not reported due to quality control result failures. See Table VI for re-analysis results.

Note: Since this summary table shows rounded results, recovery values may vary slightly from the values in the raw data.

Table VI. Matrix Spike Recovery of PFOA in SPE Extracted Sediment Samples

Sample Description	Amount Spiked (ng/g)	Amt Found in Sample (ng/g) wet weight	CB Acid PFOA	
			Amount Recovered (ng/g) wet weight	Recovery (%)
DAL-SD-BPP01-0-0000 (C0172892 Spk D, 4000 ng/g Spike)	4000	588	4420	96
DAL-SD-BPP01-0-0000 (C0172892 Spk E, 40000 ng/g Spike)	40000	588	25800	83
DAL-SD-BPP02-0-0000 (C0172893 Spk G, 4000 ng/g Spike)	4000	504	3880	84
DAL-SD-BPP02-0-0000 (C0172893 Spk H, 40000 ng/g Spike)	40000	504	25100	61
DAL-SD-BPP03-0-0000 (C0172894 Spk J, 4000 ng/g Spike)	4000	1080	4140	77
DAL-SD-BPP03-0-0000 (C0172894 Spk K, 40000 ng/g Spike)	40000	1080	23200	55
Average:				73
Standard Deviation:				16

Note: Since this summary table shows rounded results, recovery values may vary slightly from the values in the raw data.

Table VII. Surrogate Spike Recovery of ^{13}C PFOA in Sediment Samples

Exygen ID	Sample Description	Amount Spiked (ng/g)	^{13}C -PFOA	
			Amount Recovered Wet Weight (ng/g)	Recovery (%)
C0172892	DAL-SD-BPP01-0-0000	4	NR	NR
C0172892 Rep	DAL-SD-BPP01-0-0000	4	NR	NR
C0172892 Spk C	DAL-SD-BPP01-0-0000	4	NR	NR
C0172892 Spk D	DAL-SD-BPP01-0-0000	400	NR	NR
C0172892 Spk E	DAL-SD-BPP01-0-0000	4000	NR	NR
C0172893	DAL-SD-BPP02-0-0000	4	NR	NR
C0172893 Rep	DAL-SD-BPP02-0-0000	4	NR	NR
C0172893 Spk F	DAL-SD-BPP02-0-0000	4	NR	NR
C0172893 Spk G	DAL-SD-BPP02-0-0000	400	NR	NR
C0172893 Spk H	DAL-SD-BPP02-0-0000	4000	NR	NR
C0172894	DAL-SD-BPP03-0-0000	4	NR	NR
C0172894 Rep	DAL-SD-BPP03-0-0000	4	NR	NR
C0172894 Spk I	DAL-SD-BPP03-0-0000	4	NR	NR
C0172894 Spk J	DAL-SD-BPP03-0-0000	400	NR	NR
C0172894 Spk K	DAL-SD-BPP03-0-0000	4000	NR	NR

Average: NR
Standard Deviation: NR

NR = Not reported due to quality control result failures. See Table VIII for re-analysis results.
Note: Since this summary table shows rounded results, recovery values may vary slightly from the values in the raw data.

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Table VIII. Surrogate Spike Recovery of ^{13}C PFOA in Direct Injection Sediment Samples

Exygen ID	Sample Description	Amount Spiked (ng/g)	^{13}C -PFOA	
			Amount Recovered Wet Weight (ng/g)	Recovery (%)
C0172892	DAL-SD-BPP01-0-0000	40	NR	NR
C0172892 Rep	DAL-SD-BPP01-0-0000	40	NR	NR
C0172892 Spk C	DAL-SD-BPP01-0-0000	40	NR	NR
C0172892 Spk D	DAL-SD-BPP01-0-0000	4000	NR	NR
C0172892 Spk E	DAL-SD-BPP01-0-0000	40000	NR	NR
C0172893	DAL-SD-BPP02-0-0000	40	NR	NR
C0172893 Rep	DAL-SD-BPP02-0-0000	40	NR	NR
C0172893 Spk F	DAL-SD-BPP02-0-0000	40	NR	NR
C0172893 Spk G	DAL-SD-BPP02-0-0000	4000	NR	NR
C0172893 Spk H	DAL-SD-BPP02-0-0000	40000	NR	NR
C0172894	DAL-SD-BPP03-0-0000	40	NR	NR
C0172894 Rep	DAL-SD-BPP03-0-0000	40	NR	NR
C0172894 Spk I	DAL-SD-BPP03-0-0000	40	NR	NR
C0172894 Spk J	DAL-SD-BPP03-0-0000	4000	NR	NR
C0172894 Spk K	DAL-SD-BPP03-0-0000	40000	NR	NR
Average:			NR	
Standard Deviation:			NR	

NR = Not reported due to quality control result failures. See Table IX for re-analysis results.
 Note: Since this summary table shows rounded results, recovery values may vary slightly from the values in the raw data.

Table IX. Surrogate Spike Recovery of ^{13}C PFOA in SPE Extracted Sediment Samples

Exygen ID	Sample Description	^{13}C -PFOA		Recovery (%)
		Amount Spiked (ng/g)	Amount Recovered Wet Weight (ng/g)	
C0172892	DAL-SD-BPP01-0-0000	40	52.0	130
C0172892 Rep	DAL-SD-BPP01-0-0000	40	51.0	128
C0172892 Spk D	DAL-SD-BPP01-0-0000	4000	3750	94
C0172892 Spk E	DAL-SD-BPP01-0-0000	40000	25800	65
C0172893	DAL-SD-BPP02-0-0000	40	47.8	120
C0172893 Rep	DAL-SD-BPP02-0-0000	40	51.8	129
C0172893 Spk G	DAL-SD-BPP02-0-0000	4000	3500	88
C0172893 Spk H	DAL-SD-BPP02-0-0000	40000	25200	63
C0172894	DAL-SD-BPP03-0-0000	40	45.2	113
C0172894 Rep	DAL-SD-BPP03-0-0000	40	40.3	101
C0172894 Spk J	DAL-SD-BPP03-0-0000	4000	3060	77
C0172894 Spk K	DAL-SD-BPP03-0-0000	40000	22100	55
Average:				97
Standard Deviation:				27

Note: Since this summary table shows rounded results, recovery values may vary slightly from the values in the raw data.

Table X. Total Percent Solids in Sediment Samples

Exygen ID	Sample Description	Total Solids (%)
C0172892	DAL-SD-BPP01-0-0000	66.54
C0172893	DAL-SD-BPP02-0-0000	54.37
C0172894	DAL-SD-BPP03-0-0000	44.46

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

FIGURES

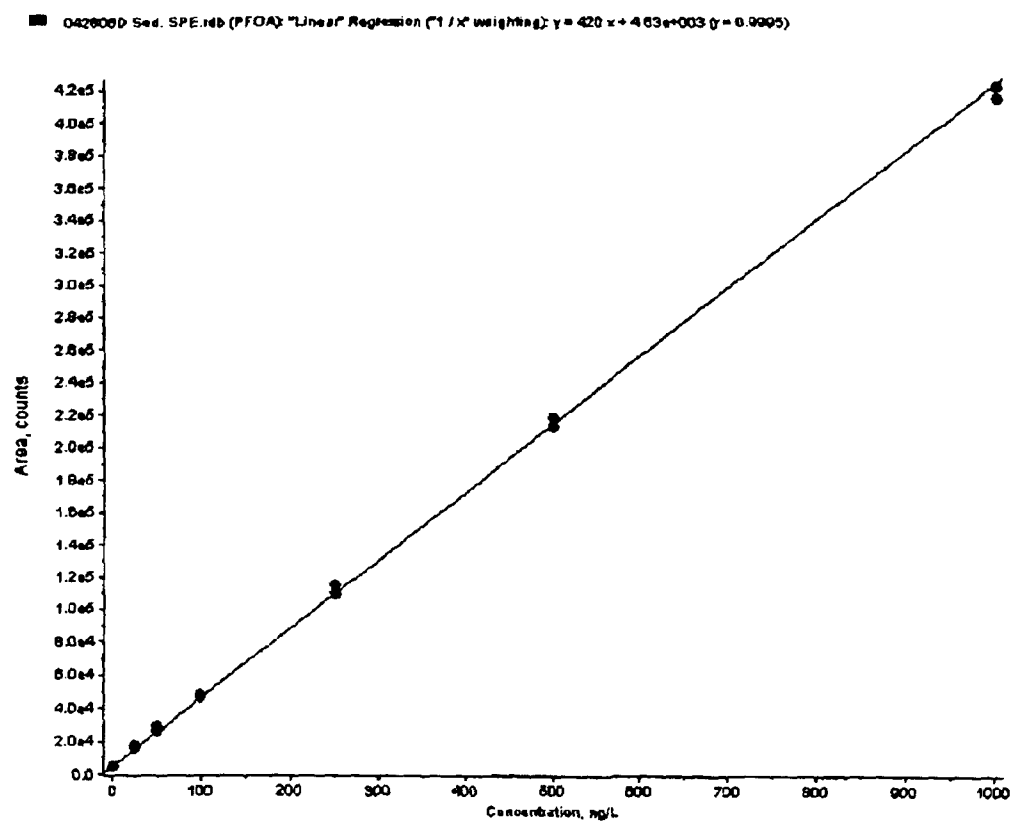
Figure 1. Typical Calibration Curve for PFOA in Methanol

Figure 2. Extracted Standards of PFOA in Methanol, 25 ng/L and 50 ng/L, Respectively

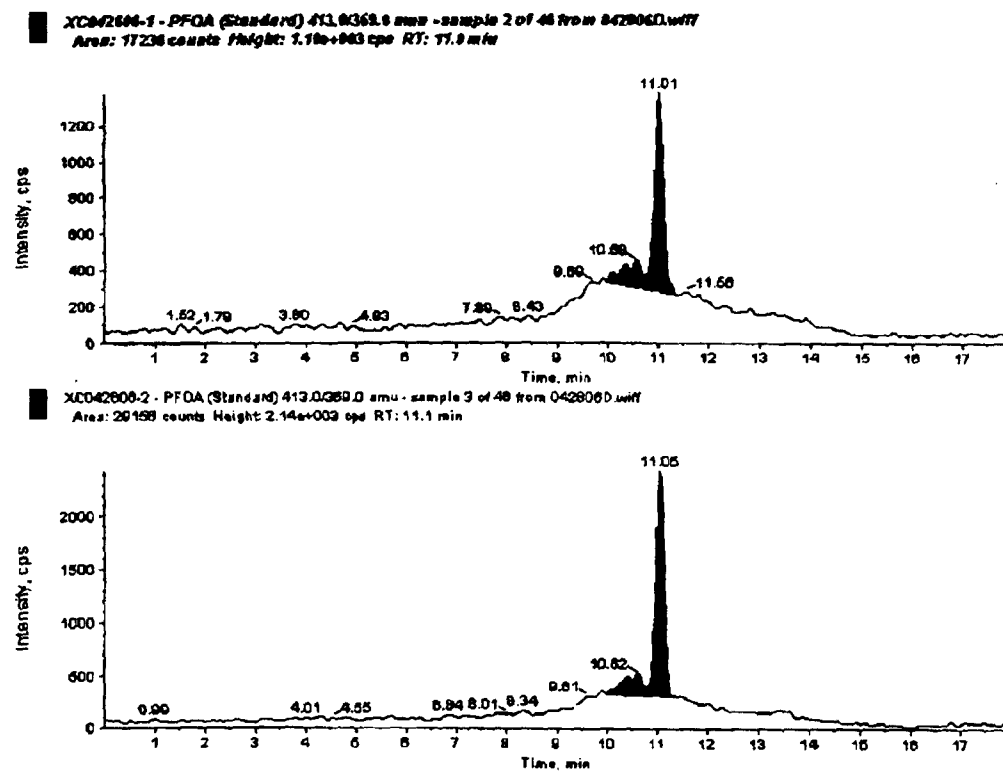


Figure 3. PFOA in Reagent Blank, 0.4 ng/mL Fortified Reagent Spk A, and 4.0 ng/mL Fortified Reagent Spk B, Respectively

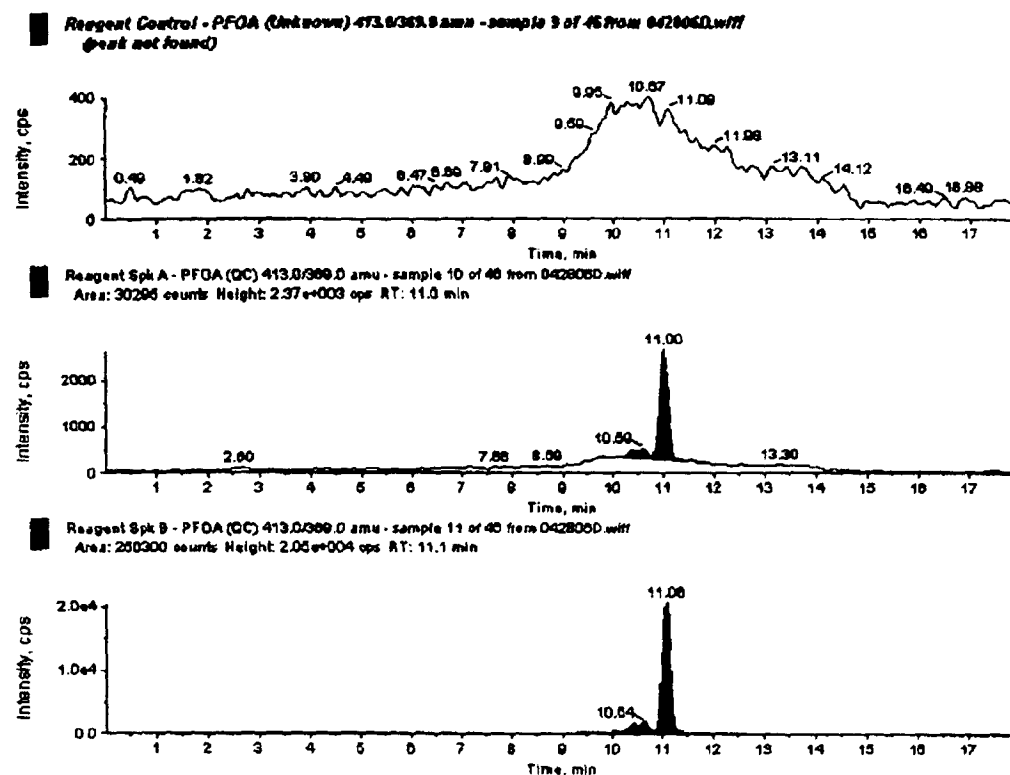


Figure 4. Chromatogram Representing a Sediment Sample Analyzed for PFOA (Exygen ID: C0172892, Data Set: 042806D)

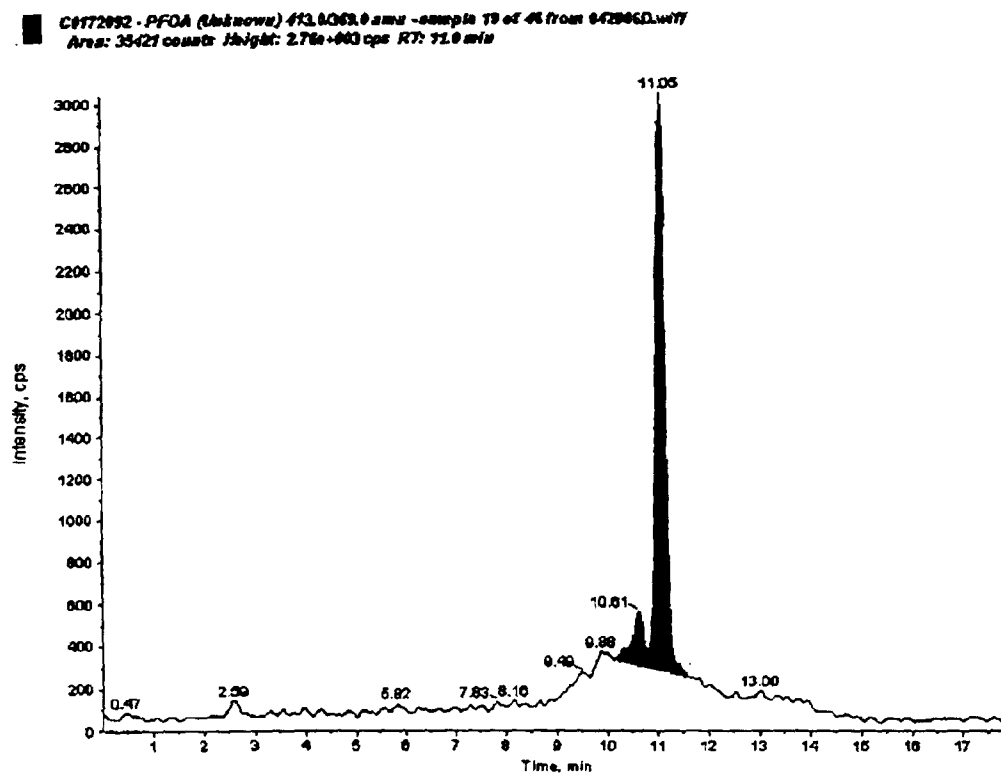


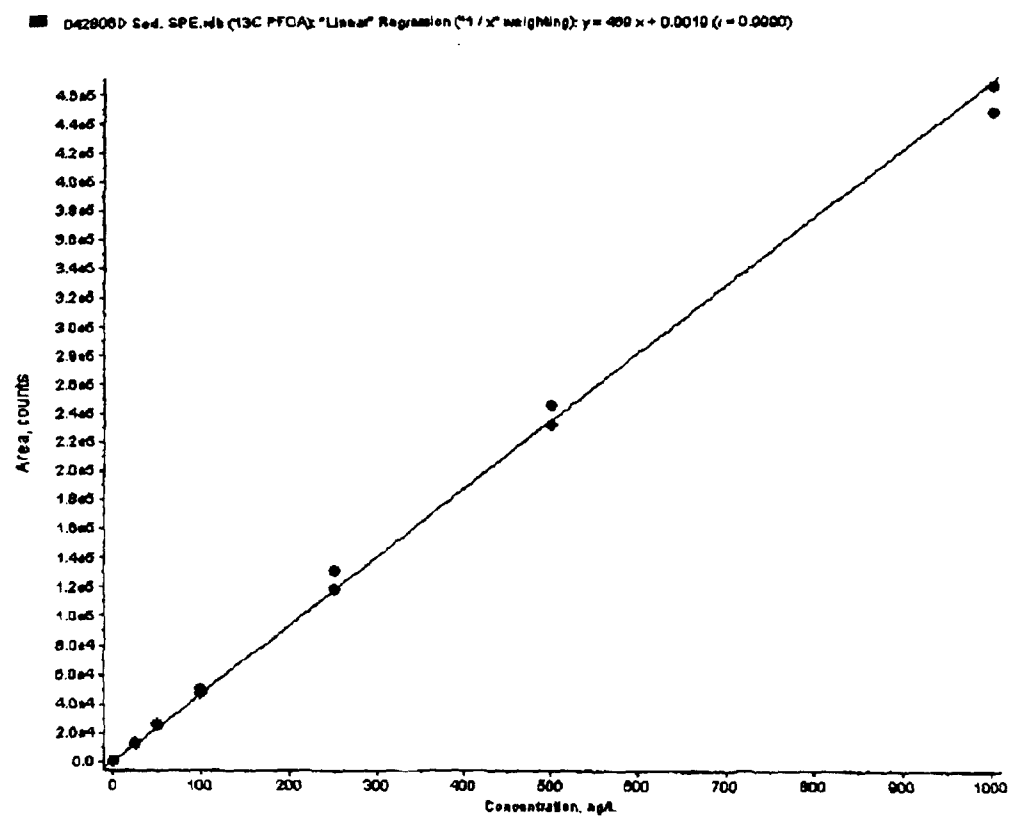
Figure 5. Typical Calibration Curve for ^{13}C PFOA in Methanol

Figure 6. Extracted Standards of ^{13}C PFOA in Methanol, 25 ng/L and 50 ng/L, Respectively

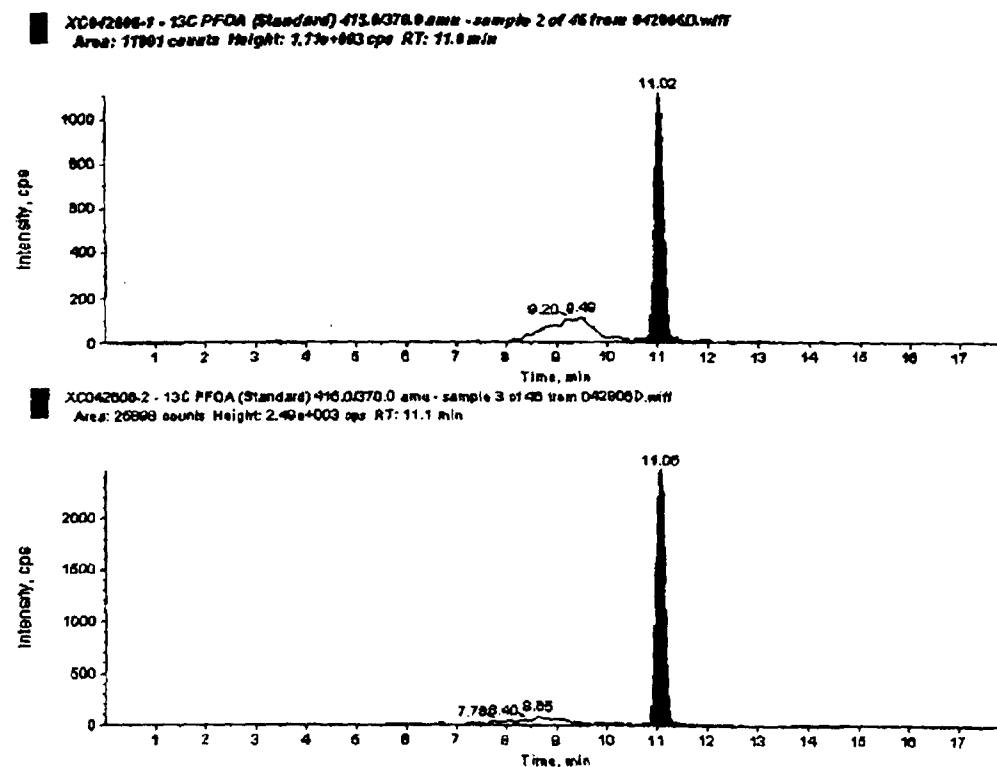
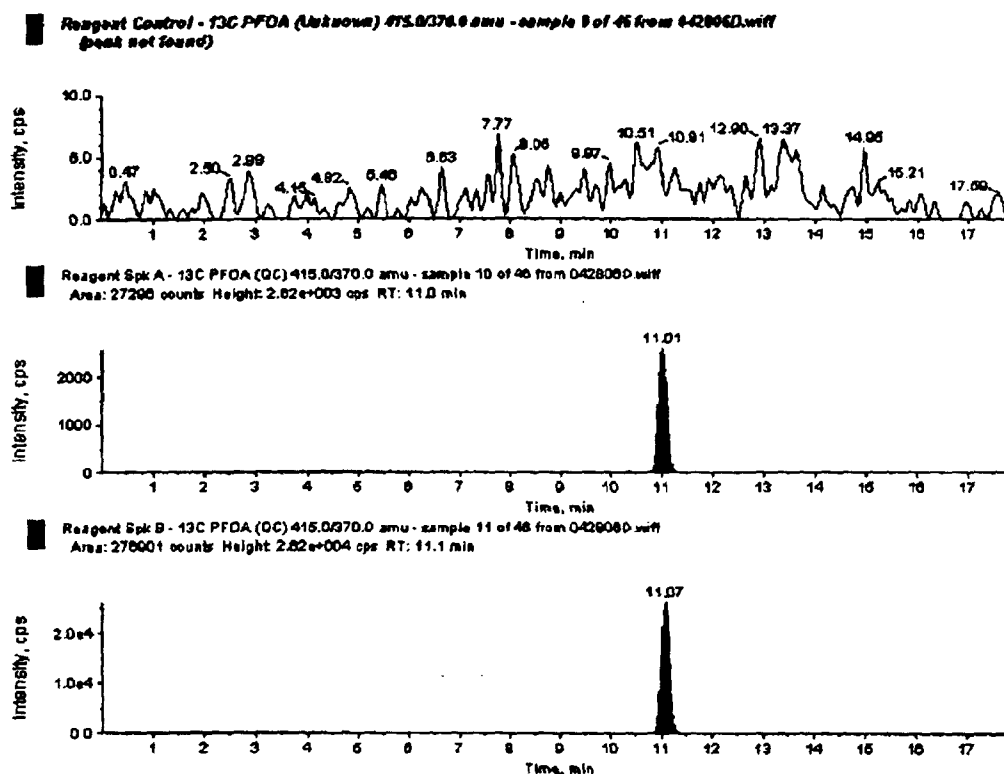


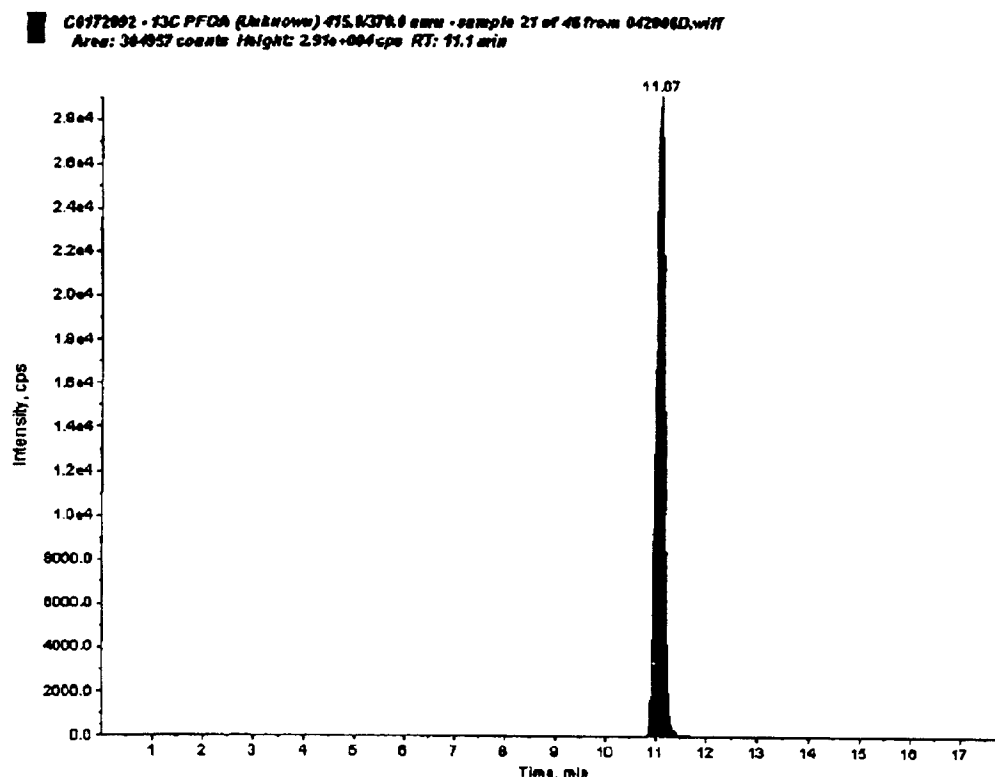
Figure 7. ^{13}C PFOA in Reagent Blank, 0.4 ng/mL Fortified Reagent Spk A, and 4.0 ng/mL Fortified Reagent Spk B, Respectively



Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Figure 8. Chromatogram Representing a Sediment Sample Analyzed for ^{13}C PFOA (Exygen ID: C0172892, Data Set: 042806D)



APPENDIX A

Study Protocol P0000760 (Exygen Study No. P0000760) with Analytical Method and Protocol Amendments

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

STUDY PROTOCOL

Study Title:

**Analysis of Perfluorooctanoic Acid (PFOA) in Water, Soil,
Sediment, Fish, Clams, Vegetation, Small Mammal Liver and
Small Mammal Serum Using LC/MS/MS for the 3M Decatur
Monitoring Program**

Exygen Protocol Number: P0000760

Performing Laboratory:

Exygen Research
3058 Research Drive
State College, PA 16801
Phone: (814) 272-1039

Sponsor Representative:

Michael A. Santoro
Director of Regulatory Affairs
3M Building 0236-01-B-10
St. Paul, MN 55144
Phone: (651) 733-6374

Page 1 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

DISTRIBUTION:

- 1) Jaisimha Kesari, Study Director, Weston Solutions
- 2) John M. Flaherty, Principal Investigator, Exygen Research
- 3) Michael A. Santoro, Sponsor Representative, 3M Company
- 4) Exygen Research Quality Assurance Unit

Page 7 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

PROTOCOL APPROVAL

Study Title: Analysis of Perfluorooctanoic Acid (PFOA) in Water, Soil, Sediment, Fish,
Clams, Vegetation, Small Mammal Livers and Small Mammal Serum
Using LC/MS/MS for the 3M Decatur Monitoring Program

Exygen Protocol Number: P0000760

APPROVALS



Jaisimha Kesari, Study Director
Weston Solutions

11/5/04

Date



Michael A. Santoro, Sponsor Representative
3M Company

11/11/04

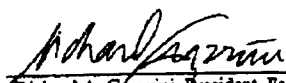
Date



John M. Flaherty, Principal Investigator
Exygen Research

10/29/04

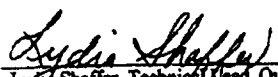
Date



Richard A. Guzzini, President, Facility Management
Exygen Research

29-0X-2004

Date



Lydia Shaffer, Technical Lead, Quality Assurance Unit
Exygen Research

10/29/04

Date

Page 3 of 63

Exygen Protocol Number: P0000760

TABLE OF CONTENTS

TITLE PAGE.....	1
DISTRIBUTION.....	2
PROTOCOL APPROVAL.....	3
TABLE OF CONTENTS.....	4
INTRODUCTION.....	5
TEST MATERIAL.....	5
SURROGATE FIELD SPIKE COMPOUND.....	5
OBJECTIVE.....	6
TESTING FACILITY.....	6
STUDY DIRECTOR.....	6
SPONSOR REPRESENTATIVE.....	7
PRINCIPAL INVESTIGATOR.....	7
PROPOSED EXPERIMENTAL START AND TERMINATION DATES.....	7
IDENTIFICATION AND JUSTIFICATION OF THE TEST SYSTEM.....	7
SAMPLE PROCUREMENT, RECEIPT AND RETENTION.....	8
SAMPLE IDENTIFICATION.....	8
ANALYTICAL PROCEDURE SUMMARY.....	9
VERIFICATION OF ANALYTICAL PROCEDURE.....	9
METHOD FOR CONTROL OF BIAS.....	11
STATISTICAL METHODS.....	11
GLP STATEMENT.....	11
REPORT.....	11
SAFETY AND HEALTH.....	12
AMENDMENTS TO PROTOCOL.....	13
DATA RECORD KEEPING.....	13
QUALITY ASSURANCE.....	14
RETENTION OF DATA AND ARCHIVING.....	14
APPENDIX I, ANALYTICAL METHODS.....	15

Page 4 of 65

Exygen Protocol Number: P0000760

INTRODUCTION

The purpose of this study is to perform analysis for perfluorooctanoic acid (PFOA) in water, soil, sediment, fish, clams, vegetation, small mammal livers and small mammal serum using LC/MS/MS for the 3M Decatur Monitoring Program.

The study will be audited for compliance with EPA TSCA Good Laboratory Practice Standards 40 CFR 792 by the Quality Assurance Unit of Exygen Research.

TEST MATERIAL

The test material is perfluorooctanoic acid (PFOA) and was purchased from Sigma-Aldrich.

PFOA

Chemical Name: Perfluorooctanoic acid

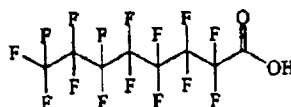
Molecular Weight: 414

Lot Number: 23116HB

Purity: 97.64%

Transitions Monitored: 413 → 369 (for quantification) and
413 → 219 (for confirmation)

Structure:



SURROGATE FIELD SPIKE COMPOUND

The surrogate field spiking compound is ^{13}C labeled perfluorooctanoic acid (^{13}C PFOA) and was purchased by 3M from Perkin-Elmer Life and Analytical Sciences.

^{13}C PFOA

Chemical Name: 1,2- ^{13}C perfluorooctanoic acid

Molecular Weight: 416

Lot Number: 3507-195

Purity: 97%

Transition Monitored: 415 → 370

Page 5 of 63

Exygen Protocol Number: P0000760

OBJECTIVE

The purpose of this study is to perform analysis for perfluorooctanoic acid (PFOA) in water, soil, sediment, fish, clams, vegetation, small mammal livers and small mammal serum for the 3M Decatur Monitoring Program using the current versions of the following Exygen analytical methods:

- V0001780: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS"
- V0001781: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by LC/MS/MS"
- V0001782: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS"
- V0001783: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Clams by LC/MS/MS"
- V0001784: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS"
- V0001785: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS"
- V0001786: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Serum by LC/MS/MS"

TESTING FACILITY

Exygen Research
3058 Research Drive
State College, PA 16801
Phone: (814) 272-1039

STUDY DIRECTOR

Jasimha Kesari P.E., DEE
Weston Solutions, Inc.
1400 Weston Way
West Chester, PA 19380
Phone: (610) 701-3761
Fax: (610) 701-7401
j.kesari@westonsolutions.com

Page 6 of 65

Exygen Protocol Number: P0000760

SPONSOR REPRESENTATIVE

Michael A. Santoro
3M Company
Director of Regulatory Affairs
3M Building 0236-01-B-10
St. Paul, MN 55144
Phone: (651) 733-6374

PRINCIPAL INVESTIGATOR

John M. Flaherty
Exygen Research
3058 Research Drive
State College, PA 16801
Phone: (814) 272-1039
john.flaherty@exygen.com

PROPOSED EXPERIMENTAL START AND TERMINATION DATES

It is proposed that the analytical portion of this study be conducted from October 01, 2004 to December 31, 2005. The actual experimental start and termination dates will be included in the final report.

IDENTIFICATION AND JUSTIFICATION OF THE TEST SYSTEM

The following are the test systems for this study:

- Water (groundwater and surface water)
- Soil
- Sediment
- Fish
- Clams
- Vegetation
- Small Mammal Liver
- Small Mammal Serum

The samples will be collected by Weston Solutions. The control samples will be purchased and prepared by the testing facility. Purchase and processing

Page 7 of 65

Exygen Protocol Number: P0000760

details for the control samples will be included in the final report associated with this study.

The test systems were chosen to assess the environmental impact of PFOA in the Decatur, Alabama area.

SAMPLE PROCUREMENT, RECEIPT AND RETENTION

Water, soil, sediment, fish, clam, vegetation, small mammal liver and small mammal serum samples will be received at Exygen directly from Weston Solutions. The details of sample procurement for this study are outlined in the 3M work plan entitled "Phase 2 Work Plan for Sampling Environmental Media." The number and types of samples collected will vary depending on availability in the field. The total number of samples received and analyzed for each matrix will be documented in the final report associated with this study.

Water, soil, and sediment samples will be used as received without further processing at Exygen. These samples will be stored refrigerated at 2°C-8°C. Fish, clam, vegetation and small mammal liver samples will be processed according to the appropriate analytical method (see Appendix I). These samples will be stored frozen at ≤ -10°C. Small mammal whole blood samples will be centrifuged in the field at the time of collection and the serum fraction will be used for the study. Small mammal serum will be stored frozen at ≤ -10°C.

The receipt and processing of the samples will be documented in the final report and raw data associated with the study.

SAMPLE IDENTIFICATION

Prior to analysis, each sample will be assigned a laboratory sample reference number. The reference number will be unique and will distinguish each laboratory sample that is processed throughout the analytical procedure. Chromatographic data will be identified by the laboratory sample reference number.

Sample storage conditions and locations will be documented throughout the study.

Page 8 of 65

Exygen Protocol Number: P0000760

ANALYTICAL PROCEDURE SUMMARY

References:

- V0001780: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS"
- V0001781: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by LC/MS/MS"
- V0001782: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS"
- V0001783: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Clams by LC/MS/MS"
- V0001784: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS"
- V0001785: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS"
- V0001786: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Serum by LC/MS/MS"

The above methods use analytical conditions capable of separating the isomers of PFOA. The final report will include the isomers summed into total PFOA found.

VERIFICATION OF ANALYTICAL PROCEDURE

A laboratory control sample will be used for the preparation of fortified control samples. The test substance will be made into solutions as per the method, and added to the matrices via a micropipette.

For water sampling, Exygen will supply one bottle per sample collected. The bottles will be 500 mL precleaned Sci/Spec Premier wide mouth HDPE bottles. These bottles have been routinely used for PFOA sample collection at the testing facility and have been shown to be free of PFOA. All containers used for water sample collection will be shipped to the sample location containing 100 ng of ¹³C PFOA (equivalent to 500 ng/L in the final sample). Samples will be added to each container to a volumetric fill line at 200 mL. A field duplicate, a low field spike and a high field spike of each sample will be collected. The low and high field spikes will contain PFOA as well as perfluorobutanesulfonate (PFBS), perfluorohexanesulfonate (PFHS), and perfluorooctanesulfonate (PFOS). These compounds are included in the solutions used to spike the samples. The results for PFBS, PFHS and PFOS will not be reported in this study. Exygen will supply one field blank (control water) and two field blank spikes (control water fortified with PFOA at a low

Page 9 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

and high level) for every twenty samples collected. At the testing facility, each water sample (excluding field duplicates and field spikes) will be extracted in duplicate and will also be fortified at a low and high concentration with PFOA and processed through the described procedure to determine method accuracy and to check for bias.

For soil, sediment, clams, and vegetation, Exygen will supply one 500 mL pre-cleaned Sci/Spec Premier wide mouth HDPE bottle per sample collected or a zip-seal bag. All containers/bags used for sample collection will be shipped to the sample location. Samples will be added to each container or bag in the field. At the testing facility, each sample will be extracted in duplicate and will also be fortified at a known concentration with PFOA at both a low and high level and processed through the described procedure to determine method accuracy and to check for bias. ^{13}C -PFOA will also be added to each sample in the laboratory prior to extraction at the level specified below.

For small mammal liver, Exygen will supply a 50 mL polypropylene centrifuge tube. For small mammal serum, Exygen will supply a collection kit for each sample containing serum separator tubes (red top), vacutainers, needle holders and needles, transfer pipettes, and polypropylene tubes. At the testing facility, each liver and serum sample will be extracted in duplicate and will also be fortified at a known concentration with PFOA at both a low and high level and processed through the described procedure to determine method accuracy and to check for bias. ^{13}C -PFOA will also be added to each sample in the laboratory prior to extraction at the level specified below.

Low and high spiking levels for each matrix are defined below:

Matrix	Low PFOA Spiking Level	High PFOA Spiking Level	^{13}C -PFOA Spiking Level
Water	500 ng/L	5000 ng/L	500 ng/L
Soil	4 ng/g	40 ng/g	40 ng/g
Sediment	4 ng/g	40 ng/g	40 ng/g
Fish	10 ng/g	100 ng/g	100 ng/g
Clams	10 ng/g	100 ng/g	100 ng/g
Vegetation	10 ng/g	100 ng/g	100 ng/g
Small Mammal Liver	10 ng/g	100 ng/g	100 ng/g
Small Mammal Serum	10 ng/g	100 ng/mL	100 ng/mL

Recoveries are anticipated to be between 70% and 130% of the fortified levels; however, the exact precision and accuracy will be determined by the analysis of the quality control samples described above. A statement of accuracy will be included in the final report.

Prior to sample analysis, a laboratory control spike study will be conducted in water to justify the use of ^{13}C PFOA results to establish precision and accuracy of PFOA samples. The following samples are to be prepared:

Page 10 of 61

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

- Two control water blanks
- Two control water samples fortified at 250 ng/L with ^{13}C PFOA
- Two control water samples fortified at 250 ng/L with ^{13}C PFOA and at 250 ng/L with PFOA.
- Two control water samples fortified at 5000 ng/L with ^{13}C PFOA and at 5000 ng/L with PFOA.

The above samples are to be extracted and analyzed according to method V0001780 entitled "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS."

METHOD FOR CONTROL OF BIAS

Control of bias will be addressed by taking representative sub-samples from a homogeneous mixture of each matrix from untreated control samples, and by analyzing at least two levels of fortifications.

STATISTICAL METHODS

Statistics will be limited to those specified in the subject methods and to the calculation of average recoveries, as applicable.

GLP STATEMENT

All aspects of this study shall be performed and reported in compliance with EPA TSCA Good Laboratory Practice Standards 40 CFR 792. The final report or data package (supplied to the Sponsor) shall contain a statement that the study was conducted in compliance with current and applicable GLP standards and will outline any deviations in the study from those standards. This statement will be signed by the Study Director and Sponsor Representative.

REPORT

A final report will be prepared by the principal investigator or their designee at the conclusion of the study. The report will include, but will not be limited to, the following:

Page 11 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

- The name and address of the Study Director, Sponsor Representative, and of the testing facility.
- A statement of GLP compliance (any related documentation, such as chain-of-custody records, must be in the study records).
- The signed and dated statement by the Exygen Research Quality Assurance Unit regarding dates of study inspections and dates findings were reported to the Study Director and Management.
- A description of the exact analytical conditions employed in the study. If the subject method was followed exactly, it is necessary to include only a copy of the analytical method. Any modifications to this method will be incorporated into the report. If the method is photo-reduced, the project number and page number must be included on each page.
- Description of the instrumentation used and operating conditions.
- All results from all sets analyzed. Control and fortified samples will be identified and the data table will include sample number and fortification level.
- Representative chromatograms for each analyte in each matrix, including chromatograms of a standard and a control sample, and a chromatogram at a fortification level. The location of the analyte peaks will be clearly identified in all chromatograms.
- All circumstances that may have affected the quality or integrity of the data will be documented in the report.
- Locations where raw data and the final report are to be archived.
- Additions or corrections to the final report shall be in the form of an amendment signed by the Study Director. The amendment shall clearly identify that part of the report that is being altered and the reasons for the alterations. The amendment will be signed and dated by the Study Director and the Sponsor Representative.

SAFETY AND HEALTH

- Laboratory personnel will practice good sanitation and health habits.
- Every reasonable precaution shall be taken to prevent inadvertent exposure of personnel and the environment to the test or reference substance(s).

Page 12 of 65

Exygen Protocol Number: P0000760

AMENDMENTS TO PROTOCOL

All significant changes to the analytical protocol outlined here will be expressed in writing, signed and dated by the Study Director and Sponsor Representative. Amendments usually will be issued prior to initiation of study plan change. However, when a change is required without sufficient time for the issue of a written amendment, that change may be effected verbally with supporting documentation signed and dated by the Study Director and followed with a written amendment as soon as possible. In this case, the effective date of the written amendment will be the date of the documented change. Copies of the signed amendments will be appended to all distributed study plan copies. The original amendment will be maintained with the original study plan. Any deviations from the study plan or from the analytical method as provided will be documented and reported promptly to the Sponsor Representative.

DATA RECORD KEEPING

Records to be maintained include the following (as appropriate):

- Sample tracking sheet(s)
- Sample receipt records, storage history, and chains of custody
- History and preparation of standards (stock, fortification, calibration)
- Description of any modifications to the method
- Instrument run sheets, bench-sheets or logs
- Analytical data tables
- All chromatographic and instrumental conditions
- Sample extraction and analysis dates
- A complete listing of study personnel, signatures and initials
- Chronological presentation of all study correspondence
- Any other documentation necessary for the reconstruction of the study

Chromatograms- All chromatograms will contain the following:

- Sample identification, injection date, arrow or other indication of the area of interest, and injection number corresponding to the run.
- Additionally, fortifications will include the amount of analyte added and the sample number of the sample that was fortified.
- Analytical standard chromatograms will additionally include the concentration (e.g., µg/mL).

Page 13 of 63

Exygen Protocol Number: P0000760

- As part of the documentation the following sheets will be included in each analytical set: a run sheet listing the samples to be run in the set, and an instrument conditions sheet describing the instrument type and operating conditions.
- All applicable requirements for reporting of study results as per 40 CFR 792.185.

QUALITY ASSURANCE

The QA Unit of Exygen Research will inspect the study at intervals adequate to assure compliance with GLP's, and will report the findings of audits to the Study Director, Exygen Management, and the Sponsor Representative.

RETENTION OF DATA AND ARCHIVING

All hard copy raw data, including, but not limited to, the original chromatograms, worksheets, correspondence, and results shall be included with the data package submitted to the Study Director. These will be archived with the original study plan, amendments, final report, and all pertinent information from the Sponsor.

The testing facility shall keep all electronic raw data and any instrument, equipment, and storage logs for the period of time specified in 40 CFR 792.195. An exact copy of the materials submitted to the study director will also be kept at Exygen Research.

Exygen will obtain permission from the study director before discarding or returning samples.

Page 14 of 63

Exygen Protocol Number: P0000760

APPENDIX I

ANALYTICAL METHODS

- V0001780: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS"
- V0001781: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by LC/MS/MS"
- V0001782: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS"
- V0001783: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Clams by LC/MS/MS"
- V0001784: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS"
- V0001785: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS"
- V0001786: "Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Serum by LC/MS/MS"

Page 13 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

ANALYTICAL METHOD

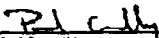
Method Number: V0001780

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water
by LC/MS/MS

Analytical Testing Facility:

Exygen Research
3055 Research Drive
State College, PA 16801

Approved By:


Paul Connolly
Technical Leader, LC-MS, Exygen Research

08/26/04
Date


John Flaherty
Vice President, Operations, Exygen Research


Date

Total Pages: 7

Page 16 of 63

Exygen Protocol Number: P0000760

Exygen Research	Method Number V0001700
ANALYTICAL METHOD	
Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS	
1.0 Scope	
This method is to be employed for the isolation and quantitation of perfluorooctanoic acid by High Performance Liquid Chromatography coupled to a tandem Mass Spectrometric Detector (LC/MS/MS) in water.	
2.0 Safety	
2.1	Always observe safe laboratory practices.
2.2	Consult the appropriate MSDS before handling any chemical for proper safety precautions.
3.0 Sample Requirement	
3.1	At least 40 mL of test sample for extraction.
3.2	No sample processing is needed for water samples.
3.3	Samples stored refrigerated should be allowed to equilibrate to room temperature.
3.4	All samples must be thoroughly mixed before being sampled for extraction.
3.5	Any samples containing particles should be centrifuged at ~3000 rpm for ~5 minutes and the supernatant used for the extraction.
3.6	Sample collection procedures will be specified in the sampling plan for this project.
4.0 Reagents and Standards	
4.1	Water – HPLC grade
4.2	Methanol – HPLC grade
4.3	Ammonium Acetate – A.C.S. Reagent Grade
4.4	Perfluorooctanoic Acid – Sigma-Aldrich
5.0 Instrument and Equipment	
5.1	A high performance liquid chromatograph capable of pumping up to 2 solvents equipped with a variable volume injector capable of injecting 5-200 µL connected to a tandem Mass Spectrometer (LC/MS/MS).
5.2	A device to collect raw data for peak integration and quantitation.
5.3	Analytical balance capable of reading to 0.00001 g.
5.4	50 mL disposable polypropylene centrifuge tubes.
5.5	15 mL disposable polypropylene centrifuge tubes.
5.6	Disposable micropipets (50-100µL, 100-200µL).
5.7	125-mL LIXPE narrow-mouth bottles.
5.8	2 mL clear HPLC vial kit.
5.9	Disposable pipettes.
5.10	Autopipettes (100-1800 µL and 10-100 µL) with disposable tips.
5.11	Waters Sep Pak Vac 6 cc (1g) C18 SPE cartridges.

Page 2 of 7

Page 17 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001700

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS

5.12 SPE vacuum manifold.

5.13 Centrifuge capable of spinning 30 mL polypropylene tubes at 3000 rpm.

6.0 Chromatographic System

6.1 Analytical Column: Poragene RP (Keyence Scientific), 2.1 mm x 50 mm, 5µ (PN: 82505-052130)

6.2 Temperature: 30°C

6.3 Mobile Phase (A): 2 mM Ammonium Acetate in Water

6.4 Mobile Phase (B): Methanol

6.5 Gradient Program:

Time (min)	% A	% B	Flow Rate (mL/min)
0.0	65	35	0.3
1.0	65	35	0.3
8.0	25	75	0.3
20.0	25	75	0.3
22.5	65	35	0.3

6.6 Injection Volume: 15 µL (can be increased to as much as 50 µL).

6.7 Quantitation: Peak Area – external standard calibration curve.

6.8 Run Time: ~ 23 minutes.

The above conditions are intended as a guide and may be changed in order to optimize the HPLC system.

7.0 MS/MS System

7.1 Mode: Electrospray Negative MRM mode, monitoring 413 → 369 m/z.

The above conditions are intended as a guide and may be changed in order to optimize the MS/MS system.

8.0 Preparation of Solutions

8.1 Mobile Phase

8.1.1 2 mM ammonium acetate in water is prepared by adding 0.154 g of ammonium acetate to 1000 mL of water.

Alternate volumes may be prepared.

Page 1 of 2

Page 18 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0081700

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS

9.0 Standard Preparation**9.1 Standard Stock/Portification Solutions**

- 9.1.1 Prepare a stock solution of ~100 µg/mL of PFOA by weighing 10 mg of analytical standard (corrected for purity) and dilute to 100 mL with methanol in a 125-mL LDPE bottle.
- 9.1.2 A 10 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 100 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.3 A 1.0 µg/mL fortification solution of PFOA is prepared by bringing 1 mL of the 10 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.4 A 0.1 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 1.0 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.5 A 0.01 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 0.1 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.6 The stock and fortification solutions are to be stored in a refrigerator at approximately 4°C and are stable for a maximum period of 6 months from the date of preparation.

9.2 Standard Calibration Solutions

- 9.2.1 LC/MS/MS calibration standards are prepared in HPLC water. The calibration standards are processed through the extraction procedure, identical to samples.
- 9.2.2 The following is a typical example; additional concentrations may be prepared as needed.

Concentration of Portification Solution (µg/L)	Portification Volume (µL)	Volume of Portified Control Sample (mL)	Final Concentration of Calibration Standard (µg/L)*	Calibration Standard ID (example)
0	0	40	0	XCalibdry-0
50	100	40	25	XCalibdry-1
10	200	40	50	XCalibdry-2
10	400	40	100	XCalibdry-3
100	100	40	250	XCalibdry-4
100	200	40	500	XCalibdry-5
100	400	40	1000	XCalibdry-6

* The extracted concentration of the calibration standard is equal to 2x its initial concentration, due to the concentration of the standard during the extraction (SPE).
XC = extracted calibration standard.

Page 4 of 7

Page 19 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0601760

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFPOA) in Water by LC/MS/MS

- 9.2.3 A new standard solution (reagent blank) must be prepared with each set of standards extracted.
- 9.2.4 Store all extracted calibration standards in 15-mL polypropylene tubes at 2°C to 6°C, up to two weeks.
- 9.2.5 Alternate volumes and concentrations of standards may be prepared as needed.
- 10.0 Batch Set Up
- 10.1 Each batch of samples extracted (typically 20 or less) must include at least one reagent control (method blank using HPLC water) and two reagent standards limited at known concentrations (lab control spike) to verify procedural recovery for the batch.
- 10.2 Requirements for field and laboratory duplicates and spikes will be specified in the quality assurance plan for this project.
- 11.0 Sample Extraction
- 11.1 Measure 40 mL of sample or a portion of sample diluted to 40 mL with water into 50 mL polypropylene centrifuge tubes (fortify as needed, replace lid and mix well).
- 11.2 Condition the C₁₈ SPE cartridge (1 g, 6 mL) by passing 10 mL methanol followed by 5 mL of HPLC water (~2 drops/sec). Do not let column run dry.
- 11.3 Load sample on conditioned C₁₈ SPE cartridge. Discard eluate.
- 11.4 Elute with ~5 mL 100% methanol. Collect 5 mL of eluate into graduated 15 mL polypropylene centrifuge tubes (final volume = 5 mL).
- 11.5 Analyze samples using electrospray LC/MS/MS.
- 12.0 Chromatography
- 12.1 Inject the same amount of each standard, sample and fortified sample into the LC/MS/MS system. A calibration standard must precede and follow all analyzed samples.
- 12.2 Standards of PFPOA corresponding to at least five or more concentration levels must be included in an analytical set.
- 12.3 An entire set of extracted calibration standards must be included at the beginning and at the end of a sample set. Extracted standards must be interspersed between every 5-10 samples. As an alternative, an entire set of extracted calibration standards may be injected at the beginning of a set followed by extracted calibration standards interspersed every 5-10 samples (in amount for a total set of extracted standards). In either case, extracted calibration standards must be the first and last injection in a sample set.
- 12.4 Use linear standard curves for quantitation. Linear standard curves are generated for the analysis by linear regression using 1/x weighting of peak area.

Page 3 of 7

Page 20 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number: U0001780

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by LC/MS/MS**

uses calibration standard concentration using Mass Lynx 3.3 (or equivalent) software system.

- 12.5 Sample response should not exceed standard responses. Any samples that exceed standard responses should be further diluted and reanalyzed.

13.0 Acceptance Criteria

- 13.1 Chromatogram must show a peak of a daughter ion at 369 m/z from a parent of 413 m/z. The 413 m/z parent corresponds to the PFOA anion, while the daughter ion (369 m/z) represents the loss of carbon dioxide.
- 13.2 Method blanks must not contain PFOA at levels greater than the LOQ. If a blank contains PFOA at levels greater than 50 ng/L, then a new blank sample must be obtained and the entire set must be re-analyzed.
- 13.3 Repetition of control spikes and matrix spikes must be between 70-130% of their known values. If a control spike falls outside the acceptable limits, the entire set of samples should be re-analyzed. Any matrix spike outside 70-130% should be evaluated by the analyst to determine if re-extraction is warranted.
- 13.4 Any calibration standard found to be a statistical outlier by using the Hugh Bow Test, may be excluded from the calculation of the calibration curve. However, the total number of extracted calibration standards that could be excluded must not exceed 20% of the total number of extracted standards injected.
- 13.5 The correlation coefficient (R^2) for calibration curves generated must be ≥ 0.992 ($R^2 \geq 0.985$). If calibration results fall outside these limits, then appropriate steps must be taken to adjust instrument operation, and the standards of the relevant set of samples should be re-analyzed.
- 13.6 Retention times between standards and samples must not drift more than $\pm 4\%$ within an analytical run. If retention time drift exceeds this limit within an analytical run then the set must be re-analyzed.

14.0 Calculations

- 14.1 Use the following equation to calculate the amount of PFOA found (in ng/L, based on peak area) using the standard curve (linear regression parameters) generated by the Mass Lynx software program:

$$\text{PFOA found (ng/L)} = \frac{(\text{Peak area} - \text{intercept})}{\text{slope}} \times \text{DF}$$

DF = factor by which the final volume was diluted, if necessary.

Page 6 of 7

Page 21 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0801780

ANALYTICAL METHODMethod of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Water by
LC/MS/MS

- 14.2 For samples fortified with known amounts of PFOA prior to extraction, use the following equation to calculate the percent recovery.

Recovery (%) =

$$\frac{(\text{total analyte found (ug/L)} - \text{analyte found in control (ug/L)})}{\text{analyte added (ug/L)}} \times 100$$

Page 7 of 7

Page 22 of 63

Interim Report #23 - Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

ANALYTICAL METHOD
Method Number: Y0001781

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by
LC/MS/MS

Analytical Testing Facility:

Exygen Research
3056 Research Drive
State College, PA 16801

Approved By:


Paul Cunniff
Technical Leader, LC-MS, Exygen Research10/26/04
Date
John Flaherty
Vice President, Operations, Exygen Research10/26/04
Date

Total Pages: 7

Page 23 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001781

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by LC/MS/MS

1.0 Scope

This method is to be employed for the isolation and quantitation of perfluorooctanoic acid by High Performance Liquid Chromatography coupled to a tandem Mass Spectrometric Detector (LC/MS/MS) in soil.

2.0 Safety

- 2.1 Always observe safe laboratory practices.
- 2.2 Consult the appropriate MSDS before handling any chemical for proper safety precautions.

3.0 Sample Requirements

- 3.1 At least 15 g of test sample for extraction.
- 3.2 No sample pretreating is needed for soil samples.
- 3.3 Samples stored refrigerated should be allowed to equilibrate to room temperature.
- 3.4 All samples must be thoroughly mixed before being sampled for extraction.
- 3.5 Sample collection procedures will be specified in the sampling plan for this project.

4.0 Reagents and Standards

- 4.1 Water – HPLC grade
- 4.2 Methanol – HPLC grade
- 4.3 Ammonium Acetate – A.C.S. Reagent Grade
- 4.4 Perfluorooctanoic Acid – Sigma-Aldrich

5.0 Instrument and Equipment

- 5.1 A high performance liquid chromatograph capable of pumping up to 2 solvents equipped with a variable volume injector capable of injecting 3-200 μ L connected to a tandem Mass Spectrometer (LC/MS/MS).
- 5.2 A device to collect raw data for peak integration and quantitation.
- 5.3 Analytical balance capable of reading to 0.00001 g.
- 5.4 50 mL disposable polypropylene centrifuge tubes.
- 5.5 15 mL disposable polypropylene centrifuge tubes.
- 5.6 Disposable micropipets (50-100 μ L, 100-200 μ L).
- 5.7 125-mL LDPE narrow-mouth bottles.
- 5.8 2 mL clear HPLC vial lot.
- 5.9 Disposable pipettes.
- 5.10 Autopipettes (100-1000 μ L and 10-100 μ L), with disposable tips.
- 5.11 Waters Sep Pak Vac 6 or (1g) μ CI8 SPE cartridges.
- 5.12 SPE vacuum manifold.
- 5.13 Ultrasonic bath.

Page 2 of 7

Page 24 of 63

Exygen Protocol Number: P0000760

Oxygen Research Method Number: V0801161

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by
LC/MS/MS

- 5.14 Wrist-action shaker.
5.15 Centrifuge capable of spinning 50 mL polypropylene tubes at 5000 rpm.

6.0 Chromatographic System

- 6.1 Analytical Column: Phenomenex RP (Knox Scientific), 2.1 mm x 50 mm, 5µ (P/N: K2505-052130)
6.2 Temperature: 30°C
6.3 Mobile Phase (A): 2 mM Ammonium Acetate in Water
6.4 Mobile Phase (B): Methanol
6.5 Gradient Program:

Time (min)	%A	%B	Flow Rate (mL/min)
0.0	65	35	0.3
1.0	65	35	0.3
8.0	25	75	0.3
20.0	25	75	0.3
22.5	65	35	0.3

- 6.6 Injection Volume: 15 µL (can be increased to as much as 50 µL).
6.7 Quantitation: Peak Area – external standard calibration curve.
6.8 Run Time: ~ 23 minutes.

The above conditions are intended as a guide and may be changed in order to optimize the HPLC system.

7.0 MS/MS System

- 7.1 Mode: Electrospray Negative MRM mode, monitoring 413 → 369 m/z for PFOA.

The above conditions are intended as a guide and may be changed in order to optimize the MS/MS system.

8.0 Preparation of Solutions**8.1 Mobile Phase**

- 8.1.1 2 mM ammonium acetate in water is prepared by adding 0.154 g of ammonium acetate to 1000 mL of water.

Alternate volumes may be prepared.

Page 3 of 7

Page 25 of 65

Oxygen Protocol Number: P0000760

Oxygen Research

Method Number V0001781

ANALYTICAL METHODMethod of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by
LCMS/MS**9.0 Standard Preparation****9.1 Standard Stock/Fortification Solution**

- 9.1.1 Prepare a stock solution of ~100 µg/mL of PFOA by weighing 10 mg of analytical standard (corrected for purity) and dilute to 100 mL with methanol in a 125-mL LDPE bottle.
- 9.1.2 A 10 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 100 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.3 A 1.0 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 10 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.4 A 0.1 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 1.0 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.5 A 0.04 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 0.1 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.6 The stock and fortification solutions are to be stored in a refrigerator at approximately 4°C and are stable for a maximum period of 6 months from the date of preparation.

9.2 Standard Calibration Solutions

- 9.2.1 LCMS/MS calibration standards are prepared in HPLC water. The calibration standards are processed through the extraction procedure, identical to samples.
- 9.2.2 The following is a typical example; additional concentrations may be prepared as needed.

Concentration of Fortification Solution (µg/L)	Fortification Volume (mL)	Volume of Fortified Control Sample (mL)	Final Concentration of Calibration Standard (µg/L)*	Calibration Standard ID (Materials)
0	0	40	0	XCstandard-0
10	100	40	25	XCstandard-1
10	200	40	50	XCstandard-2
10	400	40	100	XCstandard-3
100	100	40	250	XCstandard-4
100	200	40	500	XCstandard-5
100	400	40	1000	XCstandard-6

* The extracted concentration of the calibration standard is equal to 8x its initial concentration, due to the concentration of the standard during the extraction (SPE).
XC = extracted calibration standard.

Page 4 of 7

Page 26 of 65

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0000761

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by LC/MS/MS

- 9.2.3 A zero standard solution (reagent blank) must be prepared with each set of standards extracted.
- 9.2.4 Store all extracted calibration standards in 15-mL polypropylene tubes at 2°C to 6°C, up to two weeks.
- 9.2.5 Alternate volumes and concentrations of standards may be prepared as needed.
- 10.0 Batch Set Up
- 10.1 Each batch of samples extracted (typically 20 or less) must include at least one reagent control (method blank using 5 mL of methanol) and two reagent controls fortified at known concentrations (lab control spikes) to verify procedural recovery for the batch.
- 10.2 Requirements for field and laboratory duplicates and spikes will be specified in the quality assurance plan for this project.
- 11.0 Sample Extraction
- 11.1 Weigh 5 g of sample into 50 mL polypropylene centrifuge tubes (fortify as needed, replace lid and mix well).
- 11.2 Add 5 mL of methanol and shake on a wrist action shaker for ~15 minutes.
- 11.3 Transfer the tubes to an ultrasonic bath and sonicate for ~15 minutes.
- 11.4 Bring the volume up to 40 mL with water in the 50 mL polypropylene centrifuge tube.
- 11.5 Centrifuge for ~10 minutes at ~3000 rpm.
- 11.6 Condition the C₁₈ SPE cartridge (1 g, 6 mL) by passing 10 mL methanol followed by 5 mL of HPLC water (~2 drops/sec). Do not let column run dry.
- 11.7 Load (desorb) the sample on the conditioned C₁₈ SPE cartridge. Discard eluate.
- 11.8 Elute with ~5 mL 100% methanol. Collect 5 mL of eluate into graduated 15 mL polypropylene centrifuge tubes (final volume ~5 mL).
- 11.9 Analyze samples using electrospray LC/MS/MS.
- 12.0 Chromatography
- 12.1 Inject the same amount of each standard, sample and fortified sample into the LC/MS/MS system. A calibration standard must precede and follow all analyzed samples.
- 12.2 Standards of PFOA corresponding to at least five or more concentration levels must be included in an analytical set.
- 12.3 An entire set of extracted calibration standards must be included at the beginning and at the end of a sample set. Extracted standards must be interspersed between every 5-16 samples. As an alternative, an entire set of

Page 5 of 7

Page 27 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Records

Method Number V0001781

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by LC/MS/MS

- extracted calibration standards may be injected at the beginning of a set followed by extracted calibration standards interspersed every 5-10 samples (to account for a second set of extracted standards). In either case, extracted calibration standards must be the first and last injection in a sample set.
- 12.4 Use linear standard curves for quantitation. Linear standard curves are generated for the analyte by linear regression using 1/x weighting of peak area versus calibration standard concentration using MnasLyon 3.3 (or equivalent) software system.
- 12.5 Sample response should not exceed standard responses. Any samples that exceed standard responses should be further diluted and reanalyzed.
- 13.0 Acceptance Criteria
- 13.1 Chromatogram must show a peak of a daughter ion at 369 m/z from a parent of 413 m/z. The 413 m/z parent corresponds to the PFOA anion, while the daughter ion (369 m/z) represents the loss of carbon dioxide.
- 13.2 Method blanks must not contain PFOA at levels greater than the LOQ. If a blank contains PFOA at levels greater than 50 ng/L, then a new blank sample must be obtained and the entire set must be re-extracted.
- 13.3 Recoveries of control spikes and matrix spikes must be between 70-130% of their known values. If a control spike falls outside the acceptable limits, the entire set of samples should be re-extracted. Any matrix spikes outside 70-130% should be evaluated by the analyst to determine if re-extraction is warranted.
- 13.4 Any calibration standard found to be a statistical outlier by using the Nogue Test, may be excluded from the calculation of the calibration curve. However, the total number of extracted calibration standards that could be excluded must not exceed 20% of the total number of extracted standards injected.
- 13.5 The correlation coefficient (R^2) for calibration curves generated must be ≥ 0.992 ($R^2 \geq 0.985$). If calibration results fall outside these limits, then appropriate steps must be taken to adjust instrument operation, and the standards or the relevant set of samples should be reanalyzed.
- 13.6 Retention times between standards and samples must not drift more than $\pm 4\%$ within an analytical run. If retention time drift exceeds this limit within an analytical run then the set must be reanalyzed.

Page 6 of 7

Page 28 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0801781

ANALYTICAL METHODS

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Soil by LC/MS/MS

14.0 Calculations

- 14.1 Use the following equation to calculate the amount of PFOA found (in ng/L, based on peak area) using the standard curve (linear regression parameters) generated by the Mass Lynx software program:

$$\text{PFOA found (ng/L)} = \frac{(\text{Peak area} - \text{intercept})}{\text{slope}} \times \text{DF}$$

DF = factor by which the final volume was diluted, if necessary.

- 14.2 For samples fortified with known amounts of PFOA prior to extraction, use the following equation to calculate the percent recovery:

Recovery (%) =

$$\frac{(\text{total analyte found (ng/L)} - \text{analyte found in control (ng/L)})}{\text{analyte added (ng/L)}} \times 100$$

- 14.3 Use the following equation to convert the amount of PFOA found in ng/L to ng/g (ppb).

$$\text{PFOA found (ppb)} = \frac{(\text{PFOA found (ng/L)} \times \text{volume extracted (0.041 L)})}{\text{sample weight (5 g)}}$$

- 14.4 Use the following equation to calculate the amount of PFOA found in ppb based on dry weight.

$$\text{PFOA found (ppb) dry weight} = \text{PFOA found (ppb)} \times [100\% / \text{total solids}\%]$$

Page 7 of 7

Page 29 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

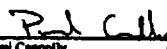
ANALYTICAL METHOD
Method Number: V0001782

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in
Sediment by LC/MS/MS

Analytical Testing Facility:

Exygen Research
3058 Research Drive
State College, PA 16801

Approved By:


Paul Connolly
Technical Leader, LC-MS, Exygen Research
Date
John Flaherty
Vice President, Operations, Exygen Research
Date

Total Pages: 7

Page 30 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research	Method Number V0001782
ANALYTICAL METHOD	
Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS	
1.0 Scope	This method is to be employed for the isolation and quantitation of perfluorooctanoic acid by High Performance Liquid Chromatography coupled to a tandem Mass Spectrometric Detector (LC/MS/MS) in sediment.
2.0 Safety	2.1 Always observe safe laboratory practices. 2.2 Consult the appropriate MSDS before handling any chemical for proper safety precautions.
3.0 Sample Requirement	3.1 At least 30 g of wet sample for extraction. 3.2 No sample processing is needed for sediment samples. 3.3 Samples stored refrigerated should be allowed to equilibrate to room temperature. 3.4 All samples must be thoroughly mixed before being sampled for extraction. 3.5 Sample collection procedures will be specified in the sampling plan for this project.
4.0 Reagents and Standards	4.1 Water – HPLC grade 4.2 Methanol – HPLC grade 4.3 Acetic Acid – Reagent grade 4.4 Ammonium Acetate – A.C.S. Reagent Grade 4.5 Perfluorooctanoic Acid – Sigma-Aldrich
5.0 Instrument and Equipment	5.1 A high performance liquid chromatograph capable of pumping up to 2 µL/min equipped with a variable volume injector capable of injecting 5-200 µL connected to a tandem mass Spectrometer (LC/MS/MS). 5.2 A device to collect raw data for peak integration and quantitation. 5.3 Analytical balance capable of reading to 0.00001 g. 5.4 50 mL disposable polypropylene centrifuge tubes. 5.5 15 mL disposable polypropylene centrifuge tubes. 5.6 Disposable micropipets (20-100µL, 100-200µL). 5.7 125-mL LDPE narrow-mouth bottles. 5.8 2 mL clear HPLC vial kit. 5.9 Disposable pipettes. 5.10 Autopipettor (100-1000 µL and 10-100 µL), with disposable tips. 5.11 Waters Sep Pak Vac 6 cc (1g) C18 SPE cartridges. 5.12 SPE vacuum manifold.

Page 2 of 7

Page 31 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0001781

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS

- 5.13 Vortexer.
 5.14 Waist-action shaker.
 5.15 Centrifuge capable of spinning 50 mL polypropylene tubes at 3000 rpm.

6.0 Chromatographic System

- 6.1 Analytical Column: Phenomenex RP (Kayson Scientific), 2.1 mm x 50 mm, 5µ (PN: 82985-032130)
 6.2 Temperature: 30°C
 6.3 Mobile Phase (A): 3 mM Ammonium Acetate in Water
 6.4 Mobile Phase (B): Methanol
 6.5 Gradient Program:

Time (min)	%A	%B	Flow Rate (mL/min)
0.0	65	35	0.3
1.0	65	35	0.3
8.0	25	75	0.3
20.0	25	75	0.3
22.5	65	35	0.3

- 6.6 Injection Volume: 15 µL (can be increased to as much as 50 µL).
 6.7 Quantitation: Peak Area – external standard calibration curve.
 6.8 Run Time: ~ 23 minutes.

The above conditions are intended as a guide and may be changed in order to optimize the HPLC system.

7.0 MS/MS System

- 7.1 Mode: Electrospray Negative MRM mode, monitoring 413 → 369 m/z for PFOA.

The above conditions are intended as a guide and may be changed in order to optimize the MS/MS system.

8.0 Preparation of Solutions**8.1 Mobile Phase**

- 8.1.1 3 mM ammonium acetate in water is prepared by adding 0.154 g of ammonium acetate to 1000 mL of water.

Page 3 of 7

Page 32 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0801782

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LOMMSMS****8.2 Extraction Solutions**

- 8.2.1 1% acetic acid in water is prepared by adding 10 mL of acetic acid to 1000 mL of water.

Alternate volumes may be prepared.

9.0 Standard Preparation**9.1 Standard Stock/Verification Solution**

- 9.1.1 Prepare a stock solution of ~100 µg/mL of PFOA by weighing 10 mg of analytical standard (uncertified for purity) and dilute to 100 mL with methanol in a 125-mL LDPE bottle.
- 9.1.2 A 10 µg/mL verification solution of PFOA is prepared by bringing 10 mL of the 100 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.3 A 1.0 µg/mL verification solution of PFOA is prepared by bringing 10 mL of the 10 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.4 A 0.1 µg/mL verification solution of PFOA is prepared by bringing 10 mL of the 1.0 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.5 A 0.01 µg/mL verification solution of PFOA is prepared by bringing 10 mL of the 0.1 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.6 The stock and verification solutions are to be stored in a refrigerator at approximately 4°C and are stable for a minimum period of 6 months from the date of preparation.

9.2 Standard Calibration Solutions

- 9.2.1 LCMS/MS calibration standards are prepared in methanol via dilution of the 0.1 µg/mL verification solution.
- 9.2.2 The following is a typical example; additional concentrations may be prepared as needed.

Concentration of Verification Solution (µg/mL)	Volume (mL)	Diluted to (mL)	Final Concentration (µg/mL)
100	10	100	10.0
100	5	100	5.0
100	2	100	2.0
10	10	100	1.0
5	10	100	0.5
2	10	100	0.2

Page 4 of 7

Page 33 of 63

Interim Report #23 - Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research	Method Number V0001762
ANALYTICAL METHOD	
Method of Analysis for the Determination of Perfluorooctanoic Acid (PFDA) in Sediment by LC/MS/MS	
9.2.3	Store all calibration standards in 125-mL LDPE narrow-mouth bottles at 2°C to 6°C, up to six months.
9.2.4	Alternate volumes and concentrations of standards may be prepared as needed.
10.0	Batch Set Up
10.1	Each batch of samples extracted (typically 20 or less) must include at least one untreated control and two untreated controls fortified at known concentrations (sub control spikes) to verify procedural recovery for the batch.
10.2	Requirements for field and laboratory duplicates and spikes will be specified in the quality assurance plan for this project.
11.0	Sample Extractions
11.1	Weigh 5 g of sample into 50 mL polypropylene centrifuge tubes (fortify as needed, replace lid and mix well).
11.2	Add 35 mL of 1% acetic acid, cap, vortex and shake on a wrist action shaker for 40 minutes.
11.3	Centrifuge the tubes at ~3000 rpm for ~30 minutes.
11.4	Condition the C ₁₈ SPE cartridges (1 g, 6 mL) by passing 10 mL methanol followed by 20 mL of HPLC water (~2 drops/sec). Do not let column run dry.
11.5	Load (descent) the sample on the conditioned C ₁₈ SPE cartridge. Discard eluate.
11.6	Add 20 mL of methanol to the sediment left in the bottom of the 50 mL centrifuge tube. Cap, vortex and shake on a wrist action shaker for ~30 minutes.
11.7	Centrifuge the tubes at ~3000 rpm for ~30 minutes.
11.8	Descent the methanol onto the same SPE cartridge. Collect the eluate.
11.9	Wash the column with 4 mL of methanol. Collect the eluate and add it to the eluate collected in step 11.3.
11.10	Condition a second C ₁₈ SPE cartridge (1 g, 6 mL) by passing 10 mL methanol followed by 20 mL of HPLC water (~2 drops/sec). Do not let column run dry.
11.11	Add the methanol to ~200 mL of water and load on the second conditioned SPE cartridge.
11.12	Elute with ~5 mL 100% methanol. Collect 5 mL of eluate into graduated 15 mL polypropylene centrifuge tubes (final volume = 5 mL).
11.13	Analyze samples using electrospray LC/MS/MS.

Page 5 of 7

Page 34 of 63

Exygen Protocol Number: P0000760

Exygen Research Method Number V0601712

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS

12.0 Chromatography

- 12.1 Inject the same amount of each standard, sample and fortified sample into the LC/MS/MS system. A calibration standard must precede and follow all sample samples.
- 12.2 Standards of PFOA corresponding to at least five or more concentration levels must be included in an analytical set.
- 12.3 An entire set of extracted calibration standards must be included at the beginning and at the end of a sample set. Standards must be interspersed between every 5-10 samples. As an alternative, an entire set of calibration standards may be injected at the beginning of a set followed by calibration standards interspersed every 5-10 samples (to account for a second set of standards). In either case, calibration standards must be the first and last injection in a sample set.
- 12.4 Use linear standard curves for quantification. Linear standard curves are generated for the analyte by linear regression using 1/x weighting of peak area versus calibration standard concentration using MassLynx 3.3 (or equivalent) software system.
- 12.5 Sample responses should not exceed standard responses. Any samples that exceed standard responses should be further diluted and reanalyzed.

13.0 Acceptance Criteria

- 13.1 Chromatograms must show a peak of a daughter ion at 369 amu from a parent of 413 amu. The 413 amu parent corresponds to the PFOA anion, while the daughter ion (369 amu) represents the loss of carbon dioxide.
- 13.2 Method blanks must not contain PFOA at levels greater than the LOQ. If a blank contains PFOA at levels greater than 0.2 ug/mL, then a new blank sample must be obtained and the entire set must be re-extracted.
- 13.3 Recoveries of control spikes and matrix spikes must be between 70-130% of their known values. If a control spike falls outside the acceptable limits, the entire set of samples should be re-extracted. Any matrix spike outside 70-130% should be evaluated by the analyst to determine if re-extraction is warranted.
- 13.4 Any calibration standard found to be a statistical outlier by using the Hugh River Test, may be excluded from the calculation of the calibration curve. However, the total number of extracted calibration standards that could be excluded must not exceed 20% of the total number of extracted standards injected.
- 13.5 The correlation coefficient (R^2) for calibration curves generated must be ≥ 0.992 ($R^2 \geq 0.985$). If calibration results fall outside these limits, then appropriate steps must be taken to adjust instrument operation, and the standards or the relevant set of samples should be reanalyzed.

Page 6 of 7

Page 35 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research Method Number V8081782

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Sediment by LC/MS/MS

13.6 Retention times between standards and samples must not drift more than $\pm 4\%$ within an analytical run. If retention time drift exceeds this limit within an analytical run then the set must be reanalyzed.

14.0 Calculations

14.1 Use the following equation to calculate the amount of PFOA found (in ng/mL, based on peak area) using the standard curve (linear regression parameters) generated by the Mini Lynx software program:

$$\text{PFOA found (ng/mL)} = \frac{(\text{Peak area} - \text{intercept})}{\text{slope}} \times \text{DF}$$

DF = factor by which the final volume was diluted, if necessary.

14.2 For samples fortified with known amounts of PFOA prior to extraction, use the following equation to calculate the percent recovery.

Recovery (%) =

$$\frac{[\text{total analyte found (ng/mL)} - \text{analyte found in control (ng/mL)}]}{\text{analyte added (ng/mL)}} \times 100$$

14.3 Use the following equation to convert the amount of PFOA found in ng/mL to ng/g (ppb).

$$\text{PFOA found (ppb)} = \frac{[\text{PFOA found (ng/mL)} \times \text{final volume (5 mL)}]}{\text{sample weight (5 g)}}$$

14.4 Use the following equation (if necessary) to calculate the amount of PFOA found in ppb based on dry weight.

$$\text{PFOA found (ppb) dry weight} = \text{PFOA found (ppb)} \times [100\% / \text{moist solids}(\%)]$$

Page 7 of 7

Page 36 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

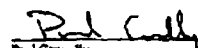
ANALYTICAL METHOD
Method Number: V0001783

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFDA) in Fish
and Clams by LC/MS/MS

Analytical Testing Facility:

Exygen Research
3048 Research Drive
State College, PA 16801

Approved By:


Paul Connolly
Technical Leader, LC-MS, Exygen Research11/26/04
Date
John Flaherty
Vice President, Operations, Exygen Research
Date

Total Pages: 8

Page 37 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0001713

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Clams by LC/MS/MS

1.0 Scope

This method is to be employed for the isolation and quantitation of perfluorooctanoic acid by High Performance Liquid Chromatography coupled to a tandem Mass Spectrometric Detector (LC/MS/MS) in fish and clams.

2.0 Safety

- 2.1 Always observe safe laboratory practices.
- 2.2 Consult the appropriate MSDS before handling any chemical for proper safety precautions.

3.0 Sample Requirements

- 3.1 At least 20 g of wet sample for extraction.
- 3.2 Samples should be preserved before extraction. Place the frozen sample in a food processor and homogenize with dry ice. Place the samples in containers and leave open in frozen storage overnight to allow for carbon dioxide sublimation. Seal and place the samples in frozen storage until time of analysis.
- 3.3 Sample collection procedures will be specified in the sampling plan for this project.

4.0 Reagents and Standards

- 4.1 Water – HPLC grade
- 4.2 Acetonitrile – HPLC grade
- 4.3 Carbon (120-400 mesh) – Reagent grade
- 4.4 Methanol – HPLC grade
- 4.5 Silica gel (60-200 mesh) – Reagent grade
- 4.6 Florisil (60-100 mesh) – Reagent grade
- 4.7 Supelco LC-NH₂ – Reagent grade
- 4.8 1-Octanol – HPLC grade
- 4.9 L-Ascorbic acid – Reagent grade
- 4.10 Dimethylsiloxane – Reagent grade
- 4.11 Toluene – Reagent grade
- 4.12 Ammonium Acetate – A.C.S. Reagent Grade
- 4.13 Perfluorooctanoic Acid – Sigma-Aldrich

5.0 Instrument and Equipment

- 5.1 A high performance liquid chromatograph capable of pumping up to 2 solvents equipped with a variable volume injector capable of injecting 5-200 µL, connected to a tandem Mass Spectrometer (LC/MS/MS).
- 5.2 A device to collect raw data for peak integration and quantitation.
- 5.3 Analytical balance capable of reading to 0.0001 g.

Page 2 of 3

Page 38 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001787

ANALYTICAL METHODMethod of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and
Clams by LC/MS/MS

- 5.4 Rotary evaporator.
- 5.5 Thermometer.
- 5.6 125 mL pear-shaped flask.
- 5.7 50 mL disposable polypropylene centrifuge tubes.
- 5.8 15 mL disposable polypropylene centrifuge tubes.
- 5.9 Disposable micropipets (50-100 μ L, 100-200 μ L).
- 5.10 125-mL LDPE narrow-mouth bottles.
- 5.11 2 mL clear HPLC vial kit.
- 5.12 Disposable pipettes.
- 5.13 Autopipettes (100-1800 μ L and 10-100 μ L), with disposable tips.
- 5.14 SPE tubes (50mL) (Supelco cat. no. NS57177).
- 5.15 Vortex action shaker.
- 5.16 Centrifuge capable of spinning 50 mL polypropylene tubes at 2000 rpm.

6.0 Chromatographic System

- 6.1 Analytical Column: Fluoropak RP (Kysentec Scientific), 2.1 mm x 50 mm, 5 μ (P/N: 82363-052120)
- 6.2 Temperature: 30°C
- 6.3 Mobile Phase (A): 2 mM Ammonium Acetate in Water
- 6.4 Mobile Phase (B): Methanol
- 6.5 Gradient Program:

Time (min)	%A	%B	Flow Rate (mL/min)
0.0	65	35	0.3
1.0	65	35	0.3
8.0	25	75	0.3
20.0	25	75	0.3
22.5	65	35	0.3

- 6.6 Injection Volume: 15 μ L (can be increased to as much as 50 μ L).
- 6.7 Quantitation: Peak Area – external standard calibration curve.
- 6.8 Run Time: ~ 23 minutes.

The above conditions are intended as a guide and may be changed in order to optimize the HPLC system.

Page 3 of 8

Page 39 of 65

Exygen Protocol Number: P0000760

Exygen Source	Method Number: V0001781
ANALYTICAL METHOD	
Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Clams by LC/MS/MS	
7.0 MS/MS System	
7.1 Mode: Electrospray Negative MRM mode, monitoring 413 → 369 m/z for PFOA.	
The above conditions are intended as a guide and may be changed in order to optimize the MS/MS system.	
8.0 Preparation of Solutions	
8.1 Mobile Phase	
8.1.1 2 mM ammonium acetate in water is prepared by adding 0.154 g of ammonium acetate to 1000 mL of water.	
8.2 Buffering Solution	
8.2.1 2% acetic acid in methanol is prepared by dissolving 2 g of acetic acid in 100 mL of methanol.	
8.2.2 30% Dimethylchlorosilane in solvent is prepared by bringing 3 mL of dimethylchlorosilane to a final volume of 10 mL with solvent.	
Alternate volumes may be prepared.	
9.0 Standard Preparation	
9.1 Standard Stock/Calibration Solution	
9.1.1 Prepare a stock solution of ~100 µg/mL of PFOA by weighing 10 mg of analytical standard (certified for purity) and dilute to 100 mL with methanol in a 125-mL LDPE bottle.	
9.1.2 A 1.0 µg/mL fortification solution of PFOA is prepared by bringing 1 mL of the 100 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.	
9.1.3 A 0.1 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 1.0 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.	
9.1.4 A 0.01 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 0.1 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.	
9.1.5 The stock and fortification solutions are to be stored in a refrigerator at approximately 4°C and are stable for a maximum period of 6 months from the date of preparation.	

Page 4 of 6

Page 10 of 63

Exygen Protocol Number: P0000760

Oxygen Research Method Number: V0001713

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Clams by LC/MS/MS

9.2 Standard Calibration Solutions

9.2.1 LC/MS/MS calibration standards are prepared in methanol via dilution of the 1.0 µg/mL certification solution.

9.2.2 The following is a typical example; additional concentrations may be prepared as needed.

Concentration of Certification Solution (µg/mL)	Volume (mL)	Diluted to (mL)	Final Concentration (µg/mL)
1.0	5.0	100	0.05
1.0	2.5	100	0.025
1.0	1.0	100	0.01
0.05	10	100	0.005
0.025	10	100	0.0025
0.1	10	100	0.001
0.005	10	100	0.0005

9.2.3 Store all calibration standards in 125-mL LDPE narrow-mouth bottles at 2°C to 4°C, up to six months.

9.2.4 Alternate volumes and concentrations of standards may be prepared as needed.

10.0 Batch Set Up

10.1 Each batch of samples extracted (typically 20 or less) must include at least one untreated control and two untreated controls fortified at known concentrations (lab control spikes) to verify procedural recovery for the batch.

10.2 Requirements for field and laboratory duplicates and spikes will be specified in the quality assurance plan for this project.

11.0 Sample Extraction

11.1 Weigh 5 g of frozen sample into 50 mL polypropylene centrifuge tubes (fluffy as needed, replace lid and mix well).

11.2 Add 30 mL of acetonitrile and shake on a wrist action shaker for ~15 minutes.

11.3 Place the tubes in a fume hood for ~1 hour.

11.4 Pack and condition the SPE tubes and attach the pear-shaped flask.

11.5 Pack the 20 mL SPE tubes in sequence with 2 g Florisil, 2 g silica gel, 2 g octadecyl, and 1 g LC-NH₂. Condition the columns with 20 mL of methanol, then 20 mL of acetonitrile. Discard all washes. Do not allow the column to dry.

11.6 Silanize the 125 mL pear-shaped flasks by rinsing with the 30% dimethylchlorosilane in toluene solution. Rinse the flask with toluene once, followed by methanol (three times). Dry the flasks completely before use, either by air-drying or with a stream of nitrogen.

Page 1 of 8

Page 41 of 65

Interim Report #23 – Analysis of Sediment Samples

Oxygen Study No.: P0000760

Oxygen Protocol Number: P0000760

Oxygen Research

Method Number V000171J

ANALYTICAL METHODMethod of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and
Clams by LC/MS/MS

- 11.7 Centrifuge the 50 mL polypropylene tubes containing sample at ~2000 rpm for ~10 minutes.
 - 11.8 Decant the extract on to a conditioned SPE column fitted inside the mouth of the pear-shaped flask. Collect the eluate in the 125 mL silanized pear-shaped flask.
 - 11.9 Add 10 mL of acetonitrile to the sample in the 50 mL centrifuge tube. Homogenize the frozen fat phase using a homogenizer for ~30 seconds and rinse the homogenizer with ~10 mL of acetonitrile into the tube.
 - 11.10 Shake the sample again for ~10 minutes on a wrist-action shaker.
 - 11.11 Place the tubes in a freezer for ~1 hour more.
 - 11.12 Centrifuge the 50 mL polypropylene tubes containing sample at ~2000 rpm for ~10 minutes.
 - 11.13 Decant the extract onto the same SPE column. Collect the eluate into the same pear-shaped flask and combine with the eluate from the initial extraction.
 - 11.14 Run 20 mL of acetonitrile through the SPE column and combine the eluate in the same pear-shaped flask.
 - 11.15 Add 3-4 drops of 1-octanol to the extract in the pear-shaped flask and evaporate at reduced pressure using a rotary evaporator (at < 40°C).
 - 11.16 Reduce the final volume, by adding 2 mL of 2% ascorbic acid in methanol to the pear-shaped flask and swirl to mix/dissolve.
 - 11.17 Transfer the extract to HPLC vials using disposable pipets.
 - 11.18 Analyze samples using electrospray LC/MS/MS.
- 12.0 Chromatography
- 12.1 Inject the same amount of each standard, sample and fortified sample into the LC/MS/MS system. A calibration standard must precede and follow all analyzed samples.
 - 12.2 Standards of PFOA corresponding to at least five or more concentration levels must be included in an analytical set.
 - 12.3 An entire set of calibration standards must be included at the beginning and at the end of a sample set. Standards must be interspersed between every 5-10 samples. As an alternative, an entire set of calibration standards may be injected at the beginning of a set followed by calibration standards interspersed every 5-10 samples (to account for a second set of standards). In either case, calibration standards must be the first and last injections in a sample set.
 - 12.4 Use linear standard curves for quantitation. Linear standard curves are generated for the analyte by linear regression using 1/x weighting of peak area versus calibration standard concentration using MAM-Lyon 3.3 (or equivalent) software system.

Page 5 of 8

Page 42 of 65

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001713

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Crusts by LC/MS/MS**

- 12.5 Sample response should not exceed standard responses. Any samples that exceed standard responses should be further diluted and reanalyzed.

13.0 Acceptance Criteria

- 13.1 Chromatograms must show a peak of a daughter ion at 369 amu from a parent of 413 amu. The 413 amu parent corresponds to the PFOA anion, while the daughter ion (369 amu) represents the loss of carbon dioxide.
- 13.2 Method blanks must not contain PFOA at levels greater than the LOQ. If a blank contains PFOA at levels greater than 0.5 ppb, then a new blank sample must be obtained and the entire set must be re-extracted.
- 13.3 Recoveries of control spikes and matrix spikes must be between 70-130% of their known values. If a control spike falls outside the acceptable limits, the entire set of samples should be re-extracted.
- 13.4 Any calibration standard found to be a statistical outlier by using the Hugs Rule Test, may be excluded from the calculation of the calibration curve. However, the total number of calibration standards that could be excluded must not exceed 20% of the total number of standards injected.
- 13.5 The correlation coefficient (R) for calibration curves generated must be ≥ 0.992 ($R^2 \geq 0.985$). If calibration results fall outside these limits, then appropriate steps must be taken to adjust instrument operation, and the standards or the relevant set of samples should be reanalyzed.
- 13.6 Retention times between standards and samples must not drift more than $\pm 0.4\%$ within an analytical run. If retention time drift exceeds this limit within an analytical run then the run must be reanalyzed.

14.0 Calculations

- 14.1 Use the following equation to calculate the amount of PFOA found (in ng/mL, based on peak area) using the standard curve (linear regression parameters) generated by the Mass Lynx software program:

$$\text{PFOA found (ng/mL)} = \frac{(\text{Peak Area} - \text{Intercept})}{\text{slope}}$$

- 14.2 Use the following equation to convert the amount of PFOA found in ng/mL to ng/g (ppb).

$$\text{PFOA found (ppb)} = \frac{(\text{PFOA found (ng/mL)} \times \text{final volume (mL)}) \times \text{DF}}{\text{sample weight (g)}}$$

DF = factor by which the final volume was diluted, if necessary.

Page 7 of 6

Page 43 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001783

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Fish and Clams by LC/MS/MS**

- 14.3 For samples fortified with known amounts of PFOA prior to extraction, use the following equation to calculate the percent recovery.

Recovery (%) =

$$\frac{(\text{total analyte found (ng/g)} - \text{analyte found in control (ng/g)})}{\text{analyte added (ng/g)}} \times 100$$

Page 5 of 8

Page 44 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

ANALYTICAL METHOD

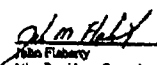
Method Number: V0001784

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in
Vegetation by LC/MS/MS

Analytical Testing Facility:

Exygen Research
3038 Research Drive
State College, PA 16801

Approved By:


Paul Connolly
Technical Leader, LC-MS, Exygen Research10/26/04
Date
John Flaherty
Vice President, Operations, Exygen Research

Date

Total Pages: 7

Page 43 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V089176

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS

1.0 Scope

This method is to be employed for the isolation and quantitation of perfluorooctanoic acid by High Performance Liquid Chromatography coupled to a tandem Mass Spectrometric Detector (LC/MS/MS) in vegetation.

2.0 Safety

- 2.1 Always observe safe laboratory practices.
- 2.2 Consult the appropriate MSDS before handling any chemical for proper safety precautions.

3.0 Sample Requirement

- 3.1 At least 20 g of test sample for extractions.
- 3.2 Samples should be processed before extraction. Place the frozen sample in a food processor and homogenize with dry ice. Place the samples in containers and leave open in frozen storage overnight to allow for carbon dioxide sublimation. Seal and place the samples in frozen storage until time of analysis.
- 3.3 Sample collection procedures will be specified in the sampling plan for this project.

4.0 Reagents and Standards

- 4.1 Water – HPLC grade
- 4.2 Acetonitrile – HPLC grade
- 4.3 Carbon (120-400 mesh) – Reagent grade
- 4.4 Methanol – HPLC grade
- 4.5 Silica gel (60-200 mesh) – Reagent grade
- 4.6 Florisil (60-100 mesh) – Reagent grade
- 4.7 Supelchem LC-NH₂ – Reagent grade
- 4.8 1-Octanol – HPLC grade
- 4.9 L-Ascorbic acid – Reagent grade
- 4.10 Dimethylsiloxane – Reagent grade
- 4.11 Toluene – Reagent grade
- 4.12 Ammonium Acetate – A.C.S. Reagent Grade
- 4.13 Perfluorooctanoic Acid – Sigma-Aldrich

5.0 Instrument and Equipment

- 5.1 A high performance liquid chromatograph capable of pumping up to 2 mL/min equipped with a variable volume injector capable of injecting 5-200 µL connected to a tandem Mass Spectrometer (LC/MS/MS).
- 5.2 A device to collect raw data for peak integration and quantitation.
- 5.3 Analytical balance capable of reading to 0.00001 g.

Page 2 of 7

Page 46 of 63

Exygen Protocol Number: P0000760

Exygen Research Method Number: V0001764

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS

- 5.4 Rotary evaporator.
- 5.5 125 mL pear-shaped flasks.
- 5.6 50 mL disposable polypropylene centrifuge tubes.
- 5.7 15 mL disposable polypropylene centrifuge tubes.
- 5.8 Disposable micropipettes (50-100µL, 100-200µL).
- 5.9 125-mL LEPN narrow-mouth bottles.
- 5.10 2 mL clear HPLC vial kit.
- 5.11 Disposable pipettes.
- 5.12 Autopipettes (100-1000 µL and 10-100 µL), with disposable tips.
- 5.13 SPE tubes (20mL) (Supelco cat. no. NS37177).
- 5.14 Wrist action shaker.
- 5.15 Centrifuge capable of spinning 50 mL polypropylene tubes at 2000 rpm.

6.0 Chromatographic System

- 6.1 Analytical Column: Fluorose RP (Keyence Scientific), 2.1 mm x 50 mm, 5µm (P/N: K2303-002134)
- 6.2 Temperature: 30°C
- 6.3 Mobile Phase (A): 2 mM Ammonium Acetate in Water
- 6.4 Mobile Phase (B): Methanol
- 6.5 Gradient Program:

Time (min)	%A	%B	Flow Rate (mL/min)
0.0	65	35	0.3
1.0	65	35	0.3
1.0	25	75	0.3
20.0	25	75	0.3
22.5	65	35	0.3

- 6.6 Injection Volume: 15 µL (can be increased to as much as 50 µL).
- 6.7 Quantitation: Peak Area – external standard calibration curve.
- 6.8 Run Time: ~ 23 minutes.

The above conditions are intended as a guide and may be changed in order to optimize the HPLC system.

7.0 MS/MS System

- 7.1 Mode: Electrospray Negative MRM mode, monitoring 413 → 369 m/z for PFOA.

Page 3 of 7

Page 47 of 55

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001754

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS

The above conditions are intended as a guide and may be changed in order to optimize the MS/MS system.

8.0 Preparation of Solutions**8.1 Mobile Phase**

8.1.1 2 mM ammonium acetate in water is prepared by adding 0.154 g of ammonium acetate to 1000 mL of water.

8.2 Extraction Solutions

8.2.1 2% acetic acid in methanol is prepared by dissolving 2 g of acetic acid in 100 mL of methanol.

8.2.2 50% Dimethylchlorosulfone in solvent is prepared by bringing 3 mL of dimethylchlorosulfone to a final volume of 10 mL with solvent.

All stock volumes may be prepared.

9.0 Standard Preparation**9.1 Standard Stock/Purification Solutions**

9.1.1 Prepare a stock solution of ~100 µg/mL of PFOA by weighing 10 mg of analytical standard (corrected for purity) and dilute to 100 mL with methanol in a 125-mL LDPE bottle.

9.1.2 A 1.0 µg/mL fortification solution of PFOA is prepared by bringing 1 mL of the 100 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.

9.1.3 A 0.1 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 1.0 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.

9.1.4 A 0.01 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 0.1 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.

9.1.5 The stock and fortification solutions are to be stored in a refrigerator at approximately 4°C and are stable for a maximum period of 6 months from the date of preparation.

9.2 Standard Calibration Solutions

9.2.1 LC/MS/MS calibration standards are prepared in methanol via dilution of the 1.0 µg/mL fortification solution.

Page 4 of 7

Page 48 of 63

Exygen Protocol Number: P0000760

Exygen Research	Method Number: V0001760
ANALYTICAL METHOD	
Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS	

9.2.2 The following is a typical example; additional concentrations may be prepared as needed.

Concentration of Fortification Solution (µg/mL)	Volume (mL)	Diluted to (mL)	Final Concentration (µg/mL)
1.0	5.0	100	0.05
1.0	2.5	100	0.025
1.0	1.0	100	0.01
0.05	10	100	0.005
0.025	10	100	0.0025
0.1	10	100	0.001
0.001	10	100	0.0005

9.2.3 Store all calibration standards in 125-mL LDPE narrow-mouth bottles at 2°C to 6°C, up to six months.

9.2.4 Alternate volumes and concentrations of standards may be prepared as needed.

10.0 Batch Set Up

- 10.1 Each batch of samples extracted (typically 20 or less) must include at least one instrument control and two measured controls (fortified at known concentrations (lab control spikes) to verify procedural recovery for the batch.
- 10.2 Requirements for field and laboratory duplicates and spikes will be specified in the quality assurance plan for this project.

11.0 Sample Extraction

- 11.1 Weigh 5 g of frozen sample into 50 mL polypropylene centrifuge tubes (fortify as needed, replace lid and mix well).
- 11.2 Add 20 mL of acetonitrile and shake on a wrist action shaker for ~15 minutes.
- 11.3 Centrifuge the 50 mL polypropylene tubes containing sample at ~2000 rpm for ~10 minutes.
- 11.4 Pack and condition the SPE tubes and alternate the pear-shaped flasks.
- 11.5 Pack the 20 mL SPE tubes in sequence with 2 g Sorbail, 2 g silica gel, 2 g carbon, and 3 g LC-NH₂. Condition the columns with 20 mL of methanol, then 20 mL of acetonitrile. Discard all washes. Do not allow the column to dry.
- 11.6 Size the 125 mL pear-shaped flasks by rinsing with the 30% dimethylchloroethane in toluene solution. Rinse the flask with toluene once, followed by methanol (three times). Dry the flasks completely before use, either by air-drying or with a stream of nitrogen.
- 11.7 Decant the extract on to a conditioned SPE column fitted inside the mouth of the pear-shaped flask. Collect the eluate in the 125 mL size the pear-shaped flask.

Page 3 of 7

Page 49 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001784

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS

- 11.8 Add 30 mL of acetonitrile to the sample in the 50 mL centrifuge tube.
 - 11.9 Shake the sample again for ~10 minutes on a wrist-action shaker.
 - 11.10 Centrifuge the 50 mL polypropylene tubes containing sample at ~2000 rpm for ~5 minutes.
 - 11.11 Decant the extract onto the same HPLC column. Collect the eluate into the same pear-shaped flask and combine with the eluent from the initial extraction.
 - 11.12 Repeat steps 11.8 through 11.11 again.
 - 11.13 Add 3-4 drops of 1-octanol to the extract in the pear-shaped flask and evaporate at reduced pressure using a rotary evaporator (at <40°C).
 - 11.14 Make the final volume, by adding 2 mL of 2% ascorbic acid in methanol to the pear-shaped flask and swirl to mix/dissolve.
 - 11.15 Transfer the extract to HPLC vials using disposable pipets.
 - 11.16 Analyze samples using electrospray LC/MS/MS.
- 12.0 Chromatography
- 12.1 Inject the same amount of each standard, sample and fortified sample into the LC/MS/MS system. A calibration standard must precede and follow all analyzed samples.
 - 12.2 Standards of PFOA corresponding to at least five or more concentration levels must be included in an analytical set.
 - 12.3 An entire set of extracted calibration standards must be included at the beginning and at the end of a sample set. Extracted standards must be interspersed between every 5-10 samples. As an alternative, an entire set of extracted calibration standards may be injected at the beginning of a set followed by extracted calibration standards interspersed every 5-10 samples (to account for a second set of extracted standards). In either case, extracted calibration standards must be the first and last injection in a sample set.
 - 12.4 Use linear standard curves for quantification. Linear standard curves are generated for the analyte by linear regression using 1/x weighting of peak area versus calibration standard concentration using MassLynx 3.3 (or equivalent) software system.
 - 12.5 Sample response should not exceed standard responses. Any samples that exceed standard responses should be further diluted and reanalyzed.
- 13.0 Acceptance Criteria
- 13.1 Chromatograms must show a peak of a daughter ion at 360 mass from a parent of 413 mass. The 413 mass parent corresponds to the PFOA anion, while the daughter ion (360 mass) represents the loss of carbon dioxide.

Page 6 of 7

Page 50 of 65

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0001784

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Vegetation by LC/MS/MS

- 13.2 Method blanks must not contain PFOA at levels greater than the LOQ. If a blank contains PFOA at levels greater than 0.5 ppb, then a new blank sample must be obtained and the entire set must be re-analyzed.
- 13.3 Recoveries of control spikes and matrix spikes must be between 70-130% of their known values. If a control spike falls outside the acceptable limits, the entire set of samples should be re-analyzed.
- 13.4 Any calibration standard found to be a statistical outlier by using the Hughes Bow Test, may be excluded from the calculation of the calibration curve. However, the total number of calibration standards that could be excluded must not exceed 20% of the total number of standards injected.
- 13.5 The correlation coefficient (R) for calibration curves generated must be ≥ 0.992 ($R^2 \geq 0.985$). If calibration results fall outside these limits, then appropriate steps must be taken to adjust instrument operation, and the standards of the relevant set of samples should be re-analyzed.
- 13.6 Retention times between standards and samples must not drift more than ± 4 % within an analytical run. If retention time drift exceeds this limit within an analytical run then the set must be re-analyzed.

14.0 Calculations

- 14.1 Use the following equation to calculate the amount of PFOA found (in ng/mL, based on peak area) using the standard curve (linear regression parameters) generated by the Mass Lynx software program:

$$\text{PFOA found (ng/mL)} = \frac{(\text{Peak Area} - \text{Intercept})}{\text{slope}}$$

- 14.2 Use the following equation to convert the amount of PFOA found in ng/mL to ng/g (ppb).

$$\text{PFOA found (ppb)} = \frac{(\text{PFOA found (ng/mL)} \times \text{final volume (mL)} \times \text{DF})}{\text{sample weight (g)}}$$

DF = factor by which the final volume was diluted, if necessary.

- 14.3 For samples fortified with known amounts of PFOA prior to extraction, use the following equation to calculate the percent recovery.

Recovery (%) =

$$\frac{(\text{total analyte found (ng/g)} - \text{analyte found in control (ng/g)})}{\text{analyte added (ng/g)}} \times 100$$

Page 7 of 7

Page 51 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

ANALYTICAL METHOD

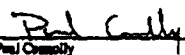
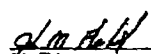
Method Number: V8001785

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS

Analytical Testing Facility:

Exygen Research
3068 Research Drive
State College, PA 16801

Approved By:


Paul Connolly
Technical Leader, LC-MS, Exygen Research10/16/04
Date
John Flaherty
Vice President, Operations, Exygen Research10/26/04
Date

Total Pages: 7

Page 52 of 65

Exygen Protocol Number: P0000760

Exygen Research	Method Number V0001781
ANALYTICAL METHOD	
Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS	
1.0	Scope
	This method is to be employed for the isolation and quantitation of perfluorooctanoic acid by High Performance Liquid Chromatography coupled to a tandem Mass Spectrometric Detector (LC/MS/MS) in small mammal liver.
2.0	Safety
2.1	Always observe safe laboratory practices.
2.2	Consult the appropriate MSDS before handling any chemical for proper safety precautions.
3.0	Sample Requirement
3.1	At least 5 g of wet sample for extraction.
3.2	Samples should be processed before extraction. Place the frozen sample in a food processor and homogenize with dry ice. Place the samples in containers and leave open in frozen storage overnight to allow for carbon dioxide sublimation. Seal and place the samples in frozen storage until time of analysis. Alternatively, if there is an insufficient amount of sample (~less than 5 g), then on processing is necessary and the sample can be used as supplied.
3.3	Sample collection procedures will be specified in the sampling plan for this project.
4.0	Reagents and Standards
4.1	Water – HPLC grade
4.2	Methanol – HPLC grade
4.3	Acetonitrile – HPLC grade
4.4	Ammonium Acetate – A.C.S. Reagent Grade
4.5	Perfluorooctanoic Acid – Sigma-Aldrich
5.0	Instrument and Equipment
5.1	A high performance liquid chromatograph capable of pumping up to 2 solvents equipped with a variable volume injector capable of injecting 1-200 µL connected to a tandem Mass Spectrometer (LC/MS/MS).
5.2	A device to collect raw data for peak integration and quantitation.
5.3	Analytical balance capable of reading to 0.00001 g.
5.4	50 mL disposable polypropylene centrifuge tubes.
5.5	15 mL disposable polypropylene centrifuge tubes.
5.6	Disposable microspigets (50-100µL, 100-200µL).
5.7	125-mL LDPE narrow-mouth bottles.
5.8	2 mL clear HPLC vial kit.

Page 1 of 3

Page 53 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0801755

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS**

- 5.9 Disposable pipettes.
- 5.10 Autopipettors (100-1000 µL and 10-100 µL), with disposable tips.
- 5.11 Waters Sep Pak Vac 6 or (1g) C18 SPE cartridges.
- 5.12 SPE vacuum manifold.
- 5.13 Thermo mixer.
- 5.14 Vortex mixer.
- 5.15 Centrifuge capable of spinning 15 mL polypropylene tubes at 3000 rpm.

6.0 Chromatographic System

- 6.1 Analytical Column: Phosphor RP (Krytox Scientific), 2.1 mm x 50 mm, 5µ (P/N: 82308-052130)
- 6.2 Temperature: 30°C
- 6.3 Mobile Phase (A): 2 mM Ammonium Acetate in Water
- 6.4 Mobile Phase (B): Methanol
- 6.5 Gradient Program:

Time (min)	% A	% B	Flow Rate (mL/min)
0.5	65	35	0.3
1.8	65	35	0.3
8.0	25	75	0.3
20.0	25	75	0.3
22.5	65	35	0.3

- 6.6 Injection Volume: 15 µL (can be increased to as much as 50 µL).
- 6.7 Quantitation: Peak Area – external standard calibration curve.
- 6.8 Run Time: ~25 minutes.

The above conditions are intended as a guide and may be changed in order to optimize the HPLC system.

7.0 MS/MS System

- 7.1 Mode: Electrospray Negative MRM mode, monitoring 413 → 369 m/z for PFOA.

The above conditions are intended as a guide and may be changed in order to optimize the MS/MS system.

Page 3 of 7

Page 54 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number V9081793

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS****8.0 Preparation of Solutions****8.1 Mobile Phase**

8.1.1 2 mM ammonium acetate in water is prepared by adding 0.154 g of ammonium acetate to 1000 mL of water.

Alternate volumes may be prepared.

9.0 Standard Preparation**9.1 Standard Stock/Fortification Solution**

9.1.1 Prepare a stock solution of ~100 µg/mL of PFOA by weighing 10 mg of analytical standard (corrected for purity) and dilute to 100 mL with methanol in a 125-mL LDPE bottle.

9.1.2 A 1.0 µg/mL fortification solution of PFOA is prepared by bringing 1 mL of the 100 µg/mL solution to a final volume of 100 with methanol in a 125-mL LDPE bottle.

9.1.3 A 0.1 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 1.0 µg/mL solution to a final volume of 100 with methanol in a 125-mL LDPE bottle.

9.1.4 The stock and fortification solutions are to be stored in a refrigerator at approximately 4°C and are stable for a maximum period of 6 months from the date of preparation.

9.2 Standard Calibration Solutions

9.2.1 LC/MS/MS calibration standards are prepared in methanol via dilution of the 0.1 µg/mL fortification solution.

9.2.2 The following is a typical example; additional concentrations may be prepared as needed.

Concentration of Fortification Solution (µg/mL)	Volume (mL)	Diluted to (mL)	Final Concentration (µg/mL)
100	1.0	100	5.0
100	2.0	100	2.0
100	1.0	100	1.0
5.0	10	100	0.5
2.0	10	100	0.2
1.0	10	100	0.1

9.2.3 Store all calibration standards in 125-mL LDPE narrow-mouth bottles at 2°C to 6°C, up to six months.

9.2.4 Alternate volumes and concentrations of standards may be prepared as needed.

Page 4 of 7

Page 55 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001763

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS

10.0 Batch Set Up

- 10.1 Each batch of samples extracted (typically 20 or less) must include at least one untreated control and two untreated controls fortified at known concentrations (lab control spikes) to verify procedural recovery for the batch.
- 10.2 Requirements for field and laboratory duplicates and spikes will be specified in the quality assurance plan for this project.

11.0 Sample Extraction

- 11.1 Weigh 1 g of sample into a 50 mL polypropylene centrifuge tubes (fortify as needed, replace lid and mix well). Note that ultimate weights of liver may be measured depending on the sample size available for use.
- 11.2 Add water to the sample for a final volume of 10 mL.
- 11.3 Homogenize sample using a tissue grinder for ~1 minute.
- 11.4 Transfer 1 mL of the sample using a disposable pipette into a 15 mL disposable centrifuge tube.
- 11.5 Add 5 mL of acetonitrile and shake for ~20 minutes on a wrist-action shaker.
- 11.6 Centrifuge the tubes at ~3000 rpm for ~5 minutes.
- 11.7 Decant the supernatant into a 50 mL disposable centrifuge tube and add 35 mL of water.
- 11.8 Condition the C₁₈ SPE cartridge (1 g, 6 mL) by passing 10 mL methanol followed by 5 mL of HPLC water (~2 drops/sec). Do not let column run dry.
- 11.9 Load the sample on conditioned C₁₈ SPE cartridge. Discard eluate.
- 11.10 Elute with ~2 mL of methanol. Collect 2 mL of eluate into a graduated 15 mL polypropylene centrifuge tube (final volume = 2 mL).
- 11.11 Analyze samples using electrospray LC/MS/MS.

12.0 Chromatography

- 12.1 Inject the same amount of each standard, sample and fortified sample into the LC/MS/MS system. A calibration standard must precede and follow all analyzed samples.
- 12.2 Standards of PFOA corresponding to at least five or more concentration levels must be included in an analytical set.
- 12.3 An entire set of calibration standards must be included at the beginning and at the end of a sample set. Standards must be interspersed between every 5-10 samples. As an alternative, an entire set of calibration standards may be injected at the beginning of a set followed by calibration standards interspersed every 5-10 samples (to account for a second set of standards). In either case, calibration standards must be the first and last injection in a sample set.
- 12.4 Use linear standard curves for quantitation. Linear standard curves are generated for the analyte by linear regression using 1/x weighting of peak area.

Page 3 of 7

Page 56 of 63

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number V800/761

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Liver by LC/MS/MS

versus calibration standard concentration using MassLynx 3.3 (or equivalent) software system.

- 12.5 Sample responses should not exceed standard responses. Any samples that exceed standard responses should be further diluted and reanalyzed.

13.0 Acceptance Criteria

- 13.1 Chromatogram must show a peak of a daughter ion at 369 m/z from a parent of 413 m/z. The 413 m/z parent corresponds to the PFOA anion, while the daughter ion (369 m/z) represents the loss of carbon dioxide.
- 13.2 Method blanks must not contain PFOA at levels greater than the LOQ. If a blank contains PFOA at levels greater than 10 ng/g, then a new blank sample must be obtained and the entire set must be re-analyzed.
- 13.3 Recovery of control spikes and matrix spikes must be between 70-130% of their known values. If a control spike falls outside the acceptable limits, the entire set of samples should be re-extracted. Any matrix spike outside 70-130% should be evaluated by the analyst to determine if re-extraction is warranted.
- 13.4 Any calibration standard found to be a statistical outlier by using the Hoge-Bever Test, may be excluded from the calculation of the calibration curve. However, the total number of calibration standards that could be excluded must not exceed 20% of the total number of standards injected.
- 13.5 The correlation coefficient (R) for calibration curves generated must be ≥ 0.992 or ≥ 0.985 . If calibration results fall outside these limits, then appropriate steps must be taken to adjust instrument operation, and the standards or the relevant set of samples should be reanalyzed.
- 13.6 Retention times between standards and samples must not drift more than ± 4 % within an analytical run. If retention time drift exceeds this limit within an analytical run then the run must be reanalyzed.

14.0 Calculations

- 14.1 Use the following equation to calculate the amount of PFOA found (in $\mu\text{g/mL}$ based on peak area) using the standard curve (linear regression parameters) generated by the Mass Lynx software program:

$$\text{PFOA found } (\mu\text{g/mL}) = \frac{(\text{Peak area} - \text{intercept})}{\text{slope}} \times \text{DF} \times \text{aliquot factor}$$

DF = factor by which the final volume was diluted, if necessary.
Aliquot factor = 10

Page 6 of 7

Page 57 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0001761

ANALYTICAL METHODMethod of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Seal
- Mammary Liver by LC/MS/MS

- 14.2 For samples fortified with known amounts of PFOA prior to extraction, use the following equation to calculate the percent recovery.

Recovery (%) =

$$\frac{[\text{total analyte found (ng/mL)} - \text{analyte found in control (ng/mL)}]}{\text{analyte added (ng/mL)}} \times 100$$

- 14.3 Use the following equation to convert the amount of PFOA found in ng/mL to ug/g (ppb).

$$\text{PFOA found (ppb)} = \frac{\text{PFOA found (ng/mL)} \times \text{final volume (mL)}}{\text{sample weight (g)}}$$

Page 7 of 7

Page 58 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

ANALYTICAL METHOD
Method Number: V0001766

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFDA) in Small
Mammal Serum by LC/MS/MS

Analytical Testing Facility:

Exygen Research
3058 Research Drive
State College, PA 16801

Approved By:


Paul Connolly
Technical Leader, LC-MS, Exygen Research10/26/01
Date
Alan Fleherty
Vice President, Operations, Exygen Research11/2/01
Date

Total Pages: 7

Page 59 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number VED01788

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Serum by LC/MS/MS

1.0 Scope

This method is to be employed for the isolation and quantitation of perfluorooctanoic acid by High Performance Liquid Chromatography coupled to a tandem Mass Spectrometric Detector (LC/MS/MS) in small mammal serum.

2.0 Safety

- 2.1 Always observe safe laboratory practices.
- 2.2 Consult the appropriate MSDS before handling any chemical for proper safety precautions.

3.0 Sample Requirements

- 3.1 At least 1 mL of test sample for extraction.
- 3.2 No sample processing is needed for serum samples. However, frozen serum samples must be allowed to completely thaw to room temperature before use.
- 3.3 Sample collection procedures will be specified in the sampling plan for this project.

4.0 Reagents and Standards

- 4.1 Water – HPLC grade
- 4.2 Methanol – HPLC grade
- 4.3 Acetonitrile – HPLC grade
- 4.4 Ammonium Acetate – A.C.S. Reagent Grade
- 4.5 Perfluorooctanoic Acid – Sigma-Aldrich

5.0 Instrument and Equipment

- 5.1 A high performance liquid chromatograph capable of purging up to 2 solvents equipped with a variable volume injector capable of injecting 5-200 μ L connected to a tandem Mass Spectrometer (LC/MS/MS).
- 5.2 A device to collect raw data for peak integration and quantitation.
- 5.3 Analytical balance capable of reading to 0.00001 g.
- 5.4 50 mL disposable polypropylene centrifuge tubes.
- 5.5 15 mL disposable polypropylene centrifuge tubes.
- 5.6 Disposable microvials (50-100 μ L, 100-200 μ L).
- 5.7 125-mL LDPE narrow-mouth bottles.
- 5.8 2 mL clear HPLC vial kit.
- 5.9 Disposable pipettes.
- 5.10 Autopipettes (100-1000 μ L and 10-100 μ L), with disposable tips.
- 5.11 Waters Sep Pak Vac 6 or (1g) C18 SPE cartridges.
- 5.12 SPE vacuum manifold.
- 5.13 Vortexer.

Page 2 of 7

Page 60 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0001786

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Serum by LC/MS/MS

- 5.14 Vial-action shaker.
5.15 Centrifuge capable of spinning 15 mL polypropylene tubes at 3000 rpm.

6.0 Chromatographic System

- 6.1 Analytical Column: Fluorophase RP (Keyence Scientific), 2.1 mm x 50 mm, 5µ (PN: K2585-052138)
6.2 Temperature: 30°C
6.3 Mobile Phase (A): 2 mM Ammonium Acetate in Water
6.4 Mobile Phase (B): Methanol
6.5 Gradient Program.

Time (min)	% A	% B	Flow Rate (mL/min)
0.0	65	35	0.3
1.0	65	35	0.3
8.0	25	75	0.3
20.0	25	75	0.3
22.5	65	35	0.3

- 6.6 Injection Volume: 15 µL (can be increased to as much as 50 µL).
6.7 Quantitation: Peak Area – external standard calibration curve.
6.8 Run Time: ~ 25 minutes.

The above conditions are intended as a guide and may be changed in order to optimize the HPLC system.

7.0 MS/MS System

- 7.1 Mode: Electrospray Negative MRM mode, monitoring 413 → 369 m/z for PFOA.

The above conditions are intended as a guide and may be changed in order to optimize the MS/MS system.

8.0 Preparation of Solutions**8.1 Mobile Phase**

- 8.1.1 2 mM ammonium acetate in water is prepared by adding 0.154 g of ammonium acetate to 1000 mL of water.

Alternate volumes may be prepared.

Page 3 of 7

Page 61 of 65

Exygen Protocol Number: P0000760

Exygen Research

Method Number V0001704

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Invertebrate Sores by LC/MS/MS****9.0 Standard Preparation****9.1 Standard Stock/Fortification Solutions**

- 9.1.1 Prepare a stock solution of ~100 µg/mL of PFOA by weighing 10 mg of analytical standard (corrected for purity) and dilute to 100 mL with methanol in a 125-mL LDPE bottle.
- 9.1.2 A 1.0 µg/mL fortification solution of PFOA is prepared by bringing 1 mL of the 100 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.3 A 0.1 µg/mL fortification solution of PFOA is prepared by bringing 10 mL of the 1.0 µg/mL solution to a final volume of 100 with methanol in a 125 mL LDPE bottle.
- 9.1.4 The stock and fortification solutions are to be stored in a refrigerator at approximately 4°C and are stable for a maximum period of 6 months from the date of preparation.

9.2 Standard Calibration Solutions

- 9.2.1 LC/MS/MS calibration standards are prepared in methanol via dilution of the 0.1 µg/mL fortification solution.
- 9.2.2 The following is a typical example; additional concentrations may be prepared as needed.

Concentration of Fortification Solution (µg/mL)	Volume (µL)	Diluted to (mL)	Final Concentration (µg/mL)
100	5.0	100	5.0
100	2.5	100	2.5
100	1.0	100	1.0
5.0	10	100	0.5
2.5	10	100	0.2
1.0	10	100	0.1

- 9.2.3 Store all calibration standards in 125-mL LDPE narrow-mouth bottles at 2°C to 6°C, up to six months.
- 9.2.4 Alternate volumes and concentrations of standards may be prepared as needed.

10.0 Batch Set Up

- 10.1 Each batch of samples extracted (typically 20 or less) must include at least one unextracted control and two extracted controls fortified at known concentrations (lab control spikes) to verify procedural recovery for the batch.
- 10.2 Requirements for field and laboratory duplicates and spikes will be specified in the quality assurance plan for this project.

Page 4 of 7

Page 62 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number: M0001786

ANALYTICAL METHOD**Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Solid Matrix Samples by LC/MS/MS****11.0 Sample Extraction**

- 11.1 Measure 1 mL of sample into a 50 mL polypropylene centrifuge tube (fortify as needed, replace lid and mix well). Note that alternate volumes of serum may be measured depending on the sample size available for use.
- 11.2 Add water to the sample for a final volume of 20 mL. Cap tightly.
- 11.3 Vortex for ~1 minute.
- 11.4 Transfer 1 mL of the sample using a disposable pipette into a 15 mL disposable centrifuge tube.
- 11.5 Add 5 mL of acetonitrile and shake for ~20 minutes on a wrist-action shaker.
- 11.6 Centrifuge the tubes at ~3000 rpm for ~5 minutes.
- 11.7 Decant the supernatant into a 50 mL disposable centrifuge tube and add 15 mL of water.
- 11.8 Condition the C₁₈ SPE cartridge (1 g, 6 mL) by passing 10 mL methanol followed by 5 mL of HPLC water (~2 drops/sec). Do not let column run dry.
- 11.9 Load the sample on conditioned C₁₈ SPE cartridge. Discard eluent.
- 11.10 Elute with ~2 mL of methanol. Collect 2 mL of eluate into a graduated 15 mL polypropylene centrifuge tube (final volume = 2 mL).
- 11.11 Analyze samples using electrospray LC/MS/MS.

12.0 Chromatography

- 12.1 Inject the same amount of each standard, sample and fortified sample into the LC/MS/MS system. A calibration standard must precede and follow all analyzed samples.
- 12.2 Standards of PFOA corresponding to at least five or more concentration levels must be included in an analytical set.
- 12.3 An entire set of calibration standards must be included at the beginning and at the end of a sample set. Standards must be interspersed between every 5-10 samples. As an alternative, an entire set of calibration standards may be injected at the beginning of a set followed by calibration standards interspersed every 5-10 samples (to account for a second set of standards). In either case, calibration standards must be the first and last injection in a sample set.
- 12.4 Use linear standard curves for quantitation. Linear standard curves are generated for the analyte by linear regression using 1/x weighting of peak area versus calibration standard concentration using MassLynx 3.3 (or equivalent) software system.
- 12.5 Sample response should not exceed standard response. Any samples that exceed standard response should be further diluted and reanalyzed.

Page 3 of 4

Page 63 of 63

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V880179a

ANALYTICAL METHOD

Method of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small Mammal Serum by LC/MS/MS

13.0 Acceptance Criteria

- 13.1 Chromatogram must show a peak of a daughter ion at 369 *m/z* from a parent of 413 *m/z*. The 413 *m/z* parent corresponds to the PFOA anion, while the daughter ion (369 *m/z*) represents the loss of trifluoromethyl.
- 13.2 Method blanks must not contain PFOA at levels greater than the LOQ. If a blank contains PFOA at levels greater than 10 ng/mL, then a new blank sample must be obtained and the entire set must be re-extracted.
- 13.3 Recovery of control spikes and matrix spikes must be between 70-130% of their known values. If a control spike falls outside the acceptable limits, the entire set of samples should be re-extracted. Any matrix spikes outside 70-130% should be evaluated by the analyst to determine if re-extraction is warranted.
- 13.4 Any calibration standard found to be a statistical outlier by using the Hoge Error Test, may be excluded from the calculation of the calibration curve. However, the total number of calibration standards that could be excluded must not exceed 20% of the total number of standards injected.
- 13.5 The correlation coefficient (R) for calibration curves generated must be ≥ 0.992 ($R^2 \geq 0.985$). If calibration results fall outside these limits, then appropriate steps must be taken to adjust instrument operation, and the standards or the relevant set of samples should be re-analyzed.
- 13.6 Retention times between standards and samples must not drift more than $\pm 4\%$ within an analytical run. If retention time drift exceeds this limit within an analytical run then the set must be re-analyzed.

14.0 Calculations

- 14.1 Use the following equation to calculate the amount of PFOA found (in ng/mL, based on peak area) using the standard curve (linear regression parameters) generated by the Mass Lynx software program:

$$\text{PFOA found (ng/mL)} = \frac{(\text{Peak area} - \text{intercept})}{\text{slope}} \times \text{DF} \times \text{aliquot factor}$$

DF = factor by which the final volume was diluted, if necessary.
 Aliquot factor = 20

- 14.2 For samples fortified with known amounts of PFOA prior to extraction, use the following equation to calculate the percent recovery.

Recovery (%) =

$$\frac{[\text{total analyte found (ng/mL)} - \text{analyte found in control (ng/mL)}]}{\text{analyte added (ng/mL)}} \times 100$$

Page 6 of 7

Page 64 of 65

Exygen Protocol Number: P0000760

Exygen Research

Method Number: V0001786

ANALYTICAL METHODMethod of Analysis for the Determination of Perfluorooctanoic Acid (PFOA) in Small
Marine Mammal Serum by LC/MS/MS

- 14.3 Use the following equation to convert the amount of PFOA found in ng/mL to
ppb.

$$\text{PFOA found (ppb)} = \frac{\text{PFOA found (ng/mL)} \times \text{final volume (}\mu\text{L)}}{\text{sample volume (}\mu\text{L)}}$$

Page 7 of 7

Page 65 of 65

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

3058 Research Drive
State College, PA 16801Phone: 814-272-1039
Fax: 814-231-1580**PROTOCOL AMENDMENT**Amendment Number: 1Effective Date: 01/19/05Exygen Study Number: P0001131 Client Study Number: None

Page 1 of 1

DESCRIPTION OF AMENDED SECTION

- 1) Analytical Procedure Summary V0001780: Section 9.1
- 2) Verification of Analytical Procedure

AMENDED TO

- 1) Add to Section 9.1: Section 9.1.8, Alternate weights of standards may be used to prepare alternate concentrations of stock solutions as necessary. Alternate levels of fortification solutions may also be prepared.
- 2) Low and high spiking levels of the analytes for each matrix may be altered depending on sample size available for extraction and/or to cover analyte concentrations expected in the samples.

RATIONALE

- 1) Higher concentrations of standards need to be prepared in order to spike the sample bottles at higher levels.
- 2) The sample size available for small mammal liver and serum was smaller than expected. Spiking at the pre-determined levels in the protocol puts the spiked concentration lower than the detection limit. Also, the analyte levels in the ground water samples are expected to greatly exceed the pre-determined spiking levels listed in the protocol. When the levels in the samples greatly exceed the spiking levels, an accurate recovery value cannot be calculated for the QC sample. Higher spiking levels in the bottles will cover the analyte concentrations expected in the water samples.

IMPACT ON STUDY

The LOQ is 100 ng/g for a 0.1 g sample of small mammal liver and is 1000 ng/mL for a 0.01 mL sample of small mammal serum. Higher levels of spiking for the water samples will ensure that more QC recovery data can be used.

[Signature]
Study Director Signature

[Signature]
Principal Investigator Signature

[Signature]
Safety Director/Management Signature

[Signature]
Sponsor Signature (if required)

[Signature]
Date 3/9/05

[Signature]
Date 2/16/05

[Signature]
Date N/A

[Signature]
Date 2/18/05

Exygen QAU Review LJS 02/16/05

LIBRARY ID: VV0001226-6

ADMINISTRATIVE FORM

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

3058 Research Drive
State College, PA 16801Phone: 814-272-1039
Fax: 814-231-1580

PROTOCOL AMENDMENT

Amendment Number: 2 Page 1 of 1
 Effective Date: 03/07/05
 Exygen Study Number P0001131 Client Study Number: None

DESCRIPTION OF AMENDED SECTION

- 1) Report, page 11 of 85
- 2) Test Materials, page 6 of 65: PFOS transition monitored 499 → 99.

AMENDED TO

- 1) Instead of one final report, interim reports will be issued.
- 2) PFOS transition monitored may also be 499 → 80.

RATIONALE

- 1) Due to the excessive sizes of the data sets, interim reports will be issued to allow the client to receive data in a timelier manner *meanwhile* *rather than*
- 2) The API 4000 LC/MS/MS systems detect the 499 → 80 PFOS transition with greater sensitivity than the 499 → 99 transition.

IMPACT ON STUDY

- 1) The client will be able to receive and review the data more quickly.
- 2) The 499 → 80 transition can be detected with greater sensitivity, therefore, giving better chromatography.

Study Director Signature

Date 3/9/05

Principal Investigator Signature

Date 3/9/05

Study Director Management Signature

Date 8-MAR-05

Exygen Management Signature

Date 3/10/05

Sponsor Signature (if required)

Date

Exygen CAU Review: 1/23 03/10/05

LIBRARY ID: V0001228-0

ADMINISTRATIVE FORM

Interim Report #23 – Analysis of Sediment Samples

Exygen Study No.: P0000760

3058 Research Drive
State College, PA 16801Phone: 814-272-1039
Fax: 814-231-1580

DEVIATION FORM

Page 1 of 2

General:	
☒ Project Specific Deviation	☐ Facility Deviation
Date of Occurrence: 04/28/08	
Exygen Project #: P760/P1131	Deviation #: 5
Client Project #: NA	Reference #: 06-076
Regulatory Driver:	Deviation Type: (Include V# for methods and SOPs)
☐ GMP ☒ GLP ☐ Other ☐ None	☒ Protocol ☐ Method ☐ SOP V#: NA Notebook reference: NA
Sample Description:	
Login #: L0008191	Container #: NA Lot #: NA
Summary of Deviation:	
The three sediment samples in L8191 (C0172892 – C0172894) were originally extracted using the sediment method V0001782. Poor recoveries were obtained for PFOS, PFOA and ¹³C PFOA. Because of this, the study sponsor requested the use of an alternative extraction for these compounds, as follows: Direct Injection Method: Before the samples were weighed for the extraction, they were mixed thoroughly by vigorously shaking the container. A one-gram portion of sediment was weighed into a 15-milliliter centrifuge tube for the extraction. Ten milliliters of 1% acetic acid in methanol was added to each sample. The samples were then shaken by hand, vortexed, and sonicated for thirty minutes. The samples were then centrifuged for ~10 minutes at ~3000 rpm. Each sample was analyzed by LC/MS/MS electrospray. Using this method acceptable data was obtained for PFOS, but the recoveries for PFOA and ¹³C PFOA were still poor. Another alternative method was then used for PFOA and ¹³C PFOA, as follows: Alternative SPE Method: The samples that were prepared in 1% acetic acid for the direct injection method were used for this extraction. Five milliliters of each sample was aliquoted into a 50-mL polypropylene centrifuge tube and the volume was taken to 40 mL with water. The samples were then centrifuged for ~10 minutes at ~3000 rpm. The supernatant was then loaded onto a C₁₈ SPE cartridge conditioned with 10 mL of methanol and 5 mL of water. The eluate was discarded. Approximately five milliliters of methanol was added to the cartridge. Five milliliters of eluate was collected into a graduated 15 mL polypropylene centrifuge tube. Each sample was analyzed by LC/MS/MS electrospray.	
Cause: ☐ Preparation ☐ Analysis ☐ Instrument ☐ Client Request ☒ Other	
Impact:	
More usable data was obtained.	

LIBRARY ID: V0001840-6

ADMINISTRATIVE FORM

Interim Report #23 - Analysis of Sediment Samples

Exygen Study No.: P0000760

CHEM EHSR 236 1B 651 733 1958 05/01 '06 13:46 NO.106 03/03
 05-01-2008 09:08pm From: BOSTON SOLUTIONS + 1-800 P.003/003 P-819

Exygen

RESEARCH

2058 Research Drive
 State College, PA 16801

Phone: 814-272-1039
 Fax: 814-272-1980

DEVIATION FORM

Section 1 ☒ Project Specific Deviation ☐ Facility Deviation		Date of Occurrence: <u>05/01/06</u>	
Exygen Project #: <u>P76071191</u>		Deviation #: <u>1</u>	
Client Project #: <u>NA</u>		Reference #: <u> </u>	
Section 2 Description of Deviation: Deviation Issued:			
Section 3 Principal Investigator: <u>[Signature]</u> Date: <u>5/1/06</u>		Exygen Manager: <u>[Signature]</u> Date: <u>5/1/06</u>	
Study Director: <u>[Signature]</u> Date: <u>5/1/06</u>		Report Representative: <u>[Signature]</u> Date: <u>5/1/06</u>	
Quality Assurance: <u>[Signature]</u> Date: <u>5/1/06</u>		QA: <u>[Signature]</u> Date: <u>5/1/06</u>	

LIBRARY ID: V30014404

ADMINISTRATIVE FORM