

AR226-0432

CORNING Hazleton

MUTAGENICITY TEST WITH

T-6564

IN THE *SALMONELLA* - *ESCHERICHIA COLI*/MAMMALIAN-MICROSOME REVERSE
MUTATION ASSAY WITH A CONFIRMATORY ASSAY


FINAL REPORT

AUTHOR

Timothy E. Lawlor, M.A.

PERFORMING LABORATORY

Corning Hazleton Inc. (CHV)
9200 Leesburg Pike
Vienna, Virginia 22182

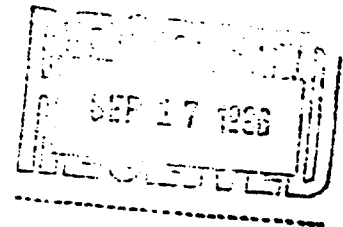

LABORATORY PROJECT ID

CHV Study No.: 17750-0-409R

SUBMITTED TO

3M Corporation
Building 220-2E-02 3M Center
St. Paul, MN 55144-1000

STUDY COMPLETION DATE

September 13, 1996

CHV Study No.: 17750-0-409R

1 of 34

000232

QUALITY ASSURANCE STATEMENT


STUDY TITLE: *Salmonella - Escherichia coli*/Mammalian-Microsome Reverse Mutation Assay with a Confirmatory Assay

ASSAY NO.: 17750-0-409R

PROTOCOL NO.: 409R, Edition 4, Modified for 3M

Quality Assurance inspections of the study and review of the final report of the above referenced project were conducted according to the Standard Operating Procedures of the Quality Assurance Unit and according to the general requirements of the appropriate Good Laboratory Practice regulations. Findings from the inspections and final report review were reported to management and to the study director on the following dates:


<u>Inspection/Date</u>	<u>Findings Reported</u>	<u>Auditor</u>
Scoring Plates - 06/28/96	06/28/96	P. Conley
Draft Report Review - 08/29/96	08/30/96	C. Orantes
Final Report Review - 09/13/96	09/13/96	C. Orantes

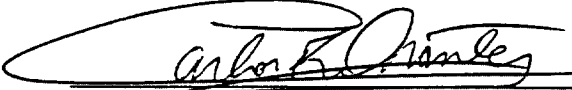
Quality Assurance Unit 9/13/96
Date Released

STUDY COMPLIANCE AND CERTIFICATION

The study was conducted in compliance with the Good Laboratory Practice regulations as set forth by the Food and Drug Administration (FDA) in Title 21 of the U.S. Code of Federal Regulations Part 58, issued December 22, 1978, (effective June 20, 1979) with any applicable amendments. There were no deviations from the aforementioned regulations or the signed protocol that would affect the integrity of the study or the interpretation of the test results. The raw data have been reviewed by the Study Director, who certifies that the evaluation of the test article as presented herein represents an appropriate conclusion within the context of the study design and evaluation criteria.

Study Director:

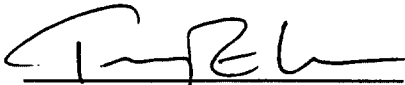

Timothy E. Lawlor, M.A.
Bacterial Mutagenesis
Genetic and Cellular Toxicology

9.13.96

Study Completion Date

TABLE OF CONTENTS

	Page No.
I. SUMMARY	5
II. STUDY INFORMATION	7
III. MATERIALS AND METHODS	9
IV. RESULTS AND CONCLUSIONS	24
V. DATA TABLES	27

CORNING Hazleton

SECTION I. SUMMARY

INTRODUCTION AND CONCLUSIONS

SUMMARY

A. Introduction

At the request of 3M Corporation, Corning Hazleton Inc. investigated T-6564 for mutagenic activity in the *Salmonella* - *Escherichia coli*/Mammalian-Microsome Reverse Mutation Assay with a Confirmatory Assay. This assay evaluated the test article and/or its metabolites for their ability to induce reverse mutations at the histidine locus in the genome of specific *Salmonella typhimurium* tester strains and at the tryptophan locus in an *Escherichia coli* tester strain both in the presence and absence of an exogenous metabolic activation system of mammalian microsomal enzymes derived from Aroclor™-induced rat liver (S9).

The doses tested in the mutagenicity assay were selected based on the results of a dose rangefinding study using tester strains TA100 and WP2uvrA and ten doses of test article ranging from 5,000 to 6.67 µg per plate, one plate per dose, both in the presence and absence of S9 mix.

The tester strains used in the mutagenicity assay were *Salmonella typhimurium* tester strains TA98, TA100, TA1535, TA1537, and *Escherichia coli* tester strain WP2uvrA. The assay was conducted with five doses of test article in both the presence and absence of S9 mix along with concurrent vehicle and positive controls using three plates per dose. The doses tested were 5,000, 3,330, 1,000, 333, and 100 µg per plate in both the presence and absence of S9 mix. The results of the initial mutagenicity assay were confirmed in an independent experiment.

B. Conclusions

The results of the *Salmonella* - *Escherichia coli*/Mammalian-Microsome Reverse Mutation Assay with a Confirmatory Assay indicate that, under the conditions of this study, 3M Corporation's test article, T-6564, did not cause a positive increase in the number of revertants per plate of any of the tester strains either in the presence or absence of microsomal enzymes prepared from Aroclor™-induced rat liver (S9).

CORNING Hazleton

SECTION II. STUDY INFORMATION

CHV Study No.: 17750-0-409R

7

000238

STUDY INFORMATION

- A. Sponsor: **3M Corporation**
- B. Test Article: **T-6564**
L-13167 FC1015-X
1. Physical Description: **clear colorless liquid**
 2. Date Received: **05/30/96**
- C. Type of Assay: *Salmonella-Escherichia coli*/Mammalian-Microsome Reverse Mutation Assay with a Confirmatory Assay
1. Protocol Number: CHV Protocol 409R, Edition 4
 2. CHV Study Number: **17750-0-409R**
- D. Study Dates
1. Study Initiation Date: **06/04/96**
 2. Experimental Start Date: **06/12/96**
 3. Experimental Termination Date: **08/27/96**
- E. Study Supervisory Personnel
- | | |
|------------------------|-------------------------|
| Study Director: | Timothy E. Lawlor, M.A. |
| Laboratory Supervisor: | Michael S. Mecchi, B.S. |

CORNING Hazleton

SECTION III. MATERIALS AND METHODS

CHV Study No.: 17750-0-409R

9

000240

MATERIALS AND METHODS

The experimental materials, methods and procedures are based on those described by Ames *et al* (1975) and Green and Muriel (1976).

MATERIALS

A. Tester Strains1. *Salmonella typhimurium*

The tester strains used were the *Salmonella typhimurium* histidine auxotrophs TA98, TA100, TA1535, and TA1537 as described by Ames *et al* (1975). The specific genotypes of these strains are shown in Table I.

TABLE I. TESTER STRAIN GENOTYPES

Histidine Mutation			Additional Mutations		
<i>his</i> G46	<i>his</i> C3076	<i>his</i> D3052	LPS	Repair	R Factor
TA1535	TA1537		<i>rfa</i>	<i>uvr</i> B	-
TA100		TA98	<i>rfa</i>	<i>uvr</i> B	+R

In addition to a mutation in the histidine operon, the tester strains contain two additional mutations which enhance their sensitivity to some mutagenic compounds. The *rfa* wall mutation results in the loss of one of the enzymes responsible for the synthesis of part of the lipopolysaccharide barrier that forms the surface of the bacterial cell wall. The resulting cell wall deficiency increases permeability to certain classes of chemicals such as those containing large ring systems (i.e. benzo(a)pyrene) that would otherwise be excluded by a normal intact cell wall.

The second mutation, a deletion of the *uvr*B gene, results in a deficient DNA excision repair system which greatly enhances the sensitivity of these strains to some mutagens. Since the *uvr*B deletion extends through the *bio* gene, all of the tester strains containing this deletion require the vitamin biotin for growth.

Strains TA98 and TA100 also contain the R-factor plasmid, pKM101, which further increases the sensitivity of these strains to some mutagens. The mechanism by which this plasmid increases sensitivity to mutagens has been suggested to be by modifying an existing bacterial DNA repair polymerase complex involved with the mismatch-repair process.

Tester strains TA98 and TA1537 are reverted from histidine dependence (auxotrophy) to histidine independence (prototrophy) by frameshift mutagens. TA1535 is reverted by base substitution mutagens and TA100 is reverted by mutagens which cause both frameshifts and base substitution mutations.

2. *Escherichia coli*

The tester strain used was the tryptophan auxotroph WP2*uvrA* as described by Green and Muriel (1976).

In addition to a mutation in the tryptophan operon, the tester strain contains a *uvrA* DNA repair deficiency which enhances its sensitivity to some mutagenic compounds. This deficiency allows the strain to show enhanced mutability since the *uvrA* repair system would normally act to remove the damaged part of the DNA molecule and accurately repair it afterwards.

Tester strain WP2*uvrA* is reverted from tryptophan dependence (auxotrophy) to tryptophan independence (prototrophy) by base substitution mutagens.

3. Source of Tester Strains

a. *Salmonella typhimurium*

The tester strains in use at CHV were received directly from Dr. Bruce Ames, Department of Biochemistry, University of California, Berkeley.

b. *Escherichia coli*

The tester strain, WP2*uvrA*, in use at CHV was received from the National Collection of Industrial Bacteria, Torrey Research Station, Scotland (United Kingdom).

4. Storage of the Tester Strains

a. Frozen Permanent Stocks

Frozen permanent stocks were prepared by growing fresh overnight cultures, adding DMSO (0.09 ml/ml of culture) and freezing small aliquots (0.5-1.5 ml) at $\leq -70^{\circ}\text{C}$.

b. Master Plates

Master plates were prepared by streaking each tester strain from a frozen permanent stock onto minimal agar appropriately supplemented with 1) for *Salmonella typhimurium*, an excess of histidine, and biotin, and for tester strains TA98 and TA100, ampicillin (25 µg/ml), to ensure the stable maintenance of the pKM101 plasmid; and 2) for *Escherichia coli*, an excess of tryptophan. Tester strain master plates were stored at $5 \pm 3^{\circ}\text{C}$.

5. Preparation of Overnight Cultures

a. Inoculation

Overnight cultures for use in all testing procedures were inoculated by transferring a colony from the appropriate master plate to a flask containing culture medium. Inoculated flasks were placed in a shaker/incubator which was programmed to begin operation (shaking, 125 ± 25 rpm; incubation, $37 \pm 2^{\circ}\text{C}$) so that the overnight cultures were in log phase or late log phase when turbidity monitoring began.

b. Harvest

To ensure that cultures were harvested in late log phase, the length of incubation was determined by spectrophotometric monitoring of culture turbidity. Cultures were harvested once a predetermined turbidity was reached as determined by a percent transmittance (%T) reading on a spectrophotometer. This target turbidity ensures that cultures have reached a density of at least 0.5×10^9 cells per ml and that the cultures have not overgrown. Overgrown (stationary) cultures may exhibit decreased sensitivity to some mutagens. Cultures were removed from incubation when the target %T was reached and were placed at $5 \pm 3^{\circ}\text{C}$.

6. Confirmation of Tester Strain Genotypes

Tester strain cultures were checked for the following genetic markers on the day of their use in the mutagenicity assay:

a. *Salmonella typhimurium*

1) *rfa* Wall Mutation

The presence of the *rfa* wall mutation was confirmed by demonstration of the sensitivity of the culture to crystal violet. An aliquot of an overnight

culture of each strain was overlaid onto plates containing selective media and an antibiotic sensitivity disk containing 10 µg of crystal violet was added. Sensitivity was demonstrated by inhibition of bacterial growth in a zone immediately surrounding the disk.

2) pKM101 Plasmid R-factor

The presence of the pKM101 plasmid was confirmed for tester strains TA98 and TA100 by demonstration of resistance to ampicillin. An aliquot of an overnight culture of each strain was overlaid onto plates containing selective media and an antibiotic sensitivity disk containing 10 µg of ampicillin was added. Resistance was demonstrated by bacterial growth in the zone immediately surrounding the disk.

3) Characteristic Number of Spontaneous Revertants

The mean number of spontaneous revertants per plate in the vehicle controls that are characteristic of the respective strains were demonstrated by plating 100 µl aliquots of the culture along with the appropriate vehicle on selective media.

b. *Escherichia coli*

1) Characteristic Number of Spontaneous Revertants

The mean number of spontaneous revertants per plate in the vehicle controls that are characteristic of the respective strains were demonstrated by plating 100 µl aliquots of the WP2*uvrA* culture along with the appropriate vehicle on selective media.

7. Tester Strain Media

a. Culturing Broth

The broth used to grow overnight cultures of the tester strains was Vogel-Bonner salt solution (Vogel and Bonner, 1956) supplemented with 2.5% (w/v) Oxoid Nutrient Broth No. 2 (dry powder).

b. Agar Plates

Bottom agar (25 ml per 15 x 100 mm petri dish) was Vogel-Bonner minimal medium E (Vogel and Bonner, 1956), supplemented with 1.5% (w/v) agar and 0.2% (w/v) glucose.

c. Overlay Agar for Selection of Revertants

Overlay (top) agar was prepared with 0.7% agar (w/v) and 0.5% NaCl (w/v) and was supplemented with 10 ml of 1) 0.5 mM histidine/biotin solution per 100 ml agar for selection of histidine revertants, or 2) 0.5 mM tryptophan solution per 100 ml of agar for selection of tryptophan revertants. When S9 mix was required, 2.0 ml of the supplemented top agar was used in the overlay. However, when S9 mix was not required, water was added to the supplemented top agar (0.5 ml of water per 2 ml of supplemented top agar) and the resulting 2.5 ml of diluted supplemented top agar was used for the overlay. This dilution ensured that the final top agar and amino acid supplement concentrations remained the same both in the presence and absence of S9 mix.

B. Liver Microsomal Enzyme Reaction Mixture (S9 Mix)1. S9 Homogenate

Liver microsomal enzymes (S9 homogenate) were purchased from Molecular Toxicology, Inc., Annapolis, MD 21401, Batch 0642 (39.5 mg of protein per ml) and Batch 0646 (43.4 mg of protein per ml). The homogenate was prepared from male Sprague-Dawley rats that had been injected (i.p.) with Aroclor™ 1254 (200 mg per ml in corn oil) at 500 mg/kg as described by Ames *et al*, 1975.

2. S9 Mix

The S9 mix was prepared immediately prior to its use in any experimental procedure. The S9 mix contained the components indicated in Table II.

TABLE II. S9 MIX COMPONENTS

H ₂ O	0.70 ml
1M NaH ₂ PO ₄ /Na ₂ HPO ₄ , pH 7.4	0.10 ml
0.25M Glucose-6-phosphate	0.02 ml
0.10M NADP	0.04 ml
0.825M KCl/0.2M MgCl ₂	0.04 ml
S9 Homogenate	<u>0.10 ml</u>
	1.00 ml

C. Controls1. Vehicle Controls

Vehicle controls were plated for all tester strains both in the presence and absence of S9 mix. The vehicle control was plated, using a 50 µl aliquot of vehicle (equal to the maximum aliquot of test article dilution plated), along with a 100 µl aliquot of the appropriate tester strain and a 500 µl aliquot of S9 mix (when necessary), on selective agar.

2. Positive Controls

The combinations of positive controls, activation condition and tester strains plated concurrently with the assay are indicated in Table III.

TABLE III. POSITIVE CONTROLS

Tester Strain	S9 Mix	Positive Control	Conc per plate
TA98	+	2-aminoanthracene	2.5 µg
TA98	-	2-nitrofluorene	1.0 µg
TA100	+	2-aminoanthracene	2.5 µg
TA100	-	sodium azide	2.0 µg
TA1535	+	2-aminoanthracene	2.5 µg
TA1535	-	sodium azide	2.0 µg
TA1537	+	2-aminoanthracene	2.5 µg
TA1537	-	ICR-191	2.0 µg
WP2uvrA	+	2-aminoanthracene	25.0 µg
WP2uvrA	-	4-nitroquinoline-N-oxide	1.0 µg

a. Source and Grade of Positive Control Articles

2-aminoanthracene (CAS #613-13-8), Sigma Chemical Co., purity ≥ 97.5%; 2-nitrofluorene (CAS #607-57-8), Aldrich Chemical Co., purity 98%; sodium azide (CAS #26628-22-8), Sigma Chemical Co., purity >98%; ICR-191 (CAS #1707-45-0), Polysciences Inc., purity >95% and Sigma Chemical Co., purity 98%; 4-nitroquinoline-N-oxide (CAS #56-57-5), Sigma Chemical Co., purity >99%.

3. Sterility Controls

a. Test Article

The most concentrated test article dilution was checked for sterility by plating a 50 µl aliquot (the same volume used in the assay) on selective agar.

b. S9 Mix

The S9 mix was checked for sterility by plating 0.5 ml on selective agar.

METHODS

A. Dose Rangefinding Study

The growth inhibitory effect (cytotoxicity) of the test article to the test system was determined in order to allow the selection of appropriate doses to be tested in the mutagenicity assay.

1. Design

The dose rangefinding study was performed using tester strains TA100 and WP2uvrA both in the presence and absence of S9 mix. Ten doses of test article were tested at one plate per dose. The test article was checked for cytotoxicity up to a maximum concentration of 5 mg per plate.

a. Rationale

The cytotoxicity of the test article observed on tester strain TA100 is generally representative of that observed on the other tester strains and because of TA100's comparatively high number of spontaneous revertants per plate, gradations of cytotoxicity can be readily discerned from routine experimental variation. The *Escherichia coli* tester strain WP2uvrA does not possess the *rfa* wall mutation that the *Salmonella typhimurium* strains have and thus, a different range of cytotoxicity may be observed. Also, the cytotoxicity induced by a test article in the presence of S9 mix may vary greatly from that observed in the absence of S9 mix. Therefore, this would require that different test article dose ranges be tested in the mutagenicity assay based on the presence or absence of the S9 mix.

2. Evaluation of the Dose Ranging Study

Cytotoxicity is detectable as a decrease in the number of revertant colonies per plate and/or by a thinning or disappearance of the bacterial background lawn.

3. Selection of the Maximum Dose for the Mutagenicity Assay

a. No Cytotoxicity Observed

Since no cytotoxicity was observed in the dose ranging study, the highest dose level of test article used in the mutagenicity assay was the same as that tested in the ranging study.

B. Mutagenicity Assay

1. Design

The assay was performed using tester strains TA98, TA100, TA1535, TA1537, and WP2uvrA both in the presence and absence of S9 mix. Five doses of test article were tested along with the appropriate vehicle and positive controls. The doses of test article were selected based on the results of the dose ranging study. The results of the initial mutagenicity assay were confirmed in an independent experiment.

2. Frequency and Route of Administration

The tester strains were exposed to the test article via the plate incorporation methodology originally described by Ames *et al* (1975) and Maron and Ames (1983). This methodology has been shown to detect a wide range of classes of chemical mutagens. In the plate incorporation methodology, the test article, the tester strain and the S9 mix (where appropriate) were combined in molten agar which was overlaid onto a minimal agar plate. Following incubation at $37 \pm 2^\circ\text{C}$ for 48 ± 8 hr, revertant colonies were counted. All doses of the test article, the vehicle controls and the positive controls were plated in triplicate.

C. Plating Procedures

These procedures were used in both the dose ranging study and the mutagenicity assay.

Each plate was labeled with a code which identified the test article, test phase, tester strain, activation condition and dose level. The S9 mix and dilutions of the test article were prepared immediately prior to their use.

When S9 mix was not required, 100 μ l of tester strain and 50 μ l of vehicle or test article dose was added to 2.5 ml of molten selective top agar (maintained at $45 \pm 2^\circ\text{C}$). When S9 mix was required, 500 μ l of S9 mix, 100 μ l of tester strain and 50 μ l of vehicle or test article dose was added to 2.0 ml of molten selective top agar. After the required components had been added, the mixture was vortexed and overlaid onto the surface of 25 ml of minimal bottom agar contained in a 15 x 100 mm petri dish. After the overlay had solidified, the plates were inverted and incubated for 48 ± 8 hr at $37 \pm 2^\circ\text{C}$. Positive control articles were plated using a 50 μ l plating aliquot.

D. Scoring the Plates

Plates which were not evaluated immediately following the incubation period were held at $5 \pm 3^\circ\text{C}$ until such time that colony counting and bacterial background lawn evaluation could take place.

1. Bacterial Background Lawn Evaluation

The condition of the bacterial background lawn was evaluated for evidence of cytotoxicity and test article precipitate. Evidence of cytotoxicity was scored relative to the vehicle control plate and was recorded along with the revertant counts for all plates at that dose level on the data tables using the code system presented at the end of the Materials and Methods Section.

2. Counting Revertant Colonies

The number of revertant colonies per plate for the vehicle controls and all plates containing test article were counted manually. The number of revertant colonies per plate for the positive controls were counted by automated colony counter with the exception of the positive controls for tester strain WP2uvrA in the absence of S9 mix, in Experiment 17750-B1, which were counted manually.

E. Analysis of Data

For all replicate platings, the mean revertants per plate and the standard deviation were calculated. The results of these calculations are presented in tabular form in the Data Tables Section of this report.

EVALUATION OF TEST RESULTS

Before assay data were evaluated, the criteria for a valid assay had to be met.

A. Criteria For A Valid Assay

The following criteria were used to determine a valid assay:

1. Tester Strain Integrity : *Salmonella typhimurium*

a. *rfa* Wall Mutation

To demonstrate the presence of the *rfa* wall mutation, tester strain cultures exhibited sensitivity to crystal violet.

b. pKM101 Plasmid

To demonstrate the presence of the R-factor plasmid, pKM101, cultures of tester strains TA98 and TA100 exhibited resistance to ampicillin.

c. Characteristic Number of Spontaneous Revertants

To demonstrate the requirement for histidine, the tester strain cultures exhibited a characteristic number of spontaneous revertants per plate when plated along with the vehicle under selective conditions. The acceptable ranges for the mean vehicle controls were as follows:

TA98	8 - 60
TA100	60 - 240
TA1535	4 - 45
TA1537	2 - 25

2. Tester Strain Integrity : *Escherichia coli*

a. Characteristic Number of Spontaneous Revertants

To demonstrate the requirement for tryptophan, the tester strain culture exhibited a characteristic number of spontaneous revertants per plate when plated along

with the vehicle under selective conditions. The acceptable range for the WP2uvrA mean vehicle controls was 5 to 40 revertants per plate.

3. Tester Strain Culture Density

To demonstrate that appropriate numbers of bacteria are plated, the density of tester strain cultures were greater than or equal to 0.5×10^9 bacteria per ml and/or had reached a target level of turbidity demonstrated to produce cultures with a density greater than or equal to 0.5×10^9 bacteria per ml.

4. Positive Control Values

a. Positive Control Values in the Absence of S9 Mix

To demonstrate that the tester strains were capable of identifying a mutagen, the mean value of a positive control for a respective tester strain exhibited at least a 3-fold increase over the mean value of the vehicle control for that strain.

b. Positive Control Values in the Presence of S9 Mix
(S9 Mix Integrity)

To demonstrate that the S9 mix was capable of metabolizing a promutagen to its mutagenic form(s), the mean value of the positive control for a respective tester strain in the presence of the S9 mix exhibited at least a 3-fold increase over the mean value of the vehicle control for that strain.

An acceptable positive control in the presence of S9 mix for a specific strain was evaluated as having demonstrated both the integrity of the S9 mix and the ability of the tester strain to detect a mutagen.

5. Cytotoxicity

A minimum of three non-toxic doses were required to evaluate assay data.

B. Criteria For A Positive Response

Once the criteria for a valid assay had been met, responses observed in the assay were evaluated as follows:

1. Tester Strains TA98, TA100, and WP2_{uvrA}

For a test article to be considered positive, it had to produce at least a 2-fold increase in the mean revertants per plate of at least one of these tester strains over the mean revertants per plate of the appropriate vehicle control. This increase in the mean number of revertants per plate had to be accompanied by a dose response to increasing concentrations of the test article.

2. Tester Strains TA1535 and TA1537

For a test article to be considered positive, it had to produce at least a 3-fold increase in the mean revertants per plate of at least one of these tester strains over the mean revertants per plate of the appropriate vehicle control. This increase in the mean number of revertants per plate had to be accompanied by a dose response to increasing concentrations of the test article.

RECORDS TO BE MAINTAINED

All raw data, documentation, records, the protocol, and the final report generated as a result of this study will be archived in the storage facilities of Corning Hazleton Inc. for at least one year following submission of the final report to the Sponsor. After the one year period, the Sponsor may elect to have the aforementioned materials retained in the storage facilities of Corning Hazleton Inc. for an additional period of time or sent to a storage facility designated by the Sponsor.

REFERENCES

- Ames, B.N., J. McCann, and E. Yamasaki. Methods for detecting carcinogens and mutagens with the *Salmonella*/Mammalian-Microsome Mutagenicity Test. *Mutation Research* 31:347-364 (1975).
- Brusick, D.J., V.F. Simmon, H.S. Rosenkranz, V.A. Ray, and R.S. Stafford.
An evaluation of the *Escherichia coli* WP2 and WP2_{uvrA} reverse mutation assay.
Mutation Research 76:169-190 (1980).
- Green, M.H.L. and W.J. Muriel. Mutagen testing using *trp*⁺ reversion in *Escherichia coli*.
Mutation Research 38:3-32. (1976).

CORNING Hazleton

Maron, D.M., and B. Ames. Revised methods for the *Salmonella* Mutagenicity Test. *Mutation Research* 113:173-215 (1983).

Vogel, H.J., and D.M. Bonner. Acetylornithinase of *E. coli*: Partial purification and some properties. *J. Biol. Chem.* 218:97-106 (1956).

BACTERIAL BACKGROUND LAWN EVALUATION CODE

The condition of the background bacterial lawn is evaluated both macroscopically and microscopically (using a dissecting microscope) for indications of cytotoxicity and test article precipitate as follows:

CODE	DEFINITION	<u>CHARACTERISTICS OF BACKGROUND LAWN</u>
1	Normal	A healthy microcolony lawn.
2	Slightly Reduced	A noticeable thinning of the microcolony lawn and an increase in the size of the microcolonies compared to the vehicle control plate.
3	Moderately Reduced	A marked thinning of the microcolony lawn and an increase in the size of the microcolonies compared to the vehicle control plate.
4	Extremely Reduced	An extreme thinning of the microcolony lawn and an increase in the size of the microcolonies compared to the vehicle control plate.
5	Absent	A complete lack of any microcolony lawn.
6	Obscured by Precipitate	The background bacterial lawn cannot be accurately evaluated due to microscopic and/or macroscopic test article precipitate.

Evidence of macroscopic test article precipitate on the plates is recorded by addition of the following precipitate code to the code number used to evaluate the condition of the background bacterial lawn.

sp	Slight Precipitate	Noticeable macroscopic precipitate on the plate, however, the precipitate does not influence automated counting of the plate.
mp	Moderate Precipitate	The amount of macroscopic precipitate on the plate would interfere with automated counting, thus requiring the plate to be hand counted.
hp	Heavy Precipitate	The large amount of macroscopic precipitate on the plate makes the required hand counting difficult.

Example: 4mp would indicate a plate observed to have an extremely reduced background lawn which had to be counted manually due to the marked amount of macroscopic test article precipitate.

SECTION IV. RESULTS AND CONCLUSIONS

RESULTS

A. Test Article Handling

The test article, T-6564, was stored at room temperature. Deionized water (CHV Lots 347, 349, and 352) was used as the vehicle. At 100 mg per ml, which was the most concentrated stock dilution prepared, the test article formed a clear, colorless solution. The test article remained a solution in all succeeding dilutions prepared for the mutagenicity assay.

B. Dose Rangefinding Study

Doses to be tested in the mutagenicity assay were selected based on the results of the dose rangefinding study conducted on the test article using tester strains TA100 and WP2uvrA in both the presence and absence of S9 mix (one plate per dose). Ten doses of test article, from 5,000 to 6.67 µg per plate, were tested and the results are presented in Tables 1 and 2. These data were generated in Experiment 17750-A1. No cytotoxicity was observed in either the presence or absence of S9 mix as evidenced by a normal background lawn and no decrease in the number of revertants per plate.

C. Mutagenicity Assay

The mutagenicity assay results for T-6564 are presented in Tables 3 through 7. These data were generated in Experiments 17750-B1, 17750-C1, and 17750-D1. The data are presented as mean revertants per plate ± standard deviation for each treatment and control group (Tables 4, 6, and 7) and as individual plate counts (Tables 3, 5, and 7).

The results of the dose rangefinding study were used to select five doses to be tested in the mutagenicity assay. The doses tested were 5,000, 3,330, 1,000, 333, and 100 µg per plate in both the presence and absence of S9 mix.

In the initial mutagenicity assay, Experiment 17750-B1 (Tables 3 and 4), and in the confirmatory assay, Experiment 17750-C1 (Tables 5 and 6), all data were acceptable and no positive increases in the number of revertants per plate were observed with any of the tester strains either in the presence or absence of S9 mix. In Experiment 17750-C1, a non-dose responsive 5.0-fold increase which did not meet the criteria for a positive response was observed with tester strain TA1537 in the absence of S9 mix. For this reason tester strain TA1537 was retested in the absence of S9 mix in Experiment 17750-D1.

In Experiment 17750-D1 (Table 7), all data were acceptable and no positive increases in the number of revertants per plate were observed with tester strain TA1537 in the absence of S9 mix.

All criteria for a valid study were met.

CHV Study No.: 17750-0-409R

CONCLUSIONS

The results of the *Salmonella - Escherichia coli*/Mammalian-Microsome Reverse Mutation Assay with a Confirmatory Assay indicate that, under the conditions of this study, 3M Corporation's test article, T-6564, did not cause a positive increase in the number of revertants per plate of any of the tester strains either in the presence or absence of microsomal enzymes prepared from AroclorTM-induced rat liver (S9).

CORNING Hazleton

SECTION V. DATA TABLES

CHV Study No.: 17750-0-409R

27

000258

CORNING Hazleton

TABLE 1

DOSE RANGE FINDING STUDY

TEST ARTICLE ID: T-6564

EXPERIMENT ID: 17750-A1

DATE PLATED: 12-Jun-96

VEHICLE: Deionized water

DATE COUNTED: 18-Jun-96

μ g/PLATE	TA100 REVERTANTS PER PLATE			
	WITH S9		WITHOUT S9	
	REVERTANTS PER PLATE	BACKGROUND LAWN EVALUATION*	REVERTANTS PER PLATE	BACKGROUND LAWN EVALUATION*
0.00 (Vehicle) (50 μ l)	85	1	78	1
Test Article 6.67	106	1	91	1
10.0	128	1	96	1
33.3	120	1	78	1
66.7	122	1	105	1
100	113	1	84	1
333	126	1	113	1
667	121	1	96	1
1000	110	1	86	1
3330	112	1	81	1
5000	118	1	99	1

* Background Lawn Evaluation Codes:

1 = normal

4 = extremely reduced

sp = slight precipitate

2 = slightly reduced

5 = absent

mp = moderate precipitate
(requires hand count)

3 = moderately reduced

6 = obscured by precipitate

hp = heavy precipitate
(requires hand count)

CORNING Hazleton

TABLE 2

DOSE RANGEFINDING STUDY

TEST ARTICLE ID: T-6564

EXPERIMENT ID: 17750-A1

DATE PLATED: 12-Jun-96

VEHICLE: Deionized water

DATE COUNTED: 18-Jun-96

µg/PLATE	WP2uvrA REVERTANTS PER PLATE			
	WITH S9		WITHOUT S9	
	REVERTANTS PER PLATE	BACKGROUND LAWN EVALUATION*	REVERTANTS PER PLATE	BACKGROUND LAWN EVALUATION*
0.00 (Vehicle) (50 µl)	10	1	12	1
Test Article 6.67	13	1	14	1
10.0	12	1	11	1
33.3	16	1	16	1
66.7	9	1	22	1
100	15	1	10	1
333	16	1	8	1
667	10	1	20	1
1000	9	1	14	1
3330	11	1	6	1
5000	10	1	5	1

* Background Lawn Evaluation Codes:

1 = normal

4 = extremely reduced

sp = slight precipitate

2 = slightly reduced

5 = absent

mp = moderate precipitate
(requires hand count)

3 = moderately reduced

6 = obscured by precipitate

hp = heavy precipitate
(requires hand count)

TABLE 3
MUTAGENICITY ASSAY RESULTS
INDIVIDUAL PLATE COUNTS

TEST ARTICLE ID: T-6564

EXPERIMENT ID: 17750-B1

DATE PLATED: 20-Jun-96

VEHICLE: Deionized water

DATE COUNTED: 28-Jun-96

PLATING ALIQUOT: 50 µl

REVERTANTS PER PLATE															BACKGROUND LAWN*		
DOSE/PLATE		TA98			TA100			TA1535			TA1537			WP2uvrA			
		1	2	3	1	2	3	1	2	3	1	2	3	1	2	3	
MICROSOMES: Rat Liver																	
VEHICLE CONTROL		25	22	13	118	106	103	15	17	9	9	10	5	15	12	12	1
TEST ARTICLE																	
100 µg		21	29	18	85	109	102	13	13	10	10	11	5	11	10	19	1
333 µg		16	27	25	101	104	96	8	11	15	5	6	2	20	24	9	1
1000 µg		27	16	29	106	100	101	10	12	15	8	5	12	14	14	14	1
3330 µg		25	17	19	114	110	125	18	10	13	12	8	10	11	16	8	1
5000 µg		23	25	25	117	99	109	12	12	10	7	7	5	13	13	5	1
POSITIVE CONTROL **		860	852	970	936	892	810	141	134	124	240	202	190	314	279	270	1
MICROSOMES: None																	
VEHICLE CONTROL		12	12	12	91	74	85	14	10	6	5	4	7	15	14	10	1
TEST ARTICLE																	
100 µg		14	10	13	81	104	83	11	9	5	6	7	1	8	7	13	1
333 µg		13	11	4	101	95	86	10	13	5	6	6	6	12	12	11	1
1000 µg		16	11	12	97	98	92	11	10	7	6	6	5	8	12	9	1
3330 µg		15	10	2	108	93	75	12	10	6	4	5	1	12	8	11	1
5000 µg		11	11	13	106	108	104	16	14	12	9	5	4	13	13	9	1
POSITIVE CONTROL ***		105	99	86	470	419	460	460	427	400	1233	918	1110	135	102	85	1

** TA98 2-aminoanthracene 2.5 µg/plate
TA100 2-aminoanthracene 2.5 µg/plate
TA1535 2-aminoanthracene 2.5 µg/plate
TA1537 2-aminoanthracene 2.5 µg/plate
WP2uvrA 2-aminoanthracene 25.0 µg/plate

*** TA98 2-nitrofluorene 1.0 µg/plate
TA100 sodium azide 2.0 µg/plate
TA1535 sodium azide 2.0 µg/plate
TA1537 ICR-191 2.0 µg/plate
WP2uvrA 4-nitroquinoline-N-oxide 1.0 µg/plate

* Background Lawn Evaluation Codes:

1 = normal
2 = slightly reduced
3 = moderately reduced
4 = extremely reduced
5 = absent
6 = obscured by precipitate
sp = slight precipitate
mp = moderate precipitate
(requires hand count)
hp = heavy precipitate
(requires hand count)

TABLE 4
MUTAGENICITY ASSAY RESULTS
SUMMARY

TEST ARTICLE ID: T-6564

EXPERIMENT ID: 17750-B1

DATE PLATED: 20-Jun-96

VEHICLE: Deionized water

DATE COUNTED: 28-Jun-96

PLATING ALIQUOT: 50 µl

MEAN REVERTANTS PER PLATE WITH STANDARD DEVIATION												BACKGROUND LAWN*
	DOSE/PLATE	TA98		TA100		TA1535		TA1537		WP2uvrA		
		MEAN	S.D.	MEAN	S.D.	MEAN	S.D.	MEAN	S.D.	MEAN	S.D.	
MICROSOMES: Rat Liver												
VEHICLE CONTROL		20	6	109	8	14	4	8	3	13	2	1
TEST ARTICLE	100 µg	23	6	99	12	12	2	9	3	13	5	1
	333 µg	23	6	100	4	11	4	4	2	18	8	1
	1000 µg	24	7	102	3	12	3	8	4	14	0	1
	3330 µg	20	4	116	8	14	4	10	2	12	4	1
	5000 µg	24	1	108	9	11	1	6	1	10	5	1
POSITIVE CONTROL **		894	66	879	64	133	9	211	26	288	23	1
MICROSOMES: None												
VEHICLE CONTROL		12	0	83	9	10	4	5	2	13	3	1
TEST ARTICLE	100 µg	12	2	89	13	8	3	5	3	9	3	1
	333 µg	9	5	94	8	9	4	6	0	12	1	1
	1000 µg	13	3	96	3	9	2	6	1	10	2	1
	3330 µg	9	7	92	17	9	3	3	2	10	2	1
	5000 µg	12	1	106	2	14	2	6	3	12	2	1
POSITIVE CONTROL ***		97	10	450	27	429	30	1087	159	107	25	1

** TA98 2-aminoanthracene 2.5 µg/plate
TA100 2-aminoanthracene 2.5 µg/plate
TA1535 2-aminoanthracene 2.5 µg/plate
TA1537 2-aminoanthracene 2.5 µg/plate
WP2uvrA 2-aminoanthracene 25.0 µg/plate

*** TA98 2-nitrofluorene 1.0 µg/plate
TA100 sodium azide 2.0 µg/plate
TA1535 sodium azide 2.0 µg/plate
TA1537 ICR-191 2.0 µg/plate
WP2uvrA 4-nitroquinoline-N-oxide 1.0 µg/plate

* Background Lawn Evaluation Codes:

1 = normal
2 = slightly reduced
3 = moderately reduced
4 = extremely reduced
5 = absent
6 = obscured by precipitate
sp = slight precipitate
mp = moderate precipitate
(requires hand count)
hp = heavy precipitate
(requires hand count)

TABLE 5
MUTAGENICITY ASSAY RESULTS
INDIVIDUAL PLATE COUNTS

TEST ARTICLE ID: T-6564

EXPERIMENT ID: 17750-C1

DATE PLATED: 25-Jul-96

VEHICLE: Deionized water

DATE COUNTED: 30-Jul-96, 31-Jul-96

PLATING ALIQUOT: 50 µl

		REVERTANTS PER PLATE												BACKGROUND LAWN*				
DOSE/PLATE		TA98			TA100			TA1535			TA1537			WP2uvrA				
		1	2	3	1	2	3	1	2	3	1	2	3	1	2	3		
MICROSOMES: Rat Liver																		
VEHICLE CONTROL		25	25	18	94	76	103	17	14	17	4	7	6	12	14	18	1	
TEST ARTICLE		100 µg	21	21	28	89	99	95	8	11	9	11	4	11	14	16	10	1
	333 µg	27	33	25	102	105	76	16	16	9	8	6	19	10	15	9	1	
	1000 µg	33	26	29	87	115	108	12	10	8	7	7	10	24	13	14	1	
	3330 µg	24	30	27	106	96	114	9	18	12	8	4	4	24	16	13	1	
	5000 µg	25	33	25	107	87	89	8	6	17	8	7	8	21	14	20	1	
POSITIVE CONTROL **		769	624	815	949	1108	1066	146	196	136	89	129	91	151	90	108	1	
MICROSOMES: None																		
VEHICLE CONTROL		13	18	15	79	89	70	16	9	19	1	4	2	13	16	23	1	
TEST ARTICLE		100 µg	16	17	17	71	105	76	16	17	17	6	5	4	15	10	10	1
	333 µg	21	8	16	74	89	73	14	9	8	6	6	8	15	15	19	1	
	1000 µg	20	9	8	74	71	88	9	9	8	1	6	7	6	13	19	1	
	3330 µg	21	15	11	78	75	78	15	14	10	16	10	3	13	9	6	1	
	5000 µg	10	17	17	99	80	77	16	13	16	8	7	1	9	17	10	1	
POSITIVE CONTROL ***		207	208	172	548	597	595	644	619	663	462	477	463	187	207	160	1	

** TA98 2-aminoanthracene 2.5 µg/plate
TA100 2-aminoanthracene 2.5 µg/plate
TA1535 2-aminoanthracene 2.5 µg/plate
TA1537 2-aminoanthracene 2.5 µg/plate
WP2uvrA 2-aminoanthracene 25.0 µg/plate

*** TA98 2-nitrofluorene 1.0 µg/plate
TA100 sodium azide 2.0 µg/plate
TA1535 sodium azide 2.0 µg/plate
TA1537 ICR-191 2.0 µg/plate
WP2uvrA 4-nitroquinoline-N-oxide 1.0 µg/plate

* Background Lawn Evaluation Codes:

1 = normal
2 = slightly reduced
4 = extremely reduced
5 = absent
sp = slight precipitate
mp = moderate precipitate
(requires hand count)

3 = moderately reduced
6 = obscured by precipitate
hp = heavy precipitate
(requires hand count)

TABLE 6
MUTAGENICITY ASSAY RESULTS
SUMMARY

TEST ARTICLE ID: T-6564

EXPERIMENT ID: 17750-C1

DATE PLATED: 25-Jul-96

VEHICLE: Deionized water

DATE COUNTED: 30-Jul-96, 31-Jul-96

PLATING ALIQUOT: 50 µl

MEAN REVERTANTS PER PLATE WITH STANDARD DEVIATION											BACKGROUND LAWN*	
	DOSE/PLATE	TA98		TA100		TA1535		TA1537		WP2uvrA		
		MEAN	S.D.	MEAN	S.D.	MEAN	S.D.	MEAN	S.D.	MEAN	S.D.	
MICROSOMES: Rat Liver												
VEHICLE CONTROL		23	4	91	14	16	2	6	2	15	3	1
TEST ARTICLE	100 µg	23	4	94	5	9	2	9	4	13	3	1
	333 µg	28	4	94	16	14	4	11	7	11	3	1
	1000 µg	29	4	103	15	10	2	8	2	17	6	1
	3330 µg	27	3	105	9	13	5	5	2	18	6	1
	5000 µg	28	5	94	11	10	6	8	1	18	4	1
POSITIVE CONTROL **		736	100	1041	82	160	32	103	23	116	31	1
MICROSOMES: None												
VEHICLE CONTROL		15	3	79	10	15	5	2	2	17	5	1
TEST ARTICLE	100 µg	17	1	84	18	17	1	5	1	12	3	1
	333 µg	15	7	79	9	10	3	7	1	16	2	1
	1000 µg	12	7	78	9	9	1	5	3	13	7	1
	3330 µg	16	5	77	2	13	3	10	7	9	4	1
	5000 µg	15	4	85	12	15	2	5	4	12	4	1
POSITIVE CONTROL ***		196	21	580	28	642	22	467	8	185	24	1

** TA98 2-aminoanthracene 2.5 µg/plate
TA100 2-aminoanthracene 2.5 µg/plate
TA1535 2-aminoanthracene 2.5 µg/plate
TA1537 2-aminoanthracene 2.5 µg/plate
WP2uvrA 2-aminoanthracene 25.0 µg/plate

*** TA98 2-nitrofluorene 1.0 µg/plate
TA100 sodium azide 2.0 µg/plate
TA1535 sodium azide 2.0 µg/plate
TA1537 ICR-191 2.0 µg/plate
WP2uvrA 4-nitroquinoline-N-oxide 1.0 µg/plate

* Background Lawn Evaluation Codes:

1 = normal
2 = slightly reduced
4 = extremely reduced
5 = absent
sp = slight precipitate
mp = moderate precipitate
(requires hand count)

3 = moderately reduced
6 = obscured by precipitate
hp = heavy precipitate
(requires hand count)

TABLE 7
MUTAGENICITY ASSAY RESULTS
INDIVIDUAL PLATE COUNTS AND SUMMARY

TEST ARTICLE ID: T-6564

EXPERIMENT ID: 17750-D1

DATE PLATED: 25-Aug-96

VEHICLE: Deionized water

DATE COUNTED: 27-Aug-96

PLATING ALIQUOT: 50 μ l

	DOSE/PLATE	REVERTANTS PER PLATE TA1537			MEAN REVERTANTS PER PLATE WITH STANDARD DEVIATION TA1537		BACKGROUND LAWN*
		1	2	3	MEAN	S.D.	
MICROSOMES: None							
VEHICLE CONTROL		10	2	9	7	4	1
TEST ARTICLE	100 μ g	9	6	5	7	2	1
	333 μ g	9	5	4	6	3	1
	1000 μ g	5	6	8	6	2	1
	3330 μ g	6	8	5	6	2	1
	5000 μ g	3	1	9	4	4	1
POSITIVE CONTROL ***		128	63	60	84	38	1

*** TA1537 ICR-191

2.0 μ g/plate

* Background Lawn Evaluation Codes:

1 = normal	2 = slightly reduced	3 = moderately reduced
4 = extremely reduced	5 = absent	6 = obscured by precipitate
sp = slight precipitate	mp = moderate precipitate (requires hand count)	hp = heavy precipitate (requires hand count)