

FILE NAME: Colgate (COL)

DATE: 1975 Aug

DOC#: COL037

DOCUMENT DESCRIPTION: National Safety News Article - Microscopical
Identification of Target Health Hazards

NATIONAL SAFETY COUNCIL

A NON-GOVERNMENTAL, NON-PROFIT, PUBLIC SERVICE ORGANIZATION

Headquarters Office: 425 N. Michigan Ave.
Chicago, IL 60611

Trustees' Fund Office: 856 Lincoln Bldg.
60 E. 42nd St.
New York, NY 10017



National SAFETY NEWS

A NATIONAL SAFETY COUNCIL PUBLICATION

CHAIRMAN, BOARD OF DIRECTORS
DR. DEWEY F. BARICH, President
Detroit Institute of Technology, Detroit, MI 48201

PRESIDENT
VINCENT L. TOFANY, National Safety Council

VICE-PRESIDENT FOR INDUSTRY
MICHAEL F. BIANCARDI, Assistant Vice-President
Employers Insurance of Wausau, Wausau, WI 54401

ACTING EXECUTIVE VICE-PRESIDENT
FRED POTENZA, National Safety Council

SECRETARY AND TREASURER
H. W. CHAMPLIN, National Safety Council

PRESIDENT EMERITUS
HOWARD PYLE, National Safety Council



DEPARTMENTS

- AUGUST 1975 VOLUME 112, NUMBER 2
55,000 copies of this issue were printed
- ### FEATURED THIS MONTH
- 62 Public Hearings on Proposed Noise Standard
 - 64 Programming Personal Protection—
To Hear or Not To Hear (It's That Simple)
 - 67 Noise Control . . . Diesel Compressors
 - 72 Industrial Hygiene Instrumentation—Vinyl Chloride Monomer
 - 80 OSHA Issues Guidelines for Vinyl Chloride Inspections
 - 81 Vinyl Chloride—Monitoring Exposures
 - 84 Microscopical Identification of Target Health Hazards
 - 93 Bulk First Aid Kits (Data Sheet 651)
 - 96 Specialized Safety Manual
 - 98 NSC Proclams . . . With Safety For All
 - 100 NSC Safety Training Institute Schedule for the Bicentennial
 - 106 Emphasis on: Occupational Health and Protection Before
Compliance

Future Article Preview	4
Officially OSHA	6
Voice of the Reader	8
Safetybriefs	16
Coming Events	28
In Focus	36
CRISP	42
Standards	48
Product Safety Memo	52
With Safety For All	56
Congress Registration	57
OSHA-Gram	59
Ideas That Worked	86
Off the Job	88
Safetyclips	91
Safety Literature	108
Keeping Posted	111
Change of Address Card	115
Wire from Washington	150
Safety Library	170
Newscope	178
Obituary	182
Personnel Directory	183
New Safety Equipment	184
OSHC	186
Reprints	193
Advertisers' Index	196
Reader Service Card	198

ON THE COVER

A model of the human ear, used at Northwestern University's Diagnostic Center for Hearing Impaired, was photographed partially assembled by James B. Lehman, NSC staff photographer. It calls attention to "Hearing Protection"—what the recent federal hearings on the proposed noise standard were all about. (See page 62.)

EDITORIAL STAFF

editor	ROY FISHER
associate editors	JAMES W. LAHEY SUSAN M. SHERERDIK EDWARD D. DIONNE
contributing editors	DOROTHEA GUTHRIE RUTH HAMMERSMITH TOM HIRSH CHARLES PRICE HARRY N. ROSENFELD ROBERT WATERBURY
data sheet	FRANK E. McELROY JAMES WEGING
technical director	ROBERT BELKNAP
editorial director	ROBERT L. MEYER
publications manager	JACK HORNER
art director	FRANK WASZAK
photographer	JAMES B. LEHMAN
advertising manager	JOHN C. LUNDGREN
production manager	PATRICIA MAZZUCCA
reader card inquiries	VICKI BACSKAI

STATEMENTS AND OPINIONS advanced in articles and advertisements are personal expressions of the authors and advertisers, and they are not necessarily those of the National Safety Council . . . unless they are presented as policy statements and are so stated.

INFORMATION AND RECOMMENDATIONS contained in this publication have been compiled from sources believed to be reliable and to be representative of the best current opinion on the subject. No warranty, guarantee, nor representation is made by the National Safety Council as to the absolute correctness or sufficiency of any representation contained in this and other Council publications, and the National Safety Council assumes no responsibility in connection therewith. Likewise, it can not be assumed that all acceptable safety measures are contained in this and in other Council publications, nor that other or additional measures may not be required under particular or exceptional conditions or circumstances.

NATIONAL SAFETY NEWS is published monthly by the National Safety Council, and it is printed in the U.S.A. Copyright © 1975 by the National Safety Council. No material from any editorial feature or department, unless otherwise stated, may be reprinted without written permission from the editor.

SUBSCRIPTION RATES: U. S., \$9.25 a year (single copies, \$1.42); Foreign: no additional charge for Canada; add \$1.16 per subscription for postage to Latin American countries, and add \$1.79 per subscription for postage to all other countries. All prices (except postage) are subject to a 20 per cent discount to members of the National Safety Council, and an 8 per cent discount to U.S. Federal Government agencies. Quantity prices for yearly subscriptions and single issues available on request.

MEMBER: Audit Bureau of Circulation



NATIONAL SAFETY NEWS is entered as second class postage paid at Chicago, IL, and at additional mailing offices.

POSTMASTER: Send Form 3572 to National Safety Council, 425 N. Michigan Ave., Chicago, IL 60611

75 2848

Microscopical Identification Of Target Health Hazards

Through sophisticated techniques and instruments, scientists can readily identify asbestos, lead, silica, cotton dust and other air-borne particulates in the workplace.

INDUSTRIAL HYGIENE is concerned with the study of hazards to the health of the industrial worker. Some of these hazards are tiny particles in the air we breathe. They are hard to trap and identify, and difficult to associate directly with the failure of health because many years may elapse after exposure before adverse effects occur.

In any case, we know we have problems. Focusing upon the need to create healthful working conditions, OSHA, in January, 1972, initiated the Target Health Hazards Program. The emphasis is on five hazardous workplace substances: asbestos, lead, silica, cotton dust, and carbon monoxide.

Asbestosis resulting from prolonged inhalation of respirable air-borne asbestos dust may appear 20, 30, or 40 years after the exposure. Other particulate hazards are cotton dust, lead compounds, and quartz. Unfortunately, air-borne particles of these materials are so ubiquitous that everyone from factory workers to housewives and school children is exposed to some extent.

Exposure to particulate lead compounds is common due to the

use of tetraethyl lead in gasoline and in paint.

Asbestos is used in filters, insulation, fire-retardant textiles, and brake linings. Abrasion of these materials produces particles that find their way into the atmosphere.

Quartz, one form of which is silica, is the most common mineral on earth and few air samples are free of small particles of the type that may eventually cause silicosis.

Cotton dust, though less of a hazard than the others, is still a potent danger to industrial workers exposed for long periods of time to high cotton dust concentrations.

The exposure time to quartz and asbestos particles must be long-term, and the dosage high before serious impairment occurs. Our lungs also have a reasonably good flushing mechanism that helps to remove the billions of particles that settle every day on those surfaces. Perhaps this capacity for cleansing the lungs must be overwhelmed for some minimal time before permanent serious damage is done.

Many varieties of technical backgrounds must be involved before the solutions are found. The microscopist has one important role—



By Walter C. McCrone, D. Sc., Scientific Advisor, Walter C. McCrone Associates, Inc., Chicago.

among others—monitoring dust and other samples for the presence and amount of the hazardous substances. He is able to do this because the microscope magnifies each particle to a size such that selective



Figure 1 shows quartz (silica) particles magnified 125 times with crossed polars.

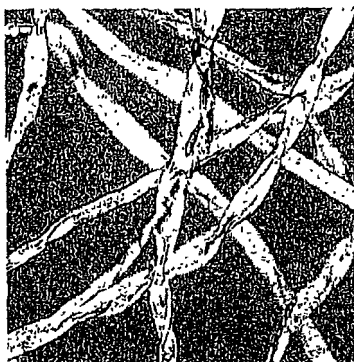


Figure 2 represents cotton fibers as they appear with crossed polars at 200 times magnification.

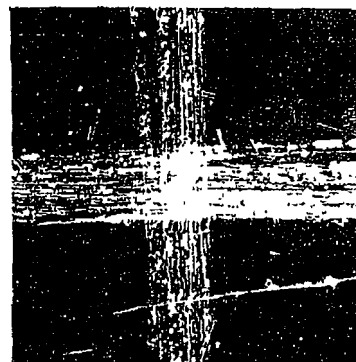


Figure 3 shows chrysotile fibers (one form of asbestos) at 100 times magnification.



Figure 4 shows electron microprobe analyzer used at McCrone Associates.

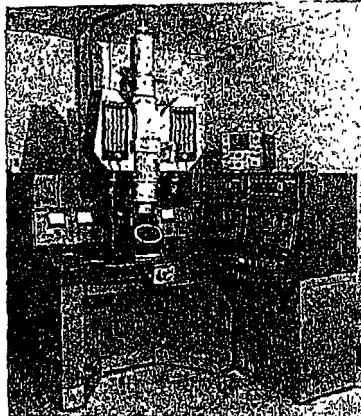


Figure 5 shows the AEI EMMA 4 electron microscope/electron probe analyzer.

and specific physical measurements can be made to identify that particle. By selecting the proper microscope and utilizing his background in crystallography the microscopist can identify many particulate substances at sight. These include most of the air-borne particles detrimental to our health—quartz, cotton fiber dust, and asbestos. Lead is a special problem because there are many different lead compounds.

Carbon monoxide is not detectable by microscopy. Refer to the October 1973 NATIONAL SAFETY News article "Target Health Hazard—Carbon Monoxide," for information regarding CO and the various means of detecting it.

First, consider the problem of identifying quartz microscopically. The shape and optical properties make this easy. Figure 1 shows quartz between two slightly uncrossed polarizing filters; the low birefringence (low numerical difference between the maximum and minimum refractive indices) and conchoidal (shell-like) fracture are obvious. To confirm his initial identification of quartz the microscopist has an optical staining procedure (called dispersion staining) that causes all quartz particles to appear blue or blue magenta. The latter procedure is excellent for detecting traces of very fine quartz in dust samples.

Next, consider cotton fibers. These, too, have characteristic optical properties and shape, and can be identified just as dependably as quartz (see Figure 2).

Asbestos presents a somewhat more complicated problem because asbestos is not a single substance, but a class of fibrous minerals of

widely varying composition and properties. Even here, however, the microscopist is able to vary his procedure so that fibers observed microscopically can be further identified as chrysotile, antigorite, anthophyllite, tremolite, amosite, or crocidolite. Each has unique optical properties observed as optical staining colors when properly mounted for examination. Figure 3 shows one of four different asbestiform minerals, chrysotile. Again, individual small fibers can be identified even when present as trace quantities in a complex mixture.

Finally, we should consider lead compounds. As suggested, these present problems because they are not a specific compound as is quartz, nor even a small group of slightly different compounds like asbestos. There is, however, one simple procedure that identifies all lead-containing particles as such, though it does not specifically identify the compound. The dust is dispersed (allowed to settle due to gravity) on a reagent-impregnated, gelatin-coated microscope slide. The reagent—sodium chromate, for example—reacts with lead to form a colored halo around each particle containing lead. One then simply counts the halos microscopically to determine the number of lead particles.

Another procedure has been used to identify the lead compounds in auto exhaust. The particles spread over the surface of a highly polished beryllium plate were examined by an electron microprobe analyzer as shown in Figure 4. This instrument, computer-controlled, systematically scans the particles

with a finely focused pencil of electrons. Each time the beam strikes a tiny particle X-rays are created and emitted. If these X-rays show the proper wavelength for lead the computer stops the scanning sequence until that lead-containing speck can be analyzed completely and quantitatively. This requires less than one minute for each particle. The scanning then proceeds to the next lead particle. In one two-week period of continuous operation the probe found and analyzed over 12,000 individual lead particles.

The problem presented by water-borne samples is one of decreased particle sizes; most large particles will have settled out. Particles smaller than about a few tenths of a micrometer (10 millionths of an inch, for example) cannot usually be identified by the methods described. The microscopist can, however, identify these small particles by electron microscopy and its auxiliary techniques. The most sophisticated of these tools is the A.E.I. Scientific Apparatus Ltd. EMMA (Electron Microscope MicroAnalyser) (See Figure 5) It is a combined electron microscope and microprobe analyzer. The identification of asbestos in water supplies is a good example of the application of an EMMA. A sample of water is filtered through the finest membrane filter. A few square millimeters of the filter are cut out and placed on a carbon-coated electron microscope grid. The filter itself is dissolved in a solvent vapor system, and the particles are left supported on the thin carbon film. The grid is then placed in the EMMA and scanned visually at about 20,000 times for fibers. Each suspected asbestos fiber is then analyzed chemically *in situ* by the microprobe and by selected area electron diffraction (identification of elements from the wavelengths of their electrons). This determines its crystal structure to confirm asbestos and identify precisely the mineral type.

Microscopy is a sophisticated technique for analyzing particles hazardous to health. A well-trained microscopist with the necessary equipment can quickly and positively identify and characterize suspect particles.—End.