



E. I. DU PONT DE NEMOURS & COMPANY
INCORPORATED
WILMINGTON, DELAWARE 19898

POLYMER PRODUCTS DEPARTMENT
EXPERIMENTAL STATION

PERSONAL AND CONFIDENTIAL

cc: A. J. Dahl - 353
B. W. Karrh - N11400
L. J. Papa - 269
Pral File
I.C.

AR 226-1578

May 20, 1981

TO: DR. C. WALKER - PPD, Chestnut Run

FROM: S. S. STAFFORD *S. Stafford*

ANALYSIS OF BLOOD SAMPLES FOR PERFLUOROOCTANOATE
(Job No. 811-610; PRAL Nos. 81-1763-1777; Notebook Nos. E22514, E26238)

As requested in your letter of 4/23/81 to L. J. Papa, the 15 blood samples submitted then have been analyzed for perfluorooctanoate (C₈). Results and sample identification are given in the attached table.

As noted there, the analyses were done using a gas chromatographic method specific for C₈ (Lab Method Number ES-567) but results have been reported as ppm F for comparison with total organic fluorine analyses. Precision is $\pm 10\%$ relative standard deviation over most of the concentration range, somewhat less at the lowest values. The lower limit for quantitation is 0.007 ppm F (0.01 ppm perfluorooctanoic acid), with a detection limit of ~ 0.004 ppm which can be distinguished from the reagent background but not well quantitated.

Please contact me (772-4440) or L. J. Papa (772-2745) if you have any questions regarding the analyses. General questions on blood sampling can be directed to J. W. Raines or L. F. Percival.

Attachment
jah

Key Words:
Perfluorooctanoic Acid
Perfluorooctanoate
Blood Analysis
GC

There's a world of things we're doing something about

EXP000018
EID713832

000081

TABLE I

CONCENTRATION OF PERFLUOROOCTANOATE IN BLOOD (a)

Sample				GC Analysis		(b)
PRAI No.	Date Sampled	S.S.No.	Name	Date Analyzed	[C ₈], $\mu\text{g F/g blood}$	
81-1763	4/24/81	195-34-9680		5/6/81	n.d.	
81-1764	4/24/81	221-26-4482		5/6/81	n.d.	
81-1765	4/24/81	221-22-3017		5/6/81	n.d.	
81-1766	4/24/81	137-24-3476		5/6/81	n.d.	
81-1767	4/24/81	222-24-4159		5/6/81	n.d.	
81-1768	4/24/81	202-18-9533		5/7/81	<.007	
81-1769	4/24/81	222-30-6425		5/7/81	n.d.	
81-1770	4/24/81	222-34-1611		5/7/81	n.d.	
81-1771	4/24/81	222-38-1332		5/7/81	n.d.	
81-1772	4/24/81	190-22-3702		5/7/81	0.009	
81-1773	4/24/81	242-14-1512		5/8/81	<.007	
81-1774	4/24/81	218-32-8049		5/8/81	<.007	
81-1775	4/24/81	214-22-6295		5/8/81	0.043	
81-1776	4/24/81	236-74-5332		5/10/81	<.007	
81-1777	4/24/81	491-58-2410		5/10/81	<.007	

(a) Analysis as described in Lab Method ES-567 ("Determination of Perfluorooctanoic Acid in Blood, Gas Chromatographic Method", S. Stafford, 4/3/81.), using the packed column GC analysis with perfluoro-n-octanoic acid as calibration standard.

(b) Although the analysis is specifically for perfluorooctanoate (acid or salts), concentrations are given in ppm fluorine for comparison with the results of total organic fluorine analyses. (ppm F = 0.688 x ppm perfluorooctanoic acid) Estimated uncertainty is $\pm 10\%$ relative standard deviation. The lower limit for quantitation is 0.007 $\mu\text{g F/g}$. The detection limit is $\sim 0.004 \mu\text{g F/g}$, but concentrations in that range cannot be well quantitated and are reported as < 0.007. None detected (n.d.) is reported for samples with $[\text{C}_8] \lesssim 0.004 \text{ ppm}$, which cannot be distinguished from reagent background.



E. I. DU PONT DE NEMOURS & COMPANY
INCORPORATED

WILMINGTON, DELAWARE 19898

POLYMER PRODUCTS DEPARTMENT
EXPERIMENTAL STATION

cc: A. J. Dahl - 353
B. W. Karrh - N11400
L. J. Papa - 269
Pral File
I.C.

AR226-1579

PERSONAL AND CONFIDENTIAL

May 20, 1981

TO: DR. J. M. HEDGES - PPD, Circleville, Ohio

FROM: S. S. STAFFORD *S.S. Stafford*

ANALYSIS OF BLOOD SAMPLES FOR PERFLUOROOCTANOATE
(Job No. 810-633; PRAL Nos. 81-1936-1950; Notebook Nos. E22514, E26238)

As requested in your letter of 5/1/81 to L. J. Papa, the 16 blood samples submitted then have been analyzed for perfluorooctanoate (C₈). Results and sample identification are given in the attached table.

As noted there, the analyses were done using a gas chromatographic method specific for C₈ (Lab Method Number ES-567) but results have been reported as ppm F for comparison with total organic fluorine analyses. Precision is $\pm 10\%$ relative standard deviation over most of the concentration range, somewhat less at the lowest values. The lower limit for quantitation is 0.007 ppm F (0.01 ppm perfluorooctanoic acid), with a detection limit of ~ 0.004 ppm which can be distinguished from the reagent background but not well quantitated.

Please contact me (772-4440) or L. J. Papa (772-2745) if you have any questions regarding the analyses. General questions on blood sampling can be directed to J. W. Raines or L. F. Percival.

Attachment
jah

Key Words:

Perfluorooctanoic Acid
Perfluorooctanoate
Blood Analysis
GC

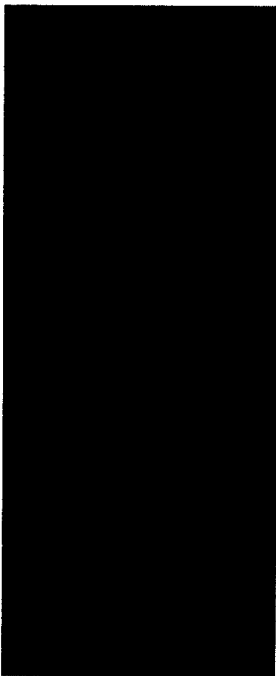
There's a world of things we're doing something about

EXP000022
EID713836

000083

TABLE I

CONCENTRATION OF PERFLUOROOCTANOATE IN BLOOD (a)

Sample				GC Analysis	
PRAL No.	Date Sampled	P.R.No.	Name	Date Analyzed	[C ₈], μ g F/g bI
81-1936	5/1/81	1485		5/12/81	0.027
81-1937	5/1/81	3518		5/12/81	n.d.
81-1938	5/1/81	2430		5/12/81	0.017
81-1939	5/1/81	3032		5/12/81	0.015
81-1940	5/1/81	5656		5/12/81	n.d.
81-1941	5/1/81	6828		5/12/81	n.d.
81-1942	5/5/81	3322		5/12/81	0.030
81-1943	5/4/81	0275		5/12/81	0.027
81-1944	5/4/81	2718		5/12/81	0.025
81-1945	5/4/81	3766		5/12/81	n.d.
81-1946	5/4/81	4748		5/13/81	0.017
81-1947	5/4/81	5646		5/13/81	0.026
81-1948	5/4/81	7350		5/13/81	0.028
81-1949	5/5/81	1505		5/13/81	0.008
81-1950	5/5/81	5645		5/13/81	0.053
81-1951	5/5/81	7452		5/13/81	0.017

- (a) Analysis as described in Lab Method ES-567 ("Determination of perfluorooctanoic Acid in Blood, Gas Chromatographic Method", S. Stafford, 4/3/81), using the packed column GC analysis with perfluoro-n-octanoic acid as calibration standard.
- (b) Although the analysis is specifically for perfluorooctanoate (acid or salts), concentrations are given in ppm fluorine for comparison with the results of total organic fluorine analyses. ($\text{ppm F} = 0.688 \times \text{ppm perfluorooctanoic acid}$) Estimated uncertainty is $\pm 10\%$ relative standard deviation. The lower limit for quantitation is $0.007 \mu\text{gF/g}$. The detection limit is $\sim 0.004 \mu\text{gF/g}$, but concentrations in that range cannot be well quantitated and are reported as < 0.007 . None detected (n.d.) is reported for samples with $[\text{C}_8] \lesssim 0.004 \text{ ppm}$, which cannot be distinguished from reagent background.

REDACTED

EXP000023
EID713837

000084